# User Guide

***i900 Mobility***

Revision 0 (July 2025)

Document No.: UG-00008

# CONTENTS

|     |             |
|-----|-------------|
| 02  | English     |
| 22  | Български   |
| 42  | 中文(简体)      |
| 62  | 中文(繁體)      |
| 82  | Hrvatski    |
| 102 | čeština     |
| 122 | Dansk       |
| 142 | Dutch       |
| 162 | Eesti keel  |
| 182 | Suomi       |
| 202 | Français    |
| 222 | Deutsch     |
| 242 | Ελληνικά    |
| 262 | magyar      |
| 282 | íslenska    |
| 302 | Italiano    |
| 322 | 한국어         |
| 342 | Latviešu    |
| 362 | Lietuvis    |
| 382 | Norsk       |
| 402 | Português   |
| 422 | Română      |
| 442 | русский     |
| 462 | Slovenčina  |
| 482 | Slovenščina |
| 502 | Español     |
| 522 | Svenska     |
| 542 | Türk        |
| 562 | українська  |

English

|                                                                  |    |
|------------------------------------------------------------------|----|
| <b>About This Guide</b>                                          | 4  |
| <b>1 Introduction and Overview</b>                               | 4  |
| 1.1 Intended Use                                                 | 4  |
| 1.2 Indication for Use                                           | 4  |
| 1.3 Contraindications                                            | 4  |
| 1.4 Qualifications of the Operating User                         | 4  |
| 1.5 Symbols                                                      | 4  |
| 1.6 i900 Mobility Components Overview                            | 5  |
| 1.6.1 Additional Components (Sold Separately)                    | 6  |
| 1.7 Setting Up the i900 Mobility System                          | 6  |
| 1.7.1 Basic Settings of the i900 Mobility System                 | 6  |
| 1.7.2 How to Use the Wireless Hub                                | 7  |
| 1.7.3 How to Use the Battery                                     | 7  |
| 1.7.4 How to Use the Battery Charger                             | 8  |
| 1.7.5 How to Place on the Desktop Cradle                         | 8  |
| 1.7.6 How to Install the Wall Mount Holder                       | 8  |
| <b>2 Medit Scan for Clinics Overview</b>                         | 9  |
| 2.1 Introduction                                                 | 9  |
| 2.2 Installation                                                 | 9  |
| 2.2.1 System Requirements                                        | 9  |
| 2.2.2 Software Installation Guide                                | 9  |
| 2.2.3 Medit Scan for Clinics User Guide                          | 9  |
| 2.3 Error Messages                                               | 10 |
| <b>3 Medit Link Express Overview</b>                             | 10 |
| 3.1 Introduction                                                 | 10 |
| 3.2 Installation                                                 | 10 |
| 3.2.1 System Requirements                                        | 10 |
| 3.2.2 Software Installation Guide                                | 10 |
| 3.2.3 Medit Link Express User Guide                              | 10 |
| 3.3 Error Messages                                               | 10 |
| <b>4 Maintenance</b>                                             | 11 |
| 4.1 Calibration                                                  | 11 |
| 4.1.1 How to Calibrate i900 Mobility                             | 11 |
| 4.1.1.1 How to Calibrate i900 Mobility in Medit Scan for Clinics | 11 |
| 4.1.1.2 How to Calibrate i900 Mobility in Medit Link Express     | 11 |
| 4.2 Cleaning, Disinfection, Sterilization Procedure              | 11 |
| 4.2.1 Reusable Tip                                               | 11 |
| 4.2.2 Mirror                                                     | 12 |
| 4.2.3 Handpiece                                                  | 12 |
| 4.2.4 Other Components                                           | 12 |
| 4.3 Disposal                                                     | 12 |
| 4.4 Battery Storage                                              | 13 |
| 4.5 The Battery Usage Precautions and Disposal Guide             | 13 |
| 4.6 Software Updates                                             | 13 |
| 4.6.1 Updates in Medit Scan for Clinics                          | 13 |
| 4.6.2 Updates in Medit Link Express                              | 13 |
| <b>5 Safety Guide</b>                                            | 14 |
| 5.1 System Basics                                                | 14 |
| 5.2 Proper Training                                              | 14 |
| 5.3 In Case of Equipment Failure                                 | 14 |
| 5.4 Hygiene                                                      | 15 |
| 5.5 Electrical Safety                                            | 15 |
| 5.6 Eye Safety                                                   | 15 |
| 5.7 Explosion Hazards                                            | 15 |
| 5.8 Pacemaker and ICD Interference Risk                          | 15 |
| 5.9 Cybersecurity Safety                                         | 15 |
| <b>6 Electro-Magnetic Compatibility Information</b>              | 16 |
| 6.1 Electromagnetic Emissions                                    | 16 |
| 6.2 Electromagnetic Immunity                                     | 16 |
| <b>7 Wireless Compliance Information</b>                         | 19 |
| 7.1 FCC Compliance Statement                                     | 19 |
| 7.2 IC Compliance Statement                                      | 19 |
| 7.3 KC Compliance Statement                                      | 19 |
| 7.4 TELEC (Japan) Compliance Statement                           | 19 |
| <b>8 Specifications</b>                                          | 20 |

About This Guide

Convention in This Guide

This guide uses various symbols to highlight important information to ensure correct usage, prevent injury to the user and others, and prevent property damage. The meanings of the symbols used are described below.

 **WARNING**

The WARNING symbol indicates information that, if ignored, could result in a medium risk of personal injury.

 **CAUTION**

The CAUTION symbol indicates safety information that, if ignored, could result in a slight risk of personal injury, property damage, or damage to the system.

 **TIPS**

The TIPS symbol indicates hints, tips, and additional information for optimal operation of the system.

1. Introduction and Overview

1.1 Intended Use

The i900 Mobility system is an intraoral 3D scanner intended to record topographical characteristics of teeth and surrounding tissues digitally. The i900 Mobility system produces 3D scans for use in computer-assisted design and manufacturing of dental restorations.

1.2 Indication for Use

The i900 Mobility system is for scanning the patient's intraoral features. Various factors (intraoral environment, operator's expertise, and laboratory workflow) may affect the final scan results when using the i900 Mobility system.

1.3 Contraindications

The i900 Mobility system is not intended to be used to create images of the internal structure of teeth or the supporting skeletal structure.

1.4 Qualifications of the Operating User

 **CAUTION**

- The i900 Mobility system is designed for use by individuals with professional knowledge in dentistry and dental laboratory technology.
- The user of the i900 Mobility system is solely responsible for determining whether or not this device is suitable for a particular patient's case and circumstances.
- The user is solely responsible for the accuracy, completeness, and adequacy of all data entered the i900 Mobility system and the provided software. The user should check the accuracy of the results and assess each individual case.
- The i900 Mobility system must be used in accordance with its accompanying User Guide.
- Improper use or handling of the i900 Mobility system will void its warranty. If you require additional information on the proper use of the i900 Mobility system, please contact your local distributor.
- The user is not allowed to modify the i900 Mobility system.

1.5 Symbols

| No. | Symbol                                                                              | Description                                         |
|-----|-------------------------------------------------------------------------------------|-----------------------------------------------------|
| 1   |  | Serial number                                       |
| 2   |  | Medical device                                      |
| 3   |  | Date of manufacture                                 |
| 4   |  | Manufacturer                                        |
| 5   |  | Caution                                             |
| 6   |  | Warning                                             |
| 7   |  | Read the user guide                                 |
| 8   |  | The official mark of the Europe Certificate         |
| 9   |  | Authorized representative in the European community |
| 10  |  | BF type of applied part                             |
| 11  |  | WEEE mark                                           |
| 12  |  | Prescription use (U.S.A)                            |
| 13  |  | Nemko NRTL/C Mark                                   |
| 14  |  | AC                                                  |
| 15  |  | DC                                                  |
| 16  |  | Temperature limitation                              |
| 17  |  | Humidity limitation                                 |
| 18  |  | Atmospheric pressure limitation                     |

|    |                                                                                   |                                                 |
|----|-----------------------------------------------------------------------------------|-------------------------------------------------|
| 19 |  | Fragile                                         |
| 20 |  | Keep dry                                        |
| 21 |  | This way up                                     |
| 22 |  | Fourteen-layer stacking prohibited              |
| 23 |  | Consult instructions for use                    |
| 24 |  | The official mark of the United Kingdom         |
| 25 |  | Authorized representative in Switzerland        |
| 26 |  | Authorized representative in the United Kingdom |
| 27 |  | Model number                                    |
| 28 |  | Quantity                                        |
| 29 |  | Unique device identifier                        |
| 30 |  | Non-sterile                                     |
| 31 |  | Stand-by                                        |

#### 1.6 i900 Mobility Components Overview

| No. | Item (Model Name)                                                                                                                                    | Qty   | Appearance                                                                            |
|-----|------------------------------------------------------------------------------------------------------------------------------------------------------|-------|---------------------------------------------------------------------------------------|
| 1   | i900 Mobility Handpiece                                                                                                                              | 1ea   |    |
| 2   | i900 Mobility Handpiece Cover (MO1-HC1)                                                                                                              | 1ea   |    |
| 3   | Wireless Hub (MO3-WFH)                                                                                                                               | 1ea   |    |
| 4   | Wall Mount Holder (MO3-WH1)                                                                                                                          | 1ea   |   |
| 5   | Power Delivery Cable Package (2 m/30 cm) (MO3-PD4)<br>• Power Delivery Cable (30 cm) – 1ea (MO3-PD3)<br>• Power Delivery Cable (2 m) – 1ea (MO1-PD2) | 1 set |  |
| 6   | Rechargeable Battery Package (MO3-REBP)<br>• Rechargeable Battery – 3ea (MO3-REB)                                                                    | 1 set |  |
| 7   | Medit Power Adapter 7 Package (MO3-MPAP)<br>• Medit Power Adapter (CDM045A3-W200V)<br>• 7 types of connectors (MA1-MPA)                              | 1 set |  |
| 8   | Battery Charger (MO3-BC1)                                                                                                                            | 1ea   |  |
| 9   | Charging Cable (2 m) (MA1-CC1)                                                                                                                       | 1ea   |  |
| 10  | Reusable Tip Package (MO1-RTP)<br>• Reusable (Large) Tip – 2ea (MO1-RTL)<br>• Reusable (Medium) Tip – 2ea (MO1-RTM)                                  | 1 set |  |
| 11  | Calibration Tool (MO1-CT1)                                                                                                                           | 1ea   |  |
| 12  | Practice Model                                                                                                                                       | 1ea   |  |
| 13  | USB Flash Drive<br>(Installer of Medit Scan for Clinics included)                                                                                    | 1ea   |  |
| 14  | User Guide                                                                                                                                           | 1ea   |  |

1.6.1 Additional Components (Sold Separately)

| No. | Item (Model Name)                                                        | Qty   | Appearance |
|-----|--------------------------------------------------------------------------|-------|------------|
| 1   | Reusable (Large) Tip Package (MO1-RTL)<br>▪ Reusable (Large) Tip – 4ea   | 1 set |            |
| 2   | Reusable (Medium) Tip Package (MO1-RTM)<br>▪ Reusable (Medium) Tip – 4ea | 1 set |            |
| 3   | Reusable (Small) Tip Package (MO1-RTS)<br>▪ Reusable (Small) Tip – 4ea   | 1 set |            |

- All components listed in 1.6 & 1.6.1 can be purchased separately.
- Availability of items for sale may vary depending on the medical device registration status in each country or region. Please contact Medit or your local distributor to check the availability of specific items.


CAUTION

- Keep the practice model in a cool place away from direct sunlight. A discolored practice model may affect the results of the practice mode.
- Medit Scan for Clinics is included in the USB drive. This product is optimized for PC, and using other devices is not recommended. Do not use anything other than a USB port. It may cause malfunction or fire.
- Medit Link Express is available for download on the App Store. This product is optimized for iPad, and use on other devices is not recommended.

1.7 Setting Up the i900 Mobility System

1.7.1 Basic Settings of the i900 Mobility System

Connect the i900 Mobility Wireless Hub

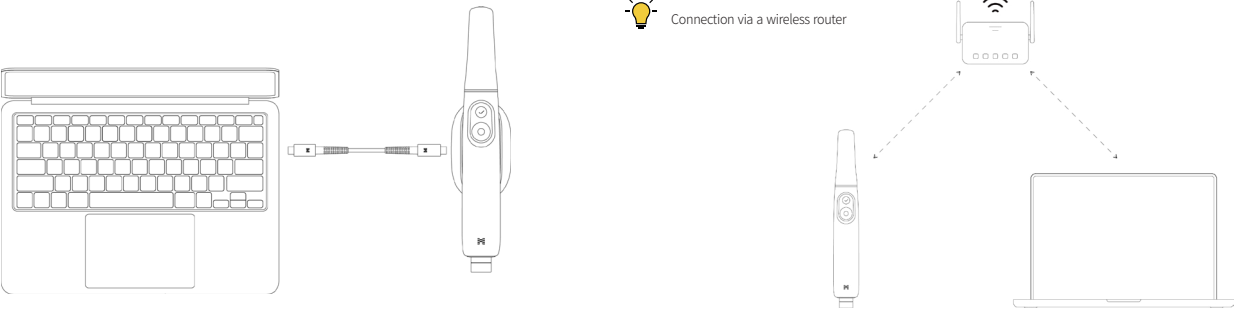- ① Connect either the 30 cm or 2 m Power Delivery Cable to the wireless hub.
- ② Connect the other end of the Power Delivery Cable to a PC.


Even if the i900 Mobility Wireless Hub is not connected to a PC, the scanner can still connect to the software via a wireless router for both Medit Scan for Clinics and Medit Link Express.

CAUTION

Medit Link Express is available exclusively for iPadOS.

Turning On the i900 Mobility

- ① Insert the battery into the i900 Mobility handpiece and press the power button once.
- ② When power is applied, the rear LED lights up.
- ③ The four LEDs on the battery's fuel gauge indicate the remaining battery level:

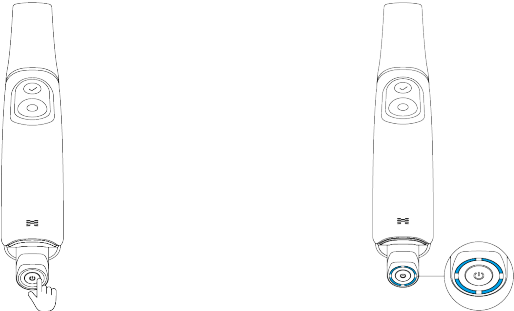
- 4 LEDs lit: 75 – 100%
- 3 LEDs lit: 50 – 74%
- 2 LEDs lit: 25 – 49%
- 1 LED lit: 12 – 24%
- 1 blinking LED: 1 – 12%

\* The indicator blinks when the charge is between 1% and 12% while the battery is inserted in the handpiece. When it's checked separately or placed in the charger, however, the indicator remains steady, with a single LED indicating the 1–12% charge range.

- No LEDs lit: 0%


Turning Off the i900 Mobility

Press and hold the power button until the rear LED is off. Removing the rechargeable battery without turning off the product may shorten the lifespan of both the handpiece and the battery.

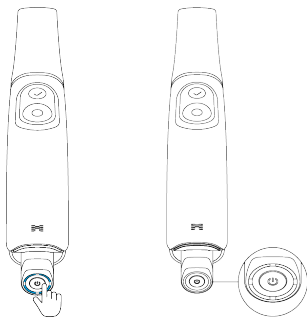1.7.2 How to Use the Wireless Hub

i900 Mobility is a wireless device that operates through a wireless module. For this reason, the i900 Mobility handpiece has a transmitter and the wireless hub has a receiver. The i900 Mobility system uses the Wi-Fi 6E frequency to transmit data and control the handpiece. Power is supplied by connecting the wireless hub to the PC using the power delivery cable. The power turns off when the cable is disconnected.


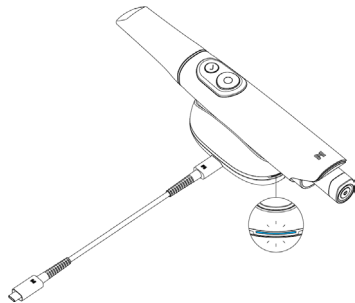
When the handpiece and wireless hub are turned on, they are ready to pair. The wireless connection is available only in Medit Scan for Clinics or Medit Link Express. An initial connection process is required, after which they will automatically connect whenever they are used again, unless the pairing is reset.

- i900 Mobility uses two antenna modules: 60 GHz for data transmission and 2.4 GHz for the control. The actual frequency depends on the local regulations.
- The actual operational range is up to 5 m, and it may vary depending on the environment.
- 5 GHz Frequency: 5210 – 5855 MHz
- 6 GHz Frequency: 6025 – 6985 MHz
- 2.4 GHz Frequency: 2402 – 2480 MHz

1.7.3 How to Use the Battery

- Rechargeable Battery
  - » Li-ion
  - » 3.6 V, 3150 mAh, 11.34 Wh
  - » Insert the battery into the i900 Mobility handpiece with the terminal aligned correctly.



- » Battery life lasts up to 1 hour. When the battery is in standalone mode (separated from the handpiece), a single click of the power button on the bottom of the battery displays the battery level via the LED.

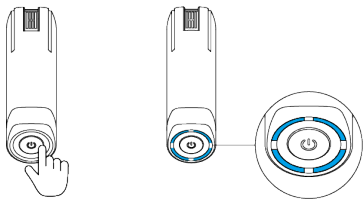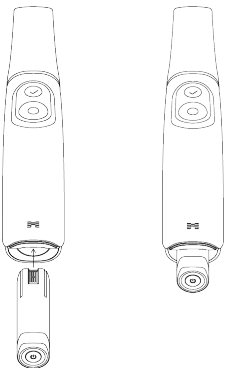Battery life may vary depending on the user environment and the battery cycle count.

1.7.4      **How to Use the Battery Charger**

- Connect the Medit Power Adapter to the battery charger's power port using the charging cable. The power turns off when the charging cable is disconnected.
- Insert the rechargeable battery into the battery charger with the charging terminals oriented correctly.

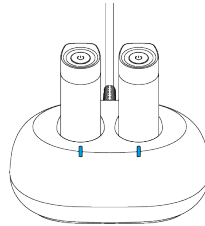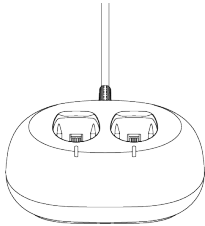
- The battery charger is only for rechargeable batteries. It takes up to 2 hours to be fully charged and the charging time may vary depending on the user environment and the battery cycle count.

 The LED on the charger lights up white when power is supplied without the battery inserted. It blinks blue while the battery is charging and remains solid blue when fully charged.


 If the battery is not correctly inserted into the battery charger, the LED on the charger will blink in red. In this case, remove the battery from the charger, check both terminals of the battery and the charger for foreign objects, wipe them gently with a soft cloth, and then insert the battery again.

1.7.5      **How to Place on the Desktop Cradle**

The wireless hub is designed to resemble a desktop cradle and functions as one.

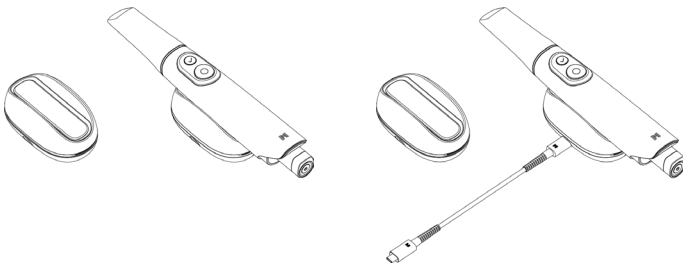1.7.6      **How to Install the Wall Mount Holder**

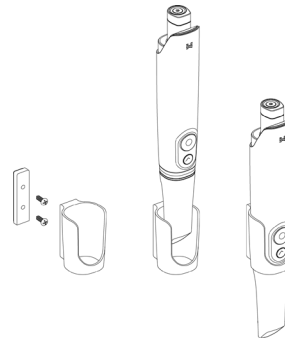2. Medit Scan for Clinics Overview

2.1 Introduction

Medit Scan for Clinics provides a user-friendly working interface to digitally record topographical characteristics of teeth and surrounding tissues using the i900 Mobility system.

2.2 Installation

2.2.1 System Requirements

Recommended System Requirements

| Windows |                                                                                                                                                                           |                                                                                            | macOS                                                                                                                                                                       |
|---------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Laptop  |                                                                                                                                                                           | Desktop                                                                                    | Laptop/Desktop                                                                                                                                                              |
| CPU     | Intel Core Ultra 7 155H<br>Intel Core i7 – 13700H<br>Intel Core i7 – 12700H<br>AMD Ryzen 7 7735H<br>AMD Ryzen 7 6800H                                                     | Intel Core i7 – 13700K<br>Intel Core i7 – 12700K<br>AMD Ryzen 7 7700X<br>AMD Ryzen 7 5800X | M1 Pro (10-core CPU, 16-core GPU)<br>M2 (8-core CPU, 10-core GPU)<br>M2 Pro (10-core CPU, 16-core GPU)<br>M3 (8-core CPU, 10-core GPU)<br>M3 Pro (11-core CPU, 14-core GPU) |
| RAM     | 32 GB                                                                                                                                                                     |                                                                                            | 24 GB                                                                                                                                                                       |
| Graphic | NVIDIA GeForce RTX 4060 (VRAM 8 GB or higher)<br>NVIDIA GeForce RTX 3070 (VRAM 8 GB or higher)<br>NVIDIA RTX A3000 (VRAM 8 GB or higher)<br>*AMD Radeon is not supported. |                                                                                            |                                                                                                                                                                             |
| OS      | Windows 10 64-bit<br>Windows 11 (recommended for Intel Core 12th Gen or newer processors)                                                                                 |                                                                                            | Monterey 12<br>Ventura 13<br>Sonoma 14<br>Sequoia 15                                                                                                                        |

Minimum System Requirements

| Windows |                                                                                                                                                                           |                                                                                        | macOS                                                      |
|---------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------------------------------------------------------------------------------|------------------------------------------------------------|
| Laptop  |                                                                                                                                                                           | Desktop                                                                                | Laptop/Desktop                                             |
| CPU     | Intel Core i5 – 13500H<br>Intel Core i5 – 12500H<br>AMD Ryzen 5 7735HS<br>AMD Ryzen 5 6600H                                                                               | Intel Core i5 – 13400<br>Intel Core i5 – 12400<br>AMD Ryzen 5 7500<br>AMD Ryzen 5 5600 | M1 (8-core CPU, 7-core GPU)<br>M2 (8-core CPU, 8-core GPU) |
| RAM     | 16 GB                                                                                                                                                                     |                                                                                        | 16 GB                                                      |
| Graphic | NVIDIA GeForce RTX 4050 (VRAM 6 GB or higher)<br>NVIDIA GeForce RTX 3060 (VRAM 6 GB or higher)<br>NVIDIA RTX A2000 (VRAM 6 GB or higher)<br>*AMD Radeon is not supported. |                                                                                        |                                                            |
| OS      | Windows 10 64-bit<br>Windows 11 (recommended for Intel Core 12th Gen or newer processors)                                                                                 |                                                                                        | Monterey 12<br>Ventura 13                                  |

Router Requirements

|          |                                                                 |
|----------|-----------------------------------------------------------------|
| Standard | Wi-Fi 6 (802.11ax) To use the 6 GHz band, Wi-Fi 6E is required. |
|----------|-----------------------------------------------------------------|

For accurate and up-to-date system requirements, please visit [www.meditlink.com](http://www.meditlink.com).

Use a PC and monitor certified to IEC 62368-1 (or IEC 60950-1), IEC 55032, and IEC 55024.

The device may not work when using cables other than the USB 3.0 cable provided by Medit. Medit is not responsible for any problems caused by cables other than the USB 3.0 cable provided by Medit. Be sure to use only the USB 3.0 cable included in the package.

2.2.2 Software Installation Guide

- ① Connect the enclosed USB flash drive to a PC.
- ② Run the installation file.
- ③ Select the setup language and click "Next."
- ④ Choose the installation path.
- ⑤ Read the "License Agreement" carefully, check "I agree to the License terms and conditions," and then click "Install."
- ⑥ The installation process may take several minutes. Please do not shut down the PC until the installation is complete.
- ⑦ After the installation is complete, restart the PC to ensure optimal program operation.

The installation will not be processed while the i900 Mobility system is connected to a PC. Please turn off the scanner before starting the installation.

2.2.3 Medit Scan for Clinics User Guide

Please refer to the User Guide of Medit Scan for Clinics: Medit Scan for Clinics > Menu > User Guide.

2.3 Error Messages

The i900 Mobility scanner system presents error messages, enabling users to intuitively identify hardware and system conditions. Some messages include solutions to aid users in self-diagnosis and problem resolution. Error messages are presented in plain text and supplemented with images when needed to enhance comprehension.

**Connection failed.**

Failed to connect to the cameras. Please follow the instructions below and try again.

- Press the "Retry" button. Device will automatically reboot and try to connect again.
- Turn off the device and disconnect all cables from the PC. Connect cables again and restart the device.
- If connection problem persists, please install the program again.
- If connection problem persists, please contact your service provider.

Cancel

Retry

**Turn Off the External Light**

 **External light was detected while scanning.**

This scan might have proceeded with an external light source, which would increase the chances of noisy and inaccurate scanned data.

 **We recommend turning the external light off before the next scan.**

Ignore

Close

If the instructions provided in the error message do not address the issue, please contact the local distributor from whom you acquired your scanner system or support@medit.com.

3. Medit Link Express Overview

3.1 Introduction

Medit Link Express offers a user-friendly interface for digitally recording the topographical characteristics of teeth and surrounding tissues using the i900 Mobility system.

3.2 Installation

3.2.1 System Requirements

| Recommended System Requirements |                               | Minimum System Requirements |                              |
|---------------------------------|-------------------------------|-----------------------------|------------------------------|
| Category                        | Details                       | Category                    | Details                      |
| Model                           | iPad Pro (7th Generation)     | Model                       | iPad Pro (6th Generation)    |
| Processor                       | M4 (10-core CPU, 10-core GPU) | Processor                   | M2 (8-core CPU, 10-core GPU) |
| RAM                             | 16 GB                         | RAM                         | 16 GB                        |
| OS                              | iPadOS 18.1.1                 | OS                          | iPadOS 18.1.1                |
| Storage                         | 1 TB                          | Storage                     | 1 TB                         |

Router Requirements

|          |                                                                 |
|----------|-----------------------------------------------------------------|
| Standard | Wi-Fi 6 (802.11ax) To use the 6 GHz band, Wi-Fi 6E is required. |
|----------|-----------------------------------------------------------------|

 For accurate and up-to-date system requirements, please visit [www.meditlink.com](http://www.meditlink.com).

 Use an iPad, and a monitor that is certified to IEC 62368-1 (or IEC 60950-1), IEC 55032, and IEC 55024.

3.2.2 Software Installation Guide

① Open the Apple App Store and search for Medit Link Express.

② Download the application.

③ Run the application.

3.2.3 Medit Link Express User Guide

Please refer to the User Guide of Medit Link Express: Medit Link Express > Settings > User Guide.

3.3 Error Messages

The i900 Mobility scanner system presents error messages, enabling users to intuitively identify hardware and system conditions. Some messages include solutions to aid users in self-diagnosis and problem resolution. Error messages are presented in plain text and supplemented with images when needed to enhance comprehension.

If the instructions provided in the error message do not address the issue, please contact the local distributor from whom you acquired your scanner system or support@medit.com.

10 User Guide

## 4. Maintenance

### CAUTION

- Equipment maintenance should only be carried out by a Medit employee or Medit-Certified company or personnel.
- In general, users are not required to perform maintenance work on the i900 Mobility system besides calibration, cleaning, and sterilization. Preventive inspections and other regular maintenance are not required.

### 4.1 Calibration

Periodic calibration is required to produce precise 3D models. You should perform calibration when:

- The quality of the 3D model is not reliable or accurate when compared to previous results.
- Environmental conditions such as temperature have changed.
- The calibration period has expired.
- You can set the calibration period in the settings:
  - Medit Scan for Clinics: Menu > Settings > Calibration Period (Days)
  - Medit Link Express: Settings > Scanner Information > Calibration Period

The calibration panel is a delicate component.  
Do not touch the panel directly. Check the calibration panel if the calibration process is not performed properly. If the calibration panel is contaminated, please contact your service provider.

If the calibration tool is long exposed to low temperatures below Medit's recommendations mentioned in the use and storage guidance, the rotation motion of the device may be compromised. In that case, forcibly turning the device may cause damage. To avoid it, keep the calibration device in the recommended temperature environment before use.

We recommend performing a calibration periodically.

- You can set the calibration period in the settings, with a default setting of 90 days:
  - Medit Scan for Clinics: Menu > Settings > Calibration Period (Days)
  - Medit Link Express: Settings > Scanner Information > Calibration Period

#### 4.1.1 How to Calibrate i900 Mobility

##### 4.1.1.1 How to Calibrate i900 Mobility in Medit Scan for Clinics

- Turn on i900 Mobility and launch Medit Scan for Clinics.
- Run the Calibration Wizard at the bottom of the main toolbar panel in Medit Scan for Clinics.
- Prepare the handpiece and calibration tool.
- Turn the dial of the calibration tool to the starting position.
- Put the handpiece into the calibration tool.
- Click "Next" to start the calibration process.
- If the calibration tool is mounted properly in the correct position, the system automatically acquires data.
- When data acquisition is complete at the starting position, turn the dial to the next position.
- Repeat the steps to the last position.
- When data acquisition is complete at the last position, the system automatically calculates and shows the calibration results.

##### 4.1.1.2 How to Calibrate i900 Mobility in Medit Link Express

- Turn on i900 Mobility and launch Medit Link Express.
- Navigate to Medit Link Express > Settings > Scanner Information > Calibration, then run the Calibration Wizard.
- Prepare the handpiece and calibration tool.
- Turn the dial of the calibration tool to the starting position.
- Put the handpiece into the calibration tool.
- Click "Calibrate Now" to start the calibration process.
- If the calibration tool is mounted properly in the correct position, the system automatically acquires data.
- When data acquisition is complete at the starting position, turn the dial to the next position.
- Repeat the steps to the last position.
- When data acquisition is complete at the last position, the system automatically calculates and shows the calibration results.

### 4.2 Cleaning, Disinfection, Sterilization Procedure

#### 4.2.1 Reusable Tip

The reusable tip is the part that is positioned in the mouth of patients during scanning, and is reusable for a limited number of times. The tip needs to be cleaned and sterilized between patient uses to avoid cross-contamination.

##### Cleaning & Disinfection

- Prepare a cleaning solution.
  - Dilute a general neutral detergent at a 1:100 ratio before use.
- Clean the reusable tips with a cleaning solution and a brush.
  - Ensure the mirror of the tip is completely clean and stain-free. If the mirror appears stained or foggy, repeat the cleaning process.

### CAUTION

- The reusable tip has a complex structure, and automatic cleaning may not achieve complete cleaning; therefore, do not clean the reusable tip in an automatic washer.
- Rinse the reusable tips three times with purified water.
- Remove moisture with a paper towel and allow them to air-dry completely at room temperature for at least 80 minutes.
- Disinfect the reusable tips using a disinfectant containing 15% or less Isopropyl Alcohol (IPA) for 1 minute. Then, ensure they are thoroughly dried at room temperature for at least 5 minutes.
  - Before using a disinfectant product, please refer to the product user manual for proper use.
  - You can find the list of recommended disinfectants in the Medit Help Center at <http://support.medit.com/hc>.

## Sterilization

- The tip should be cleaned manually using a disinfecting solution. After cleaning and disinfecting, inspect the mirror inside the tip to ensure there aren't any stains or smudges.
- Repeat the cleaning and disinfection process if necessary. Carefully dry the mirror with a paper towel.
- Insert the tip into a paper sterilization pouch and seal it, making sure that it is airtight. Use either a self-adhesive or heat-sealed pouch.
- Sterilize the wrapped tip in an autoclave with the following conditions:
  - » Sterilize in a gravity autoclave at 135°C (275°F) for 10 minutes and dry for 30 minutes.
  - » Sterilize in a pre-vacuum autoclave at 134°C (273.2°F) for 4 minutes and dry for 20 minutes.
- Use an autoclave program that dries the wrapped tip before opening the autoclave.
- Scanner tips can be re-sterilized up to 150 times. After reaching this limit, they must be disposed of according to the guidelines in the disposal section.
- Autoclave times and temperatures may vary depending on the autoclave type and manufacturer. For this reason, it may not be able to meet the maximum number of times. Please refer to the user manual of the autoclave manufacturer you are using to determine whether the required conditions are met.

## CAUTION

- The mirror of the scanner tip is a delicate optical component that should be handled with care to ensure optimal scan quality. Be careful not to scratch or smudge it as any damage or blemishes may affect the data acquired.
- Make sure to always wrap the tip before autoclaving. If you autoclave an exposed tip, it will cause stains on the mirror, which cannot be removed. Check the autoclave manual for more information.
- Tips that have been cleaned, disinfected, and sterilized must remain sterile until they are used on the patient.
- Medit is not responsible for any damage, including tip distortion, caused by cleaning, disinfection, or sterilization procedures that do not comply with the guidelines outlined above.

## 4.2.2 Mirror

The presence of impurities or smudges on the tip mirror may lead to poor scan quality and an overall poor scanning experience. In such situations, clean the mirror following the steps below:

- ① Disconnect the scanner tip from the i900 Mobility handpiece.
- ② Pour alcohol on a clean cloth or cotton-tipped swab and wipe the mirror. Make sure to use alcohol that is free of impurities or it may stain the mirror. You can use either ethanol or propanol (ethyl/propyl alcohol).
- ③ Wipe the mirror dry using a dry, lint-free cloth.
- ④ Make sure the mirror is free of dust and fibers. Repeat the cleaning process as necessary.

## 4.2.3 Handpiece

After treatment, clean and disinfect all other surfaces of the i900 Mobility handpiece except for the scanner front (optical window) and end (air vent hole). Cleaning and disinfecting must be done with the device turned off. Use the device only after it is completely dry.

The recommended cleaning and disinfecting solution is denatured alcohol (ethyl alcohol or ethanol) – typically 60 – 70% Alc/Vol.

The general cleaning and disinfecting procedures are as follows:

- ① Turn off the device using the power button.
- ② Unplug all the cables from the wireless hub.
- ③ Clean the filter on the front end of the i900 Mobility handpiece.
  - » If alcohol is poured directly into the filter, it may seep inside the i900 Mobility handpiece and cause a malfunction.
  - » Do not clean the filter by pouring alcohol or cleaning solution directly into the filter. The filter must be gently wiped with a cotton or soft cloth moistened with alcohol.
  - » Do not wipe by hand or apply excessive force.
  - » Medit is not responsible for any damage or malfunction that occurs during cleaning that does not follow the guidelines above.
- ④ After cleaning the filter, put the cover on the front of the i900 Mobility handpiece.
- ⑤ Pour the disinfectant onto a soft, lint-free, and non-abrasive cloth.
- ⑥ Wipe the scanner surface with the cloth.
- ⑦ Dry the surface with a clean, dry, lint-free, and non-abrasive cloth.

## CAUTION

- Do not clean the i900 Mobility handpiece when the device is turned on as the fluid may enter the scanner and cause malfunction.
- Use the device after it is completely dry.
- Chemical cracks may appear if improper cleaning and disinfecting solutions are used during cleaning.

## 4.2.4 Other Components

- Pour the cleaning and disinfecting solution onto a soft, lint-free, and non-abrasive cloth.
- Wipe the component surface with the cloth.
- Dry the surface with a clean, dry, lint-free, and non-abrasive cloth.

## CAUTION

- Chemical cracks may appear if improper cleaning and disinfecting solutions are used during cleaning.

## 4.3 Disposal

## CAUTION

- The scanner tip must be sterilized before disposal. Sterilize the tip as described in the section "4.2.1 Reusable Tip."
- Dispose of the scanner tip as you would any other clinical waste.
- Other components are designed to conform with the following directives:
  - RoHS, Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment. (2011/65/EU)
  - WEEE, Waste Electrical and Electronic Equipment Directive. (2012/19/EU)

4.4 Battery Storage

- Put it in a pack or a box and store indoors in a cool environment out of direct sunlight.
- Store the battery in a dry place with an ambient temperature of from -20°C to +30°C (-4°F to +86°F).
- If not used for a long time, the battery self-discharge may accelerate and slip into sleep mode. To minimize the effect of deactivation, store the battery pack in between +10°C - +30°C (+50°F to +86°F).
- When charging for the first time after long-term storage, the battery capacity may decrease due to pack deactivation. Restore the battery through several full charge/discharge cycles.
- If the battery is stored for more than 6 months, it should be charged at least once every 6 months to prevent the decrease in shelf life due to self-discharge.

CAUTIONS: Battery Safety Specifications

| Safety Specifications                 |                   |                     |
|---------------------------------------|-------------------|---------------------|
| Overcharge                            | Sense Voltage     | 4.250 V ± 0.020     |
|                                       | Steady Voltage    | 4.050 V ± 0.05      |
|                                       | Sense Delay       | 1.0 s ± 0.3         |
| Over-Discharge                        | Sense Voltage     | 2.50 V ± 0.050      |
|                                       | Steady Voltage    | 2.90 V ± 0.10       |
|                                       | Sense Delay       | 128 ms ± 38.4       |
| Overcurrent                           | Sense (Charge)    | 12.0 A + 5.0 / -4.5 |
|                                       | Sense Delay       | 16.0 ms ± 4.8       |
|                                       | Sense (Discharge) | Min. 7.5 A          |
|                                       | Sense Delay       | 16.0 ms ± 4.8       |
| Current consumption at operating mode |                   | Max. 150.0 µA       |

The possible safety specifications are determined by the propulsion control module (PCM) in the parts list.

4.5 The Battery Usage Precautions and Disposal Guide

CAUTION

- Be sure to fully understand how to replace the battery before use.
- Use a charger suitable for the specified voltage and current.
- Do not attempt to reverse the charge. Reverse charging may increase the gas pressure inside the battery and cause the battery to leak.
- Do not attempt to recharge a fully charged battery. Repeated overcharging could result in battery performance degradation and overheating.
- Charging efficiency decreases at temperatures below 0°C (32°F) and above +45°C (113°F).
- Do not short-circuit the positive (+) and negative (-) terminals with metallic objects such as metal wires, necklaces, or chains.
- To avoid malfunction or damage, do not drop or throw the battery.
- Do not deform a battery by excessive pressure.
- Do not directly solder anything onto the battery pack.
- Do not allow children to change batteries without adult supervision.
- Do not dispose of the battery as general waste and separate it from recyclables.
- Do not dispose of or throw the battery into a fire. Heat can cause a battery explosion and a fire.
- Separate them from each other when disposing of the secondary batteries with different electrochemical systems.
- Dispose of the battery by discharging it to prevent heat from short-circuiting.
- Battery disposal methods may vary by country and region. Dispose of used batteries under local laws and regulations.

4.6 Software Updates

4.6.1 Updates in Medit Scan for Clinics

Medit Scan for Clinics automatically checks for updates when the software is in operation.  
If a new version of the software is released, it will be downloaded automatically.

4.6.2 Updates in Medit Link Express

Medit Link Express automatically checks for updates when the software is in operation.  
If a new version of the software is released, it will be downloaded automatically, provided that automatic updates are enabled in the iPad settings.

## 5. Safety Guide

Please adhere to all the safety procedures as detailed in this User Guide to prevent human injury and equipment damage. This document uses the words WARNING and CAUTION when highlighting precautionary messages.

Carefully read and understand the guidelines, including all preventive messages as prefaced by the words WARNING and CAUTION. To avoid bodily injury or equipment damage, make sure to adhere strictly to the safety guidelines. All instructions and precautions as specified in the Safety Guide must be observed to ensure the system's proper functionality and personal safety.

The i900 Mobility system should only be operated by dental professionals and technicians who are trained to use the system. Using the i900 Mobility system for any purpose other than its intended usage as outlined in the section "1.1 Intended Use" may result in injury or damage to the equipment. Please handle the i900 Mobility system according to the guidelines in the safety guide.

Any serious incident involving the device should be reported to the manufacturer and the competent authority of the member state in which the user and patients are located.

### 5.1 System Basics

The i900 Mobility system is a high-precision optical medical device. Acquaint all the following safety and operating instructions before the installation, usage, and operation of the i900 Mobility system.

#### CAUTION

- The wireless hub is designed specifically for the i900 Mobility system and should not be used with any other device.
- This device has been assessed for conformity for use in business environments. When used in residential environments, there is a risk of radio wave interference.
- If the product has been stored in a cold environment, give it time to adjust to the temperature of the environment before use.
- If used immediately, condensation may occur, which may damage the electronic parts inside the unit.
- Ensure that all components provided are free from physical damage. Safety cannot be guaranteed if there is any physical damage to the unit.
- Before using the system, check that there are no issues such as physical damage or loose parts. If there is any visible damage, do not use the product and contact the manufacturer or your local representative.
- Check the i900 Mobility handpiece and its accessories for any sharp edges.
- When not in use, the i900 Mobility system should be kept mounted on a desk stand or a wall mount stand.
- Do not install the desk stand on an inclined surface.
- Do not place any object on the i900 Mobility system.
- Do not place the i900 Mobility system on any heated or wet surface.
- Do not block the air vents located at the rear of the i900 Mobility system. If the equipment overheats, the i900 Mobility system may malfunction or stop working.
- The i900 Mobility battery can only be compatible with the i900 Mobility system.
- Do not touch the charging terminals of the rechargeable battery with your hands or other instruments.
- If the charging terminal of the rechargeable battery is damaged, do not use it and contact the manufacturer or regional manager.
- If the shape of the rechargeable battery has been deformed by dropping or physical damage, never use it, and contact the manufacturer or regional manager.
- Do not spill any liquid on the i900 Mobility system.
- The i900 Mobility handpiece and other included components are made of electronic components. Do not allow any kind of liquid or foreign objects to enter.
- Do not pull or bend the cable connected to the wireless hub.
- Carefully arrange all the cables so that you or your patient do not trip or get caught in the cables. Any pulling tension on the cables may cause damage to the i900 Mobility system.
- Always place the plug of the Medit Power Adapter of the i900 Mobility system in an easily accessible location.
- Always keep an eye on the product and your patient while using the product to check for abnormalities.
- Proceed with calibration, cleaning, disinfection, and sterilization in accordance with the contents of the user guide.
- If you drop the i900 Mobility tip on the floor, do not attempt to reuse it. Discard the tip immediately as there is a risk that the mirror attached to the tip may have been dislodged.
- Due to its fragile nature, the i900 Mobility tips should be handled with care. To prevent damage to the tip and its internal mirror, be careful to avoid contact with the patient's teeth or restorations.
- If the i900 Mobility system is dropped on the floor or if the unit is impacted, it must be calibrated before use. If the instrument is unable to connect to the software, consult the manufacturer or authorized resellers.
- If the equipment fails to operate normally, such as having issues with accuracy, stop using the product, and contact the manufacturer or authorized resellers.
- Install and use only approved programs to ensure the proper functionality of the i900 Mobility system.
- In the event of a severe accident involving the i900 Mobility system, notify the manufacturer and report it to the competent national authority of the country where the user and patient reside.
- If the PC with the software installed does not have security software or if there is a risk of malicious code intrusion into the network, the PC may be breached with malware (malicious software such as viruses or worms that damage your computer).
- The software for this product must be used in compliance with medical and personal information protection laws.
- To provide electrical insulation and maintain electrical safety, a coating is applied to insulate the device except on the areas where the USB ports are located.
- This high-precision 3D scanner is particularly sensitive to impacts. Handle with care to prevent performance degradation from shocks, including accidental drops.

### 5.2 Proper Training

#### WARNING

Before using your i900 Mobility system on patients:

- You should have been trained to use the system or read and fully understand this User Guide.
- You should be familiar with the safe use of the i900 Mobility system, as detailed in this User Guide.
- Before use or after changing any settings, the user should check that the live image is displayed properly in the camera preview window of the program.

### 5.3 In Case of Equipment Failure

#### WARNING

If your i900 Mobility system is not working properly or if you suspect that there is a problem with the equipment:

- Remove the device from the patient's mouth and discontinue use immediately.
- Disconnect the device from the PC and check for errors.
- Remove the rechargeable battery from the i900 Mobility system.
- Contact the manufacturer or authorized resellers.
- Modifications to the i900 Mobility system are prohibited by law as they may compromise the safety of the user, patient, or a third party.

## 5.4 Hygiene

### WARNING

For clean working conditions and patient safety, ALWAYS wear clean surgical gloves when:

- Handling and replacing the tip.
- Using the i900 Mobility system on patients.
- Touching the i900 Mobility system.

### WARNING

The i900 Mobility system and its optical window should be always kept clean. Before using the i900 Mobility system on a patient, be sure to:

- Sterilize the i900 Mobility system as described in the section "4.2 Cleaning, Disinfection, Sterilization Procedure."
- Use a sterilized tip.

## 5.5 Electrical Safety

### WARNING

- To prevent electric shock, the i900 Mobility system must only be connected to a power source with a protective earth connection.
- If you cannot insert the i900 Mobility-supplied plug into the main outlet, contact a qualified electrician to replace the plug or outlet. Do not try to circumvent these safety guidelines.
- Do not use a grounding-type plug connected to the i900 Mobility system for any other purpose than its intended use.
- The i900 Mobility system only uses RF energy internally. The amount of RF radiation is low and does not interfere with surrounding electromagnetic radiation.
- There is a risk of electric shock if you attempt to access the inside of the i900 Mobility system. Only qualified service personnel should access the system.
- Do not connect the i900 Mobility system to a regular power strip or extension cord, as these connections are not as safe as grounded outlets.

Failure to adhere to these safety guidelines may result in the following hazards:

- » All connected equipment's total short circuit current may exceed the limit specified in EN/IEC 60601-1.
  - » The impedance of the ground connection may exceed the limit specified in EN/IEC 60601-1.
  - Do not place liquids such as beverages near the i900 Mobility system and avoid spilling any liquid on the system.
  - Never spill liquid of any kind on the i900 Mobility system.
  - Condensation due to changes in temperature or humidity can cause moisture buildup inside the i900 Mobility system, which may damage the system.
- Before connecting the i900 Mobility system to a power supply, be sure to keep the i900 Mobility system at room temperature for at least two hours to prevent condensation.
- If condensation is visible on the product surface, the i900 Mobility system should be left at room temperature for more than 8 hours.

- You should only disconnect the i900 Mobility system from the power supply via its power cord or rechargeable battery.
- When disconnecting the power cord or the rechargeable battery, hold the surface to remove it.
- Before disconnecting, make sure to turn off the power on the device using the power switch on the handpiece.
- The EMISSIONS characteristics of this equipment make it suitable for use in industrial areas and hospitals (CISPR 11 Class A).
- If used in a residential environment (for which CISPR 11 Class B is normally required), this equipment might not offer adequate protection to radio-frequency communication services.
- Use only batteries provided for use with the i900 Mobility system. Other batteries could damage the i900 Mobility system.
- Avoid pulling on the communication cables, power cables, etc. used with the i900 Mobility system.
- Use only medical adapters provided for use with the i900 Mobility system. Other adapters could damage the i900 Mobility system.
- Do not touch the connectors of the device and the patient simultaneously.

## 5.6 Eye Safety

### WARNING

- The i900 Mobility system projects a bright light from its tip during scanning.
  - The bright light projected from the tip of i900 Mobility is not harmful to the eyes. However, you should not look directly at the bright light nor aim the light beam into the eyes of others. Generally, intense light sources can cause eyes to become brittle and the likelihood of secondary exposure is high. As with other intense light source exposure, you may experience a temporary reduction in visual acuity, pain, discomfort, or visual impairment, all of which increase the risk of secondary accidents.
  - There is an LED that emits UV-C wavelengths inside the i900 Mobility handpiece. It is irradiated only inside the i900 Mobility handpiece and does not go outside.
  - The blue light visible inside the i900 Mobility handpiece is for guidance, not UV-C light. It is harmless to the human body.
  - The UV-C LED operates with a wavelength of 270 – 285 nm.
  - Disclaimer for Risks Involving Patients with Epilepsy
- The i900 Mobility system should not be used on patients that have been diagnosed with epilepsy due to the risk of seizures and injury.
- For the same reason, dental staff who have been diagnosed with epilepsy should not operate the i900 Mobility system.

## 5.7 Explosion Hazards

### WARNING

- The i900 Mobility system is not designed to be used near flammable liquids, gases, or in environments with high oxygen concentrations.
  - There is a risk of explosion if you use the i900 Mobility system near flammable anesthetics.
  - The rechargeable battery used with the i900 Mobility system is designed with safety devices included.
  - The rechargeable battery shall not be exposed to excessive heat such as sunshine or the like. Failure to observe this could result in the battery explosion.
  - Please be careful with the storage and maintenance of the battery.
  - The battery charger is designed to adjust the current supply after charging is complete.
- But if it is not going to be used for a long time, disconnect the power to the battery charger or remove the fully charged battery from the battery charger.

## 5.8 Pacemaker and ICD Interference Risk

### WARNING

- Implantable Cardioverter Defibrillators (ICDs) and pacemakers may have interference due to some devices.
- Maintain a moderate distance from the patient's ICD or pacemaker when using the i900 Mobility system.
- For more information on peripherals used with the i900 Mobility system, check the respective manufacturer's manuals.

## 5.9 Cybersecurity Safety

- If a cybersecurity incident occurs, stop using the scanner and software immediately. Power off the scanner and log out of the software.
- Promptly report the incident to our support team through email, phone, or other available means of contact. Please refer to the last page of the User Guide for contact information.
- When reporting an incident, please provide as much information as possible, including the time of occurrence and any unusual behavior you noticed. This information will assist us in resolving the issue quickly.

# 6. Electro-Magnetic Compatibility Information

## 6.1 Electromagnetic Emissions

The i900 Mobility system is intended for use in the electromagnetic environment as specified below. The customer or the user of the i900 Mobility system should ensure that it is used in such an environment.

| Guidance and Manufacturer's Declaration – Electromagnetic Emission |            |                                                                                                                                                                                                                                              |
|--------------------------------------------------------------------|------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Emission Test                                                      | Compliance | Electromagnetic Environment – Guidance                                                                                                                                                                                                       |
| RF Emissions CISPR 11                                              | Group 1    | The i900 Mobility system is suitable for use in all establishments.<br>This includes domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes. |
| RF Emissions CISPR 11                                              | Class A    |                                                                                                                                                                                                                                              |
| Harmonic Emissions IEC 61000-3-2                                   | Class A    |                                                                                                                                                                                                                                              |
| Voltage Fluctuations / Flicker Emissions                           | Complies   |                                                                                                                                                                                                                                              |

### ⚠ WARNING

The i900 Mobility system is intended for use by healthcare professionals only. This equipment/system may cause radio interference or may disrupt the operation of nearby equipment. It may be necessary to take mitigation measures, such as re-orienting or relocating the i900 Mobility system or shielding the location.

## 6.2 Electromagnetic Immunity

### • Guidance 1

The i900 Mobility system is intended for use in the electromagnetic environment as specified below.  
The customer or the user of the i900 Mobility system should ensure that it is used in such an environment.

| Guidance and Manufacturer's Declaration – Electromagnetic Immunity                                      |                                                                                                                                 |                                                                                                                                 |                                                                                                                                                                                                                                                                                                              |
|---------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Immunity Test                                                                                           | IEC 60601 Test Level                                                                                                            | Compliance Level                                                                                                                | Electromagnetic Environment – Guidance                                                                                                                                                                                                                                                                       |
| Electrostatic Discharge (ESD)<br>IEC 61000-4-2                                                          | ± 8 kV contact ± 2 kV, ± 4 kV,<br>± 8 kV, ± 15 kV air                                                                           | ± 8 kV contact ± 2 kV, ± 4 kV,<br>± 8 kV, ± 15 kV air                                                                           | Floors should be made of wood, concrete, or ceramic tiles.<br>If floors are covered with a synthetic material, relative humidity of at least 30% is recommended.                                                                                                                                             |
| Electrical Fast Transient / Burst<br>IEC 61000-4-4                                                      | ±2 kV for power supply lines<br>±1 kV for input/output lines                                                                    | ±2 kV for power supply lines<br>±1 kV for input/output lines                                                                    | The mains power quality should be that of a typical commercial or hospital environment.                                                                                                                                                                                                                      |
| Surge IEC 61000-4-5                                                                                     | ±0.5 kV, ±1 kV differential mode                                                                                                | ±0.5 kV, ±1 kV differential mode                                                                                                | The mains power quality should be that of a typical commercial or hospital environment.                                                                                                                                                                                                                      |
|                                                                                                         | ±0.5 kV, ±1 kV, ±2 kV common mode                                                                                               | ±0.5 kV, ±1 kV, ±2 kV common mode                                                                                               |                                                                                                                                                                                                                                                                                                              |
| Voltage Dips, Short Interruptions, and Voltage Variations on Power Supply Input Lines<br>IEC 61000-4-11 | 0% Ur (100% dip in Ur) for 0.5/1 cycles<br>70% Ur (30% dip in Ur) for 25/30 cycles<br>0% Ur (100% dip in Ur) for 250/300 cycles | 0% Ur (100% dip in Ur) for 0.5/1 cycles<br>70% Ur (30% dip in Ur) for 25/30 cycles<br>0% Ur (100% dip in Ur) for 250/300 cycles | The mains power quality should be that of a typical commercial or hospital environment. If the user of the i900 Mobility system requires continued operation during power mains interruptions, it is recommended that the i900 Mobility system be powered from an uninterruptible power supply or a battery. |
| Power Frequency Magnetic Fields<br>(50/60Hz) IEC 61000-4-8                                              | 30 A/m                                                                                                                          | 30 A/m                                                                                                                          | Power frequency magnetic fields should be at levels characteristic of a location in a typical commercial or hospital environment.                                                                                                                                                                            |
| Proximity Magnetic Fields in the Frequency Range 9 kHz to 13.56 MHz<br>Immunity IEC 61000-4-39          | 8 A/m<br>30 kHz CW modulation                                                                                                   | 8 A/m<br>30 kHz CW modulation                                                                                                   | Resistance to magnetic fields was tested and applied only to surfaces of enclosures or accessories accessible during intended use.                                                                                                                                                                           |
|                                                                                                         | 65 A/m<br>134.2 kHz<br>PM 2.1 kHz<br><br>7.5 A/m<br>13.56 MHz<br>PM 50 kHz                                                      | 65 A/m<br>134.2 kHz<br>PM 2.1 kHz<br><br>7.5 A/m<br>13.56 MHz<br>PM 50 kHz                                                      |                                                                                                                                                                                                                                                                                                              |

NOTE: Ur is the main voltage (AC) prior to the application of the test level.

Guidance 2

| Recommended Separation Distances Between Portable and Mobile Communication Equipment and the i900 Mobility System |                                                                   |                                         |  |
|-------------------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------|-----------------------------------------|--|
| Rated Maximum Output Power of the Transmitter [W]                                                                 | Separation Distance According to the Frequency of Transmitter [M] |                                         |  |
|                                                                                                                   | IEC 60601-1-2:2014                                                |                                         |  |
|                                                                                                                   | 150 kHz to 80 MHz<br>$d = 1.2 \sqrt{P}$                           | 80 MHz to 2.7 GHz<br>$d = 2.0 \sqrt{P}$ |  |
| 0.01                                                                                                              | 0.12                                                              | 0.20                                    |  |
| 0.1                                                                                                               | 0.38                                                              | 0.63                                    |  |
| 1                                                                                                                 | 1.2                                                               | 2.0                                     |  |
| 10                                                                                                                | 3.8                                                               | 6.3                                     |  |
| 100                                                                                                               | 12                                                                | 20                                      |  |

For transmitters rated at a maximum output power not listed above, the recommended separation distance (d) in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

NOTE 1: At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

NOTE 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.

Guidance 3

The i900 Mobility system is intended for use in the electromagnetic environment specified below. The customer or the user of the i900 Mobility system should ensure that it is used in such an environment.

| Guidance and Manufacturer's Declaration – Electromagnetic Immunity |                                                       |                  |                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                        |
|--------------------------------------------------------------------|-------------------------------------------------------|------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Immunity Test                                                      | IEC 60601 Test Level                                  | Compliance Level | Electromagnetic Environment – Guidance                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                 |
| Conducted RF<br>IEC 61000-4-6                                      | 3 Vrms 150 kHz to 80 MHz<br>Outside ISM Bands amateur | 3 Vrms           | <p>Portable and mobile RF communications equipment should not be used closer to any part of the Ultrasound System, including cables, than the recommended separation distance. This is calculated using the equation applicable to the frequency of the transmitter.</p> <p><b>Recommended Separation Distance (d):</b><br/><math>d = 1.2 \sqrt{P}</math><br/><b>IEC 60601-1-2:2007</b><br/><math>d = 1.2 \sqrt{P}</math> 80 MHz to 800 MHz<br/><math>d = 2.3 \sqrt{P}</math> 80 MHz to 2.5 GHz<br/><b>IEC 60601-1-2:2014</b><br/><math>d = 2.0 \sqrt{P}</math> 80 MHz to 2.7 GHz</p> <p>Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer, d is the recommended separation distance in meters (m).<br/>Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey should be less than the compliance level in each frequency range.<br/>Interference may occur in the vicinity of equipment marked with the following symbol:</p> |
|                                                                    | 6 Vrms 150 kHz to 80 MHz<br>In ISM Bands amateur      | 6 Vrms           |                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                        |
| Radiated RF<br>IEC 61000-4-3                                       | 3 V/m 80 MHz to 2.7 GHz                               | 3 V/m            |                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                     |

NOTE 1: At 80 MHz and 800 MHz, the higher frequency range applies.

NOTE 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.

NOTE 3: The ISM (Industrial, Scientific, and Medical) bands between 150 kHz and 80 MHz are 6.765 MHz to 6.795 MHz; 13.553 MHz to 13.567 MHz; 26.957 MHz to 27.283 MHz; and 40.66 MHz to 40.70 MHz.

Guidance 4

The i900 Mobility system is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. Portable RF communications equipment should be used no closer than 30 cm (12 inches) to any part of the i900 Mobility system. Otherwise, degradation of the performance of this equipment could result.

| Guidance and Manufacturer's Declaration – Electromagnetic Immunity |                    |                                                                 |                                        |                      |                  |
|--------------------------------------------------------------------|--------------------|-----------------------------------------------------------------|----------------------------------------|----------------------|------------------|
| Immunity Test                                                      | Band <sup>1)</sup> | Service <sup>1)</sup>                                           | Modulation                             | IEC 60601 Test Level | Compliance Level |
| Proximity Fields from RF Wireless Communications IEC 61000-4-3     | 380 – 390 MHz      | TETRA 400                                                       | Pulse Modulation 18 Hz                 | 27 V/m               | 27 V/m           |
|                                                                    | 430 – 470 MHz      | GMRS 460; FRS 460                                               | FM $\pm 5$ kHz Deviation<br>1 kHz sine | 28 V/m               | 28 V/m           |
|                                                                    | 704 – 787 MHz      | LTE Band 13, 17                                                 | Pulse Modulation 217 Hz                | 9 V/m                | 9 V/m            |
|                                                                    | 800 – 960 MHz      | GSM 800/900; TETRA 800; iDEN 820; CDMA 850; LTE Band 5          | Pulse Modulation 18 Hz                 | 28 V/m               | 28 V/m           |
|                                                                    | 1700 – 1990 MHz    | GSM 1800; CDMA 1900; GSM 1900; DECT; LTE Band 1, 3, 4, 25; UMTS | Pulse Modulation 217 Hz                | 28 V/m               | 28 V/m           |
|                                                                    | 2400 – 2570 MHz    | Bluetooth; WLAN 802.11b/g/n; RFID 2450; LTE Band 7              | Pulse Modulation 217 Hz                | 28 V/m               | 28 V/m           |
|                                                                    | 5100 – 5800 MHz    | WLAN 802.11a/n                                                  | Pulse Modulation 217 Hz                | 9 V/m                | 9 V/m            |

NOTE: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.

 WARNING

- Use of the i900 Mobility system adjacent to or on other equipment must be avoided as it may result in improper operation. If this use is necessary, it is advisable that this and the other equipment be observed to verify that they are operating normally.
- The use of accessories, transducers, and cables other than those specified or provided by the Medit of the i900 Mobility system could result in high electromagnetic emissions or reduced electromagnetic immunity from this equipment and result in improper operation.

<sup>1)</sup> For some services, only the uplink frequencies are included.

7. Wireless Compliance Information

7.1 FCC Compliance Statement

Federal Communications Commission Interference Statement

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.  
FCC Caution: Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate this equipment.  
This device and its antenna(s) must not be co-located or operation in conjunction with any other antenna or transmitter.

|                   |                |              |              |
|-------------------|----------------|--------------|--------------|
| FCC ID:           |                |              |              |
| Intraoral Scanner | 2A2QMMO3-I900M | Wireless Hub | 2A2QMMO3-WFH |

IMPORTANT NOTE:

FCC Radiation Exposure Statement:

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. The wireless hub should be installed and operated with a minimum distance of 20 cm (8 inches) between the radiator and your body. Incidental physical contact, such as during installation or operation, is permitted. The handpiece is intended for use in direct contact with the human body.

7.2 IC Compliance Statement

This Class B digital apparatus complies with Canadian ICES-003.  
This device complies with Industry Canada license-exempt RSS standard(s).  
Operation is subject to the following two conditions: (1) this device may not cause interference, and (2) this device must accept any interference, including interference that may cause undesired operation of the device.  
Any changes or modifications not expressly approved by the manufacturer could void the user's authority to operate this equipment.  
This device and its antenna(s) must not be co-located or operation in conjunction with any other antenna or transmitter.  
The device could automatically discontinue transmission in case of the absence of information to transmit or operational failure.  
Note that this is not intended to prohibit transmission of control or signaling information or the use of repetitive codes where required by the technology.

|                   |                |              |              |
|-------------------|----------------|--------------|--------------|
| IC ID:            |                |              |              |
| Intraoral Scanner | 27675-MO3I900M | Wireless Hub | 27675-MO3WFH |

IMPORTANT NOTE:

IC Radiation Exposure Statement

This equipment complies with IC RSS-102 radiation exposure limits set forth for an uncontrolled environment. The wireless hub should be installed and operated with a minimum distance of 20 cm (8 inches) between the radiator and your body. Incidental physical contact, such as during installation or operation, is permitted. The handpiece is intended for use in direct contact with the human body.

Transmit Antenna Notice

This radio transmitter has been approved by Innovation, Science and Economic Development Canada to operate with the antenna types listed below, with the maximum permissible gain indicated. Antenna types not included in this list that have a gain greater than the maximum gain indicated for any type listed are strictly prohibited from use with this device.

Antenna List

| Model           | Type | Maximum Gain (dBi)                                |
|-----------------|------|---------------------------------------------------|
| INNO-AFP-0420   | FPCB | 2.4 GHz: 4.35 dBi 5 GHz: 4.51 dBi 6 GHz: 6.86 dBi |
| INNO-AFP-0404-R | FPCB | 2.4 GHz: 2.09 dBi 5 GHz: 3.27 dBi 6 GHz: 3.69 dBi |
| INNO-AFP-0404-L | FPCB | 2.4 GHz: 2.95 dBi 5 GHz: 4.58 dBi 6 GHz: 4.79 dBi |

7.3 KC Compliance Statement

KC ID: R-C-SX0-MO3-I900M

This device has been assessed for conformity for use in business environments. When used in residential environments, there is a risk of radio wave interference.

7.4 TELEC (Japan) Compliance Statement

218-250616  
218-250617

|                   |            |              |            |
|-------------------|------------|--------------|------------|
| MIC ID:           |            |              |            |
| Intraoral Scanner | 218-250616 | Wireless Hub | 218-250617 |

8. Specifications

|                                                                      |                                               |                                                                               |
|----------------------------------------------------------------------|-----------------------------------------------|-------------------------------------------------------------------------------|
| Model Name                                                           |                                               | MO3-i900M                                                                     |
| Trade Name                                                           |                                               | i900 Mobility                                                                 |
| Packing Unit                                                         |                                               | 1 set                                                                         |
| Classifications for Protection Against Electric Shock, Applied Parts |                                               | Class II, Type BF Applied Parts                                               |
| * This product is a medical device.                                  |                                               |                                                                               |
| Handpiece                                                            |                                               |                                                                               |
| Dimension                                                            |                                               | 209 x 42 x 42 mm ± 5% (excluding reusable tip and rechargeable battery)       |
| Weight                                                               |                                               | 242 g ± 5% (excluding reusable tip and rechargeable battery)                  |
| Power Supply                                                         |                                               | Internally powered                                                            |
| Rechargeable Battery                                                 |                                               |                                                                               |
| Model Name                                                           |                                               | MO3-REB                                                                       |
| Type                                                                 |                                               | Lithium-ion                                                                   |
| Output                                                               |                                               | 3.6 Vdc, 11.34 Wh                                                             |
| Dimension                                                            |                                               | 24.7 x 22 x 77.8 mm ± 5%                                                      |
| Weight                                                               |                                               | 66 g ± 5%                                                                     |
| Capacity                                                             |                                               | 3,150 mAh                                                                     |
| Battery Charger                                                      |                                               |                                                                               |
| Model Name                                                           |                                               | MO3-BC1                                                                       |
| Input                                                                |                                               | 12 Vdc, 3 A                                                                   |
| Dimension                                                            |                                               | 110 x 66 x 50.2 mm ± 5%                                                       |
| Weight                                                               |                                               | 158 g ± 5%                                                                    |
| Wireless Hub                                                         |                                               |                                                                               |
| Model Name                                                           |                                               | MO3-WFH                                                                       |
| Input                                                                |                                               | 5 Vdc / 3 A                                                                   |
| Dimension                                                            |                                               | 109.6 x 65.6 x 51.9 mm ± 5%                                                   |
| Weight                                                               |                                               | 168 g ± 5%                                                                    |
| Medit Power Adapter                                                  |                                               |                                                                               |
| Model Name                                                           |                                               | CDM045A3-W200V                                                                |
| Dimension                                                            |                                               | 31.5 x 53.5 x 58 mm ± 5%                                                      |
| Weight                                                               |                                               | 142 g ± 5%                                                                    |
| Input Voltage                                                        |                                               | Universal 100 – 240 Vac, 50 – 60 Hz input, without any slide switch           |
| Output                                                               |                                               | 5 Vdc / 3 A<br>9 Vdc / 3 A<br>12 Vdc / 3 A<br>15 Vdc / 3 A<br>20 Vdc / 2.25 A |
| EMI                                                                  |                                               | CE / FCC Class B, Conduction & Radiation met                                  |
| Protection                                                           | OVP (Over Voltage Protection)                 |                                                                               |
|                                                                      | SCP (Short Circuit Protection)                |                                                                               |
|                                                                      | OCP (Over Current Protection)                 |                                                                               |
| Protection Against Electric Shock                                    |                                               | Class II                                                                      |
| Mode of Operation                                                    |                                               | Continuous                                                                    |
| Reusable Tip                                                         |                                               |                                                                               |
| Dimension – Full Tip                                                 | Large 36.1 x 34.1 x 90.8 mm ± 5% (W x H x L)  |                                                                               |
|                                                                      | Medium 36.1 x 34.1 x 90.4 mm ± 5% (W x H x L) |                                                                               |
|                                                                      | Small 36.1 x 34.1 x 90.3 mm ± 5% (W x H x L)  |                                                                               |
| Dimension – Tip Head                                                 | Large 26.9 x 19.7 mm ± 5% (W x H)             |                                                                               |
|                                                                      | Medium 22.4 x 16.3 mm ± 5% (W x H)            |                                                                               |
|                                                                      | Small 18.36 x 13.1 mm ± 5% (W x H)            |                                                                               |

|                                                                   |                                                        |                                             |
|-------------------------------------------------------------------|--------------------------------------------------------|---------------------------------------------|
| Weight                                                            | Large 14 g ± 5%<br>Medium 12 g ± 5%<br>Small 11 g ± 5% |                                             |
| Calibration Tool                                                  |                                                        |                                             |
| Model Name                                                        | MOI-CT1                                                |                                             |
| Dimension                                                         | 160 x 48.5 mm ± 5% (H x Ø)                             |                                             |
| Weight                                                            | 205 g ± 5%                                             |                                             |
| Wireless Module                                                   |                                                        |                                             |
| 5 GHz                                                             | Frequency Bands                                        | 5210 – 5855 MHz                             |
|                                                                   | Modes                                                  | 802.11ax 80 MHz                             |
| 6 GHz (Indoor Only)                                               | Frequency Bands                                        | 6025 – 6985 MHz                             |
|                                                                   | Modes                                                  | 802.11ax 160 MHz                            |
| 2.4 GHz (Bluetooth / Bluetooth LE)                                | Frequency Bands                                        | 2402 – 2480 MHz                             |
|                                                                   | Channels                                               | 79/40                                       |
|                                                                   | Modulation                                             | GFSK                                        |
| 2.4 GHz (Bluetooth LE)<br>Wireless Hub only.<br>Antenna included. | Frequency Bands                                        | 2402 – 2480 MHz                             |
|                                                                   | Channels                                               | 40                                          |
|                                                                   | Channel Bandwidth                                      | 1 MHz<br>2 MHz                              |
|                                                                   | e.i.r.p.                                               | 9.3 dBm                                     |
|                                                                   | Modulation                                             | GFSK                                        |
| Operating, Storage, and Transport Conditions                      |                                                        |                                             |
| Operating Condition                                               | Temperature                                            | 18 – 28°C (64.4 – 82.4°F)                   |
|                                                                   | Humidity                                               | 20 – 75% relative humidity (non-condensing) |
|                                                                   | Air Pressure                                           | 800 – 1,100 hPa                             |
| Storage Condition                                                 | Temperature                                            | -10 – 50°C (14 – 122°F)                     |
|                                                                   | Humidity                                               | 20 – 80% relative humidity (non-condensing) |
|                                                                   | Air Pressure                                           | 800 – 1,100 hPa                             |
| Transport Condition                                               | Temperature                                            | -10 – 50°C (14 – 122°F)                     |
|                                                                   | Humidity                                               | 20 – 80% relative humidity (non-condensing) |
|                                                                   | Air Pressure                                           | 620 – 1,200 hPa                             |

 CAUTION: The frequency bands listed below are mandatory for the specified regions.

| Region     | Center Frequencies (MHz)                             | Frequency Bands (MHz)                                 |
|------------|------------------------------------------------------|-------------------------------------------------------|
| Korea      | 5210, 5775, 6025, 6185, 6345, 6505, 6665, 6825, 6985 | 5150–5250, 5725–5850, 5925–6425, 6425–7125            |
| USA/Canada | 5210, 5775, 6025, 6185, 6345, 6505, 6665, 6825, 6985 | 5150–5250, 5725–5850, 5925–6425, 6425–6875, 6875–7125 |
| Europe     | 5210, 5775, 6025, 6185, 6345                         | 5150–5250, 5725–5850, 5945–6425                       |
| Japan      | 5210, 5775, 6025, 6185, 6345                         | 5150–5250, 5725–5850, 5945–6425                       |

EC REP

**Meditrial Srl**  
Via Po 9 00198, Rome Italy

CH REP

**Meditrial Europe Ltd**  
Bahnhofstrasse 23 6300 Zug, Switzerland

UK REP

**Meditrial UK Ltd**  
19 The Circle, Queen Elizabeth Street, London, SE1 2JE, United Kingdom

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**Medit Corp.**