USER GUIDE

Medit Link (Revision 10) 2025.03.18

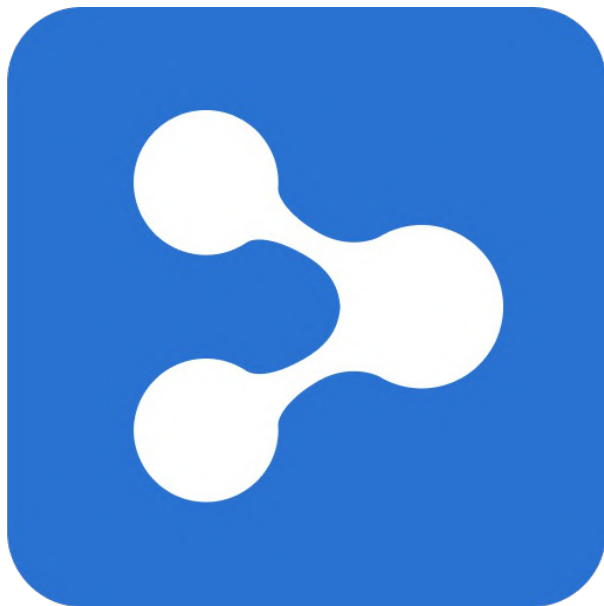
Table of Contents

Medit Link

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Administrator

Dentist, Manager

Have the authority to create a new dental clinic account and have full access to Medit Link.

- ✓ Scan/Order/Manage orders
- ✓ Manage members

Register

Lab

Administrator

Laboratory Director, Manager

Have the authority to create a new dental lab account and have full access to Medit Link.

- ✓ Scan/Manage ordered cases
- ✓ Manage members

Register

Member

Dental Hygienist, Dental Staff

Belong to the dental clinic their admin created and have limited access to Medit Link after the admin's approval.

- ✓ Scan/Order/Manage orders

Register

Member

Dental Technician, Laboratory Staff

Belong to the dental lab their admin created and have limited access to Medit Link after the admin's approval.

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Email Address

Check

Please ensure your email address is active. Medit Link will send a verification email.
Registering with a disposable and/or fake e-mail address can make you vulnerable to security issues and make it difficult to receive information from Medit Link.

Password

Password must be 8 to 16 characters long, including a combination of at least three letters, numbers, and special characters. (Passwords are accepted in English language only.)

Confirm Password

Please type in the same password as above.

Select Country/Region

☐ I agree to all of the Terms and Conditions.
 ☐ Personal Information Collection Agreement (Mandatory)

View

☐ Terms & Conditions (Mandatory)

View

☐ Agree to the Marketing and Promotion Agreement (optional)

View

Check this option to receive marketing emails, which include product updates, software news, and promotional events.

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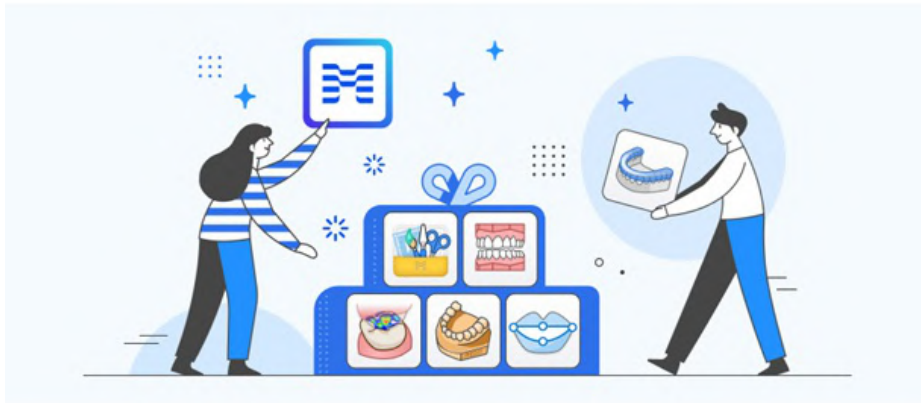
1 2 3

Utilize Various Apps

Download Medit Link to your desktop to utilize scanning programs and various apps.
All data generated is synced to the cloud and is also available on the web.

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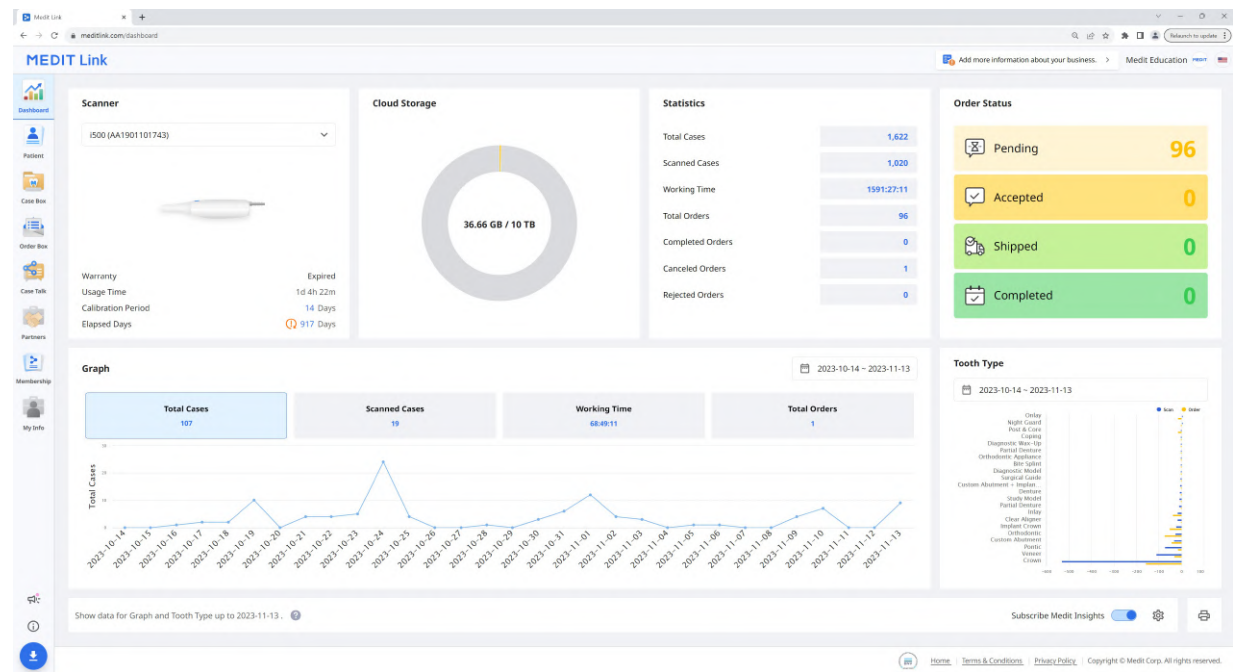
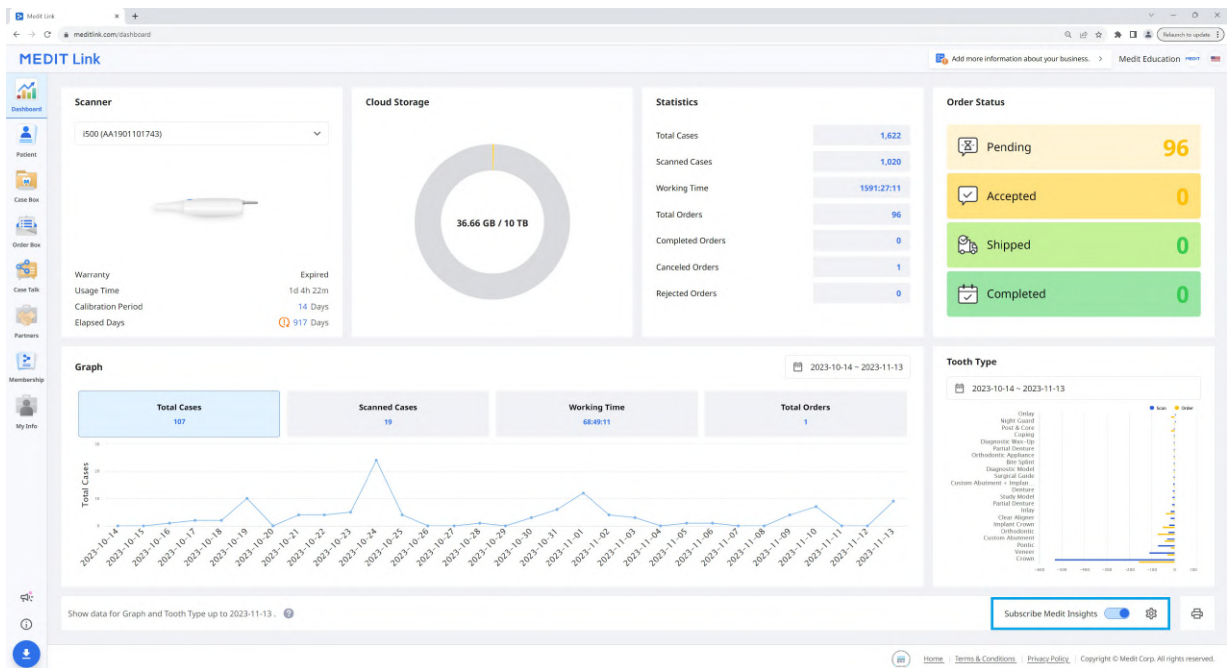
☒ Medit Link account


☐ Another email address

Enter another email address to receive the report email.

The report will be sent on the 2nd day of every month.

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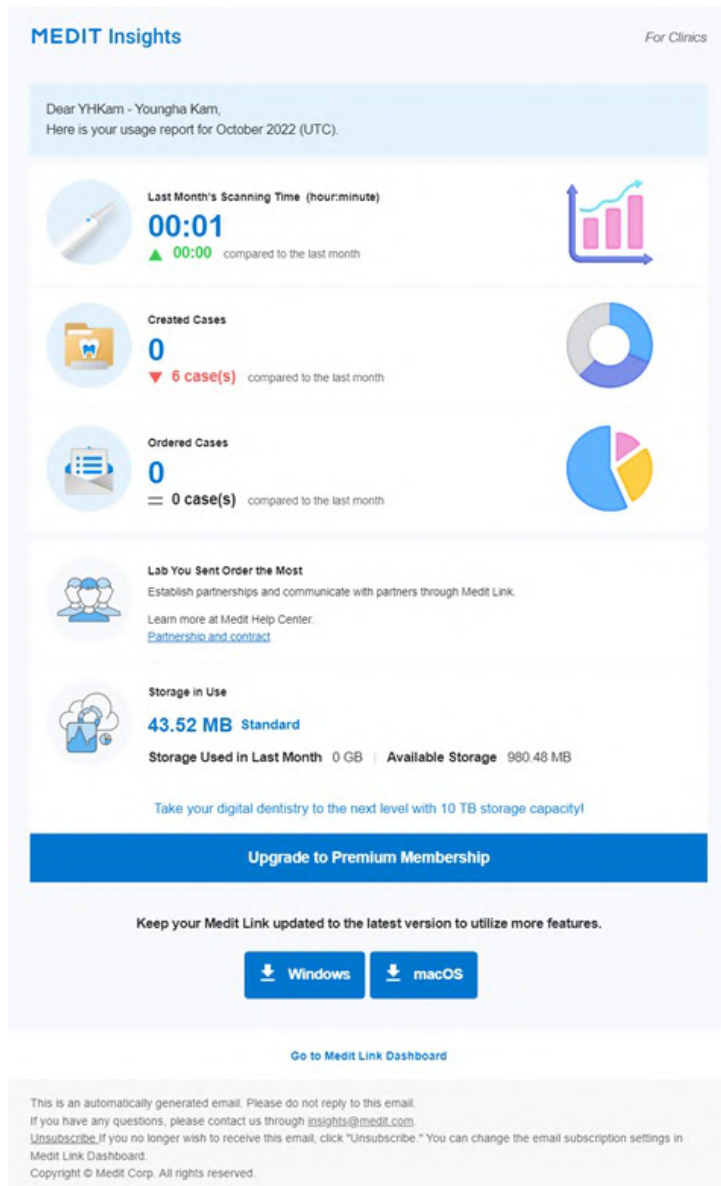
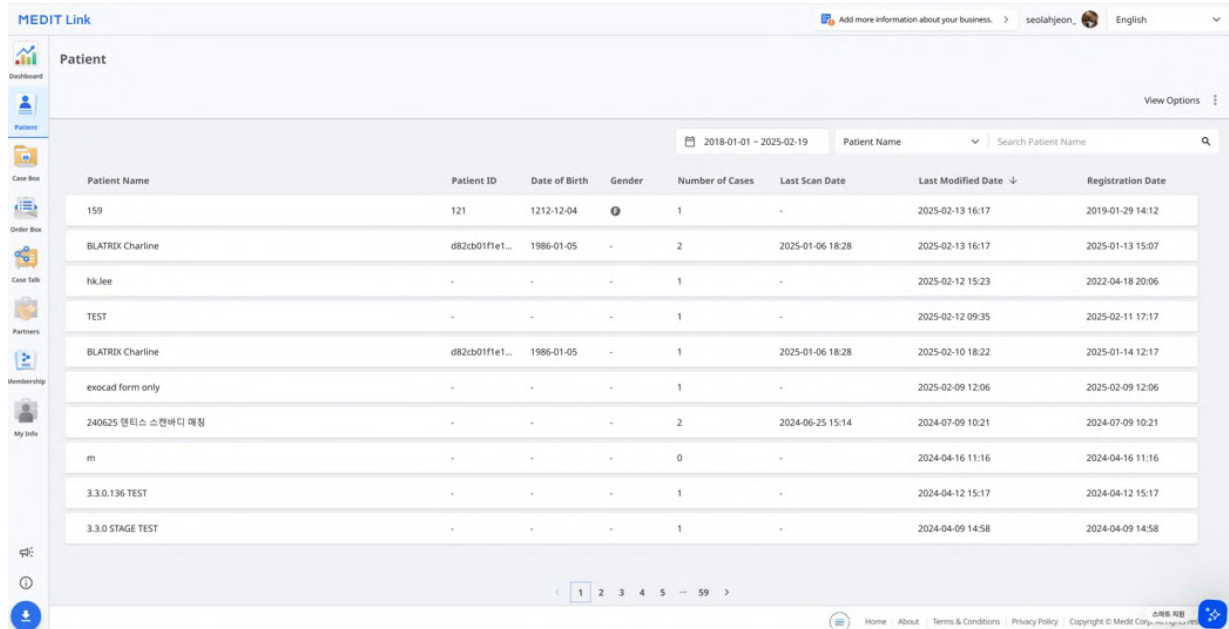
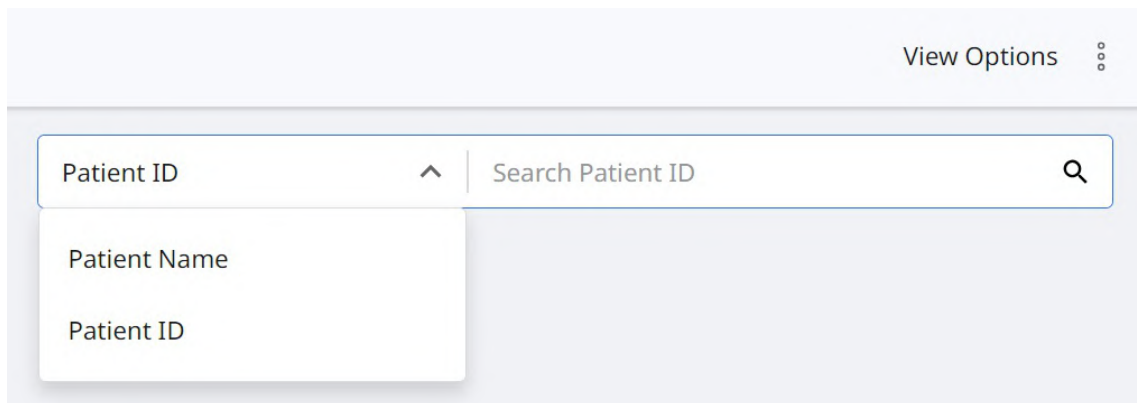
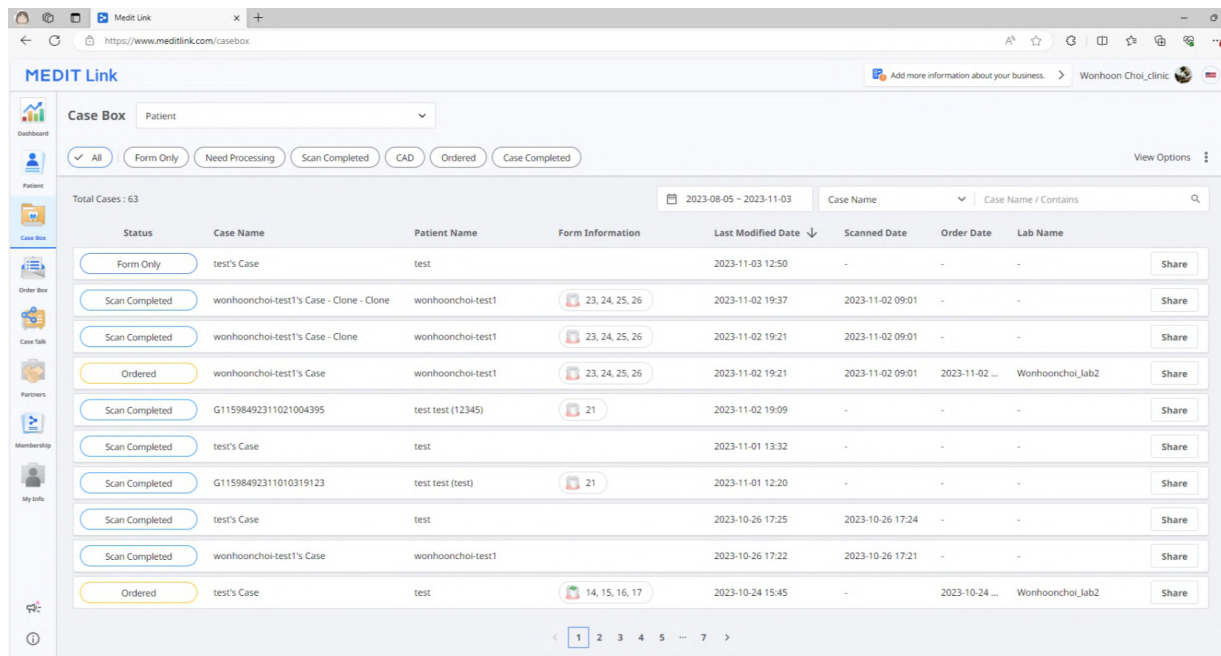
The screenshot shows the Medit Link Case Box interface. At the top, there's a 'Case Box' dropdown menu set to 'Patient'. Below it are filter buttons: 'All', 'Form Only', 'Need Processing', 'Scan Completed', 'CAD', 'Ordered', and 'Case Completed'. A 'Total Cases: 63' indicator is present. A date range filter shows '2023-08-05 - 2023-11-03'. A search bar for 'Case Name' is also visible. The main table lists cases with columns: Status, Case Name, Patient Name, Form Information, Last Modified Date, Scanned Date, Order Date, and Lab Name. Each row has a 'Share' button. The table is paginated with a '1' button highlighted.

Status	Case Name	Patient Name	Form Information	Last Modified Date	Scanned Date	Order Date	Lab Name
Form Only	test's Case	test		2023-11-03 12:50	-	-	-
Scan Completed	wonhoonchoi-test1's Case - Clone - Clone	wonhoonchoi-test1	23, 24, 25, 26	2023-11-02 19:37	2023-11-02 09:01	-	-
Scan Completed	wonhoonchoi-test1's Case - Clone	wonhoonchoi-test1	23, 24, 25, 26	2023-11-02 19:21	2023-11-02 09:01	-	-
Ordered	wonhoonchoi-test1's Case	wonhoonchoi-test1	23, 24, 25, 26	2023-11-02 19:21	2023-11-02 09:01	2023-11-02 ...	Wonhoonchoi_Lab2
Scan Completed	G11598492311021004395	test test (12345)	21	2023-11-02 19:09	-	-	-
Scan Completed	test's Case	test		2023-11-01 13:32	-	-	-
Scan Completed	G11598492311010319123	test test (test)	21	2023-11-01 12:20	-	-	-
Scan Completed	test's Case	test		2023-10-26 17:25	2023-10-26 17:24	-	-
Scan Completed	wonhoonchoi-test1's Case	wonhoonchoi-test1		2023-10-26 17:22	2023-10-26 17:21	-	-
Ordered	test's Case	test	14, 15, 16, 17	2023-10-24 15:45	-	2023-10-24 ...	Wonhoonchoi_Lab2

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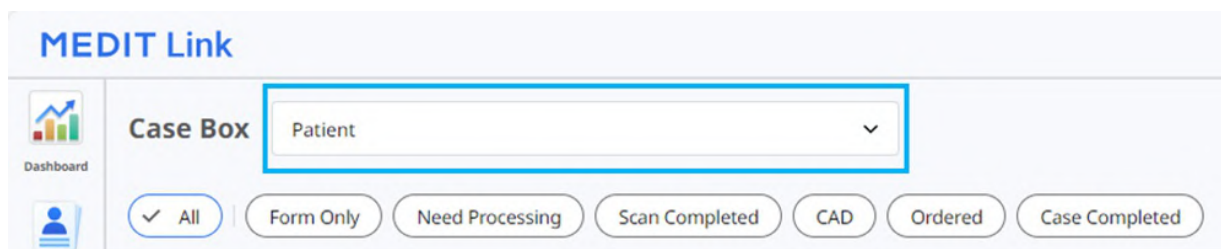The screenshot shows the Medit Link Case Box interface. The 'Case Box' dropdown menu is highlighted with a blue box, and it is set to 'Patient'. Below it are filter buttons: 'All', 'Form Only', 'Need Processing', 'Scan Completed', 'CAD', 'Ordered', and 'Case Completed'.

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+

☰

❗ Create a shared link and share it with specific people.

Expiration Date

📅 2022-06-28

❗ You can extend the period up to 90 days from today.

PIN Code

☐ PIN Code

Enter PIN Code

❗ Set the PIN code using six digits for privacy protection.

☒ Share the patient's name

☒ Grant Download Permission [Premium Feature](#)

❗ If you check this option, all users who have access to the shared link can download the files.

❗ The premium feature is provided to all users for a limited time.
Feature availability is subject to change depending on the membership plan.

☒ Allow comments

❗ If you check this option, you will be able to communicate through the shared link.

Comment

Please enter your comment.

0/1024

⚠

By pressing the "Get Link" button, you agree that you understand personal information can be collected and provided to a third party (or shared in any public domain) ONLY with the consent of the person whose images on which you have worked. You also agree the consent of that person was truthfully obtained, and Medit is not responsible for any problems which may arise there after.

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Invite Partner

Invite your desired partner to join Medit Link. When the lab accepts your invitation, a partnership is automatically formed.

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The screenshot displays the 'Medit Link' interface. The top navigation bar includes 'Add more information about your business', a user profile 'seolahjeon_', and a language dropdown set to 'English'. The main navigation menu on the left lists 'Dashboard', 'Patient', 'Case Box', 'Order Box', 'Case Talk', 'Partners', 'Membership', and 'My Info'. The 'Membership' section is active, showing tabs for 'Membership', 'Payment Details', 'Redeem', and 'My Subscription'. The 'Medit Cloud Storage Plan' is highlighted, with options for '2 GB Plan' and '10 TB Plan'. The '10 TB Plan' is selected, showing a 'Free (Promotion)' status with benefits: 'A paid plan', '10 TB of cloud storage', 'Efficient for users with many cases and various data', and 'One-month free trial included'. A 'Subscribe' button is visible. Below, '10 TB Plan Benefits' are listed, including a '10 TB Plan' with a code 'qqqq' and an expiration date of '2025-06-26'. A note states: 'The promotion has been applied. Your 10 TB Plan is now active and will remain valid until 2025-05-14. View More. Learn more about the Medit Cloud Storage Plan in our Q&A Here.' The 'Cloud Storage' section on the right shows a 'Set Period' from '2018-01-01' to '2025-02-17' with filters for 'All', '1 Week Ago', '1 Month Ago', '6 Months Ago', and '1 Year Ago'. A circular progress indicator shows '3.46 GB / 10 TB'. Below, a table lists storage items: 'Total Selected' (0.00 KB), 'Raw Data' (52.36 GB), and 'Result Files & Attachments' (3.46 GB). The footer contains links for 'Home', 'About', 'Terms & Conditions', 'Privacy Policy', and 'Copyright © Medit Co.', along with a Korean text '소속 의료기관'.

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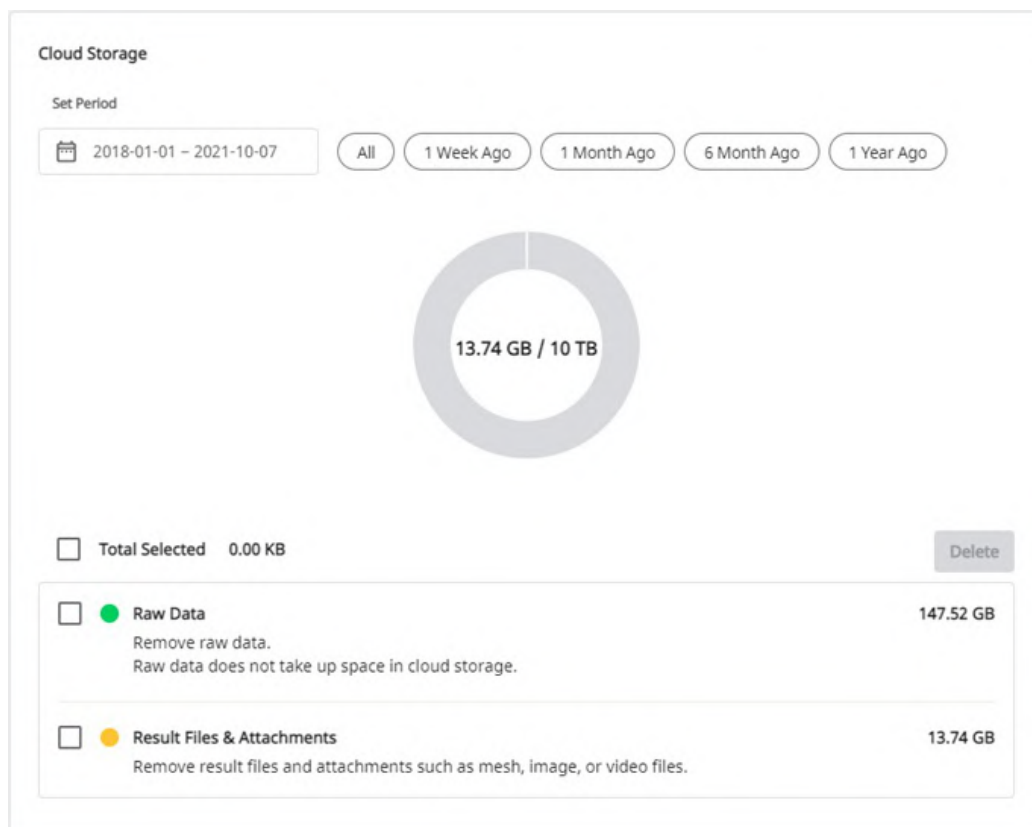
2 GB Plan **10 TB Plan**


Free (Promotion)

- ✓ A paid plan
- ✓ 10 TB of cloud storage
- ✓ Efficient for users with many cases and various data
- ✓ One-month free trial included

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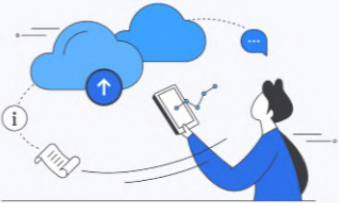10 TB Plan


9.99\$ USD / Monthly

- ✓ Subscribe now to keep your 10 TB Plan after the promotion expires.
- ✓ You will pay a subscription fee after the current promotion expires. If you don't wish to pay, you can unsubscribe before payment.

* VAT is not included.
* Learn more about the Medit Cloud Storage Plan in our Q&A [Here](#).

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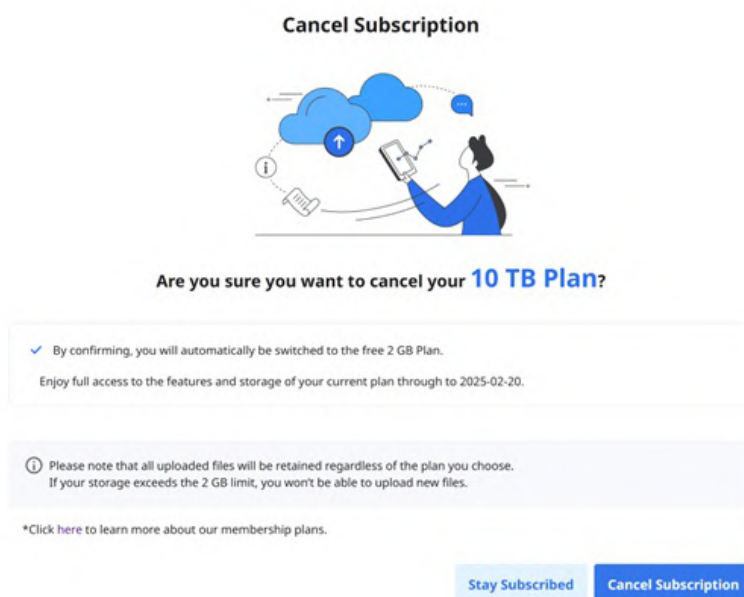Confirm Subscription Cancellation

Your subscription will be canceled on 2025-02-20.

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Confirm Subscription Cancellation

Your subscription has been canceled.


You can use your current plan until 2025-02-20. After that, the cloud capacity will be limited to 2 GB. If the used space exceeds the limit, data will not be backed up, and functions such as Order and Share will not be available.

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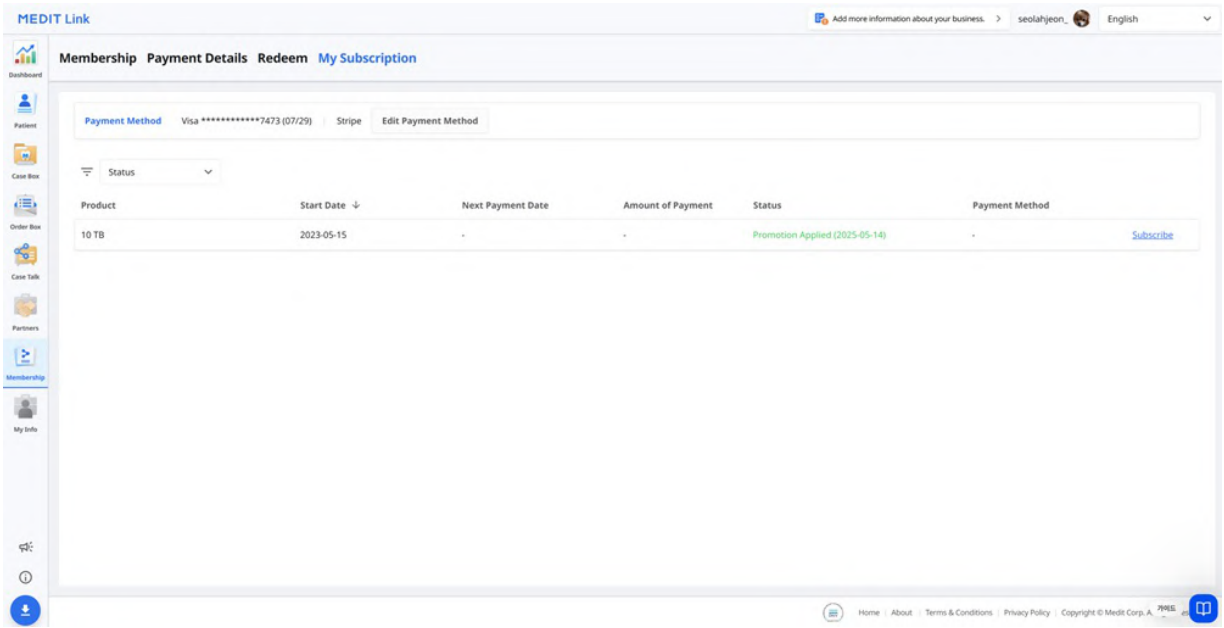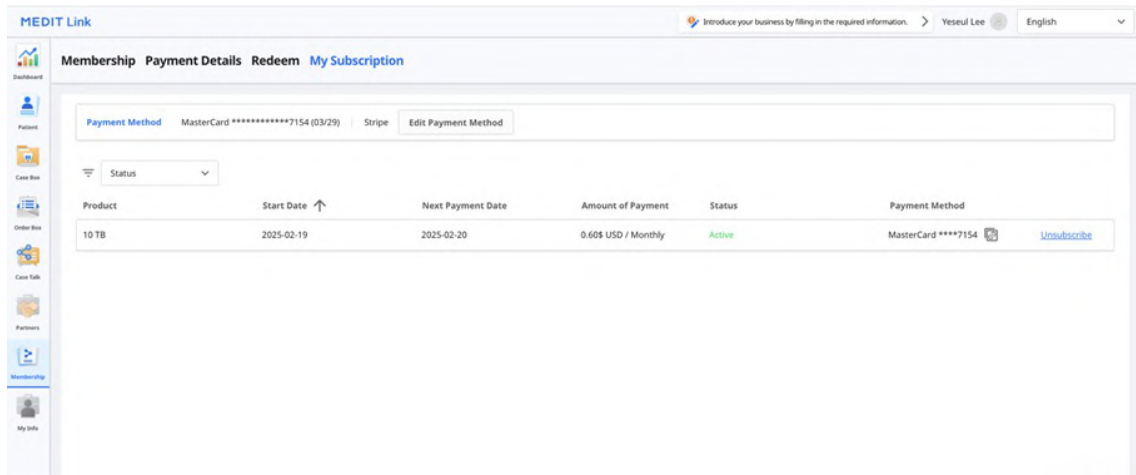
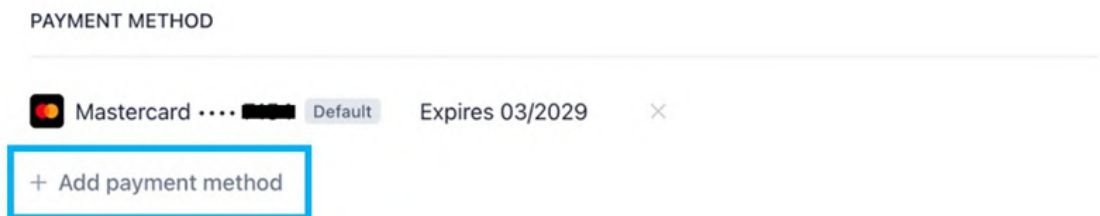
Country

South Korea

By providing your card information, you allow MEDIT USA Inc to charge your card for future payments in accordance with their terms.

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CURRENT SUBSCRIPTION

[DEV] Medit Cloud Storage 10TB Plan for Clinics
\$9.99 per day

[View details](#) ▾

Your subscription renews on February 20, 2025.

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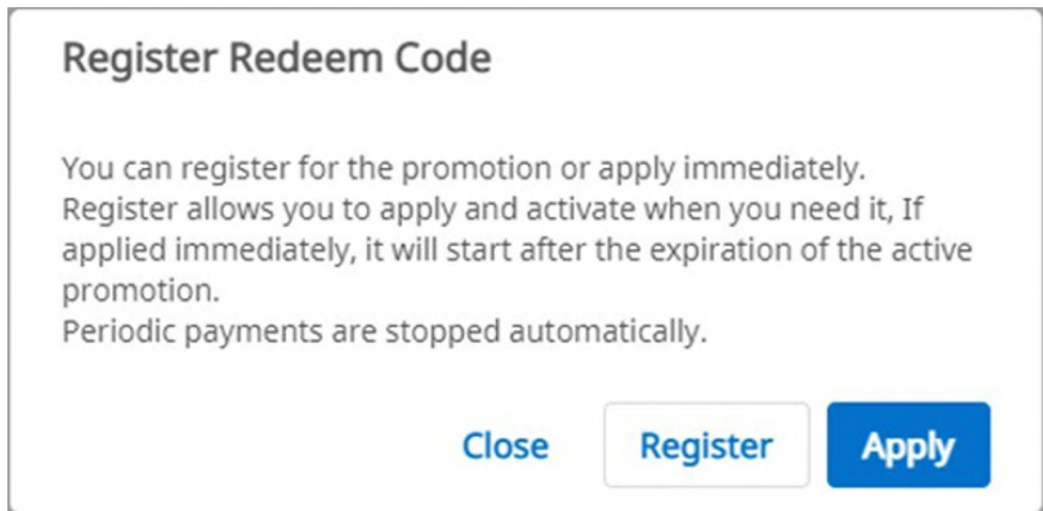
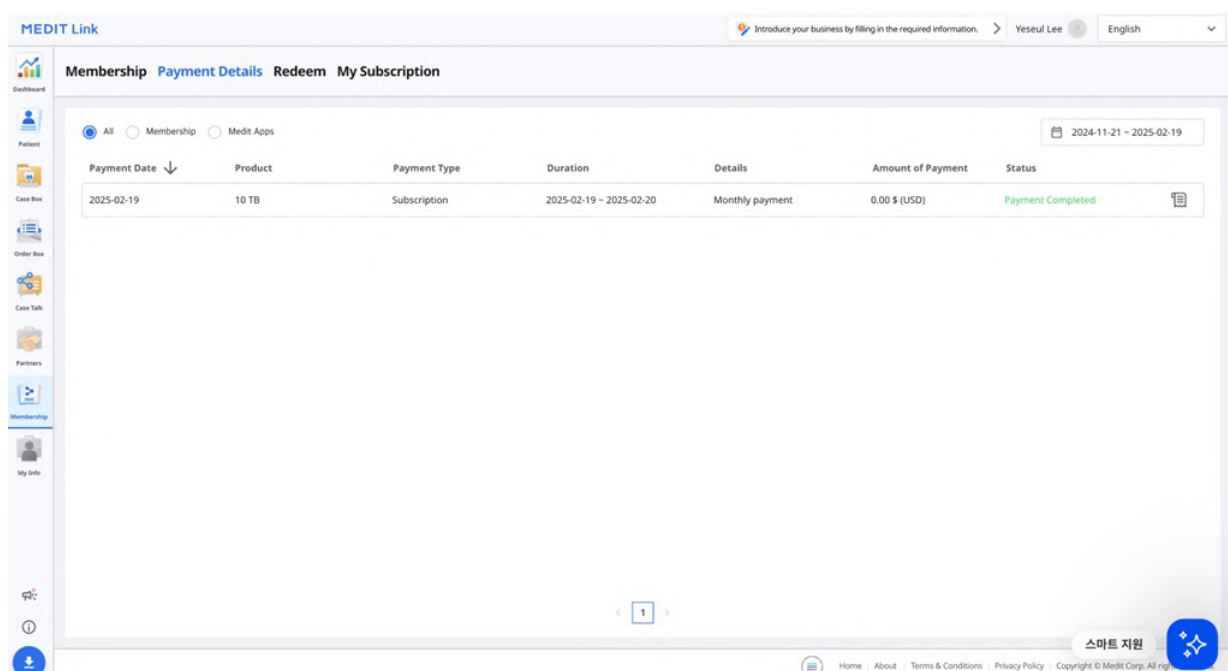
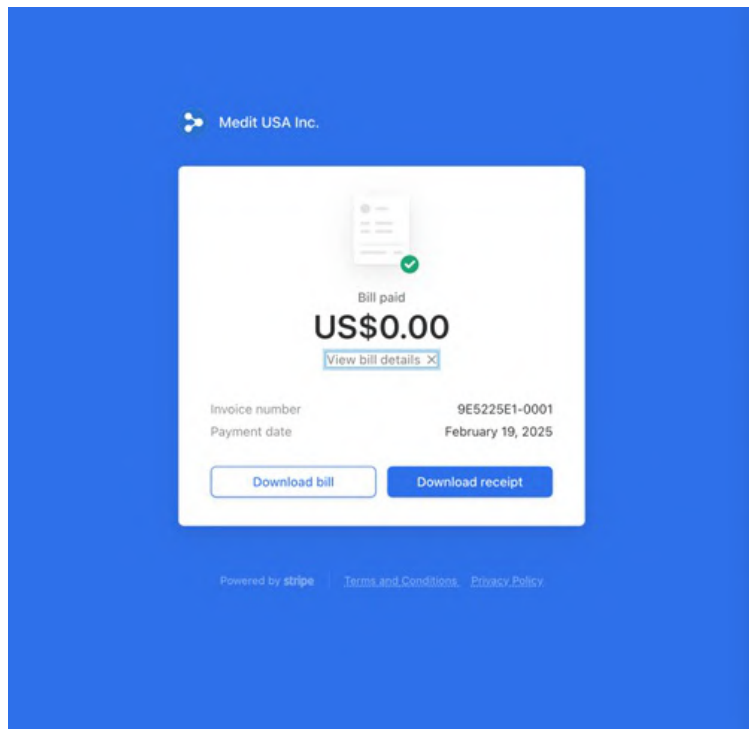
ITEM

FEBRUARY 19 - FEBRUARY 20, 2025

Medit Cloud Storage 10TB Plan for Clinics US\$0.60
Quantity 1 piece

subtotal US\$0.60

1 Month Free for New Subscribers(100% 할인) -US\$0.60

Payment amount US\$0.00

Amount paid US\$0.00

Remaining amount US\$0.00

Do you have any questions? [Contact Medit USA Inc.](#)

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Country/Region

Korea, Republic of

📘 Your chosen language is applied to the email sent by Medit Link and the notification.

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
Edit

Cover and Logo

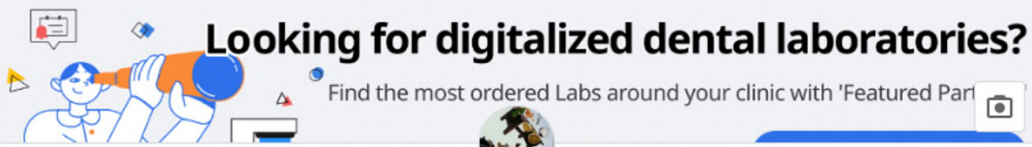Looking for digitalized dental laboratories?
Find the most ordered Labs around your clinic with 'Featured Partners'

1

2

3

Services

Select

▼

Specialties (choose up to 3)

Select

▼

 Select the services your business specializes in from the list.

About Us

e.g., Welcome to OOO Dental Clinic.
We are a licensed dental clinic situated in California specializing mainly in aesthetic restorations and implantology.

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Notice
NFT Service Termination Notice
2023-04-19

To our users,

We want to thank all Medit Link users who have supported and used our NFT service.

Unfortunately, we must inform you that this service will be terminated following the upcoming program update.

Please see the details below for more information:

- **Feature in Question**
 - Share via NFT (in Medit Link Case Talk)
- **Effective Date**
 - Starting from the update in the first half of 2023 (scheduled for April 27)

※ The update schedule is subject to change.

We'll continuously work to provide better service.

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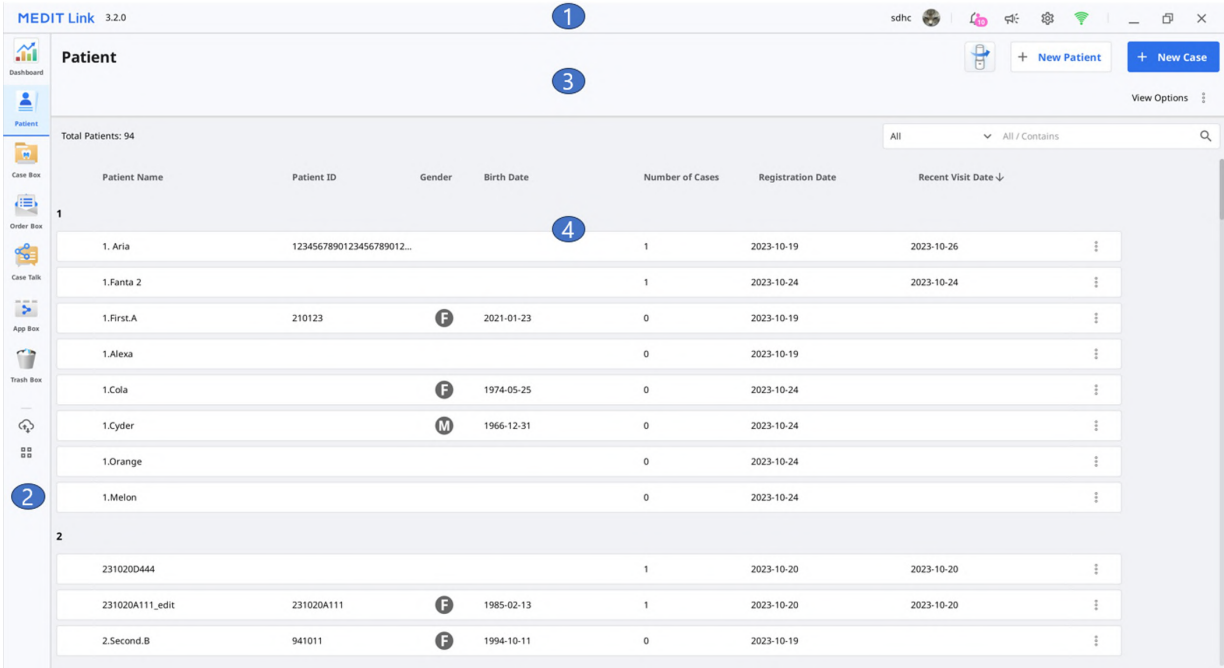
2

Delivery

Wonhoonchoi_Lab_Dev has dispatched the ordered package.

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Order

Wonhoonchoi_Lab_Dev has accepted the order.

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Order

Wonhoonchoi_Lab_Dev has accepted the order.

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Order

Wonhoonchoi_Lab_Dev has accepted the order.

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Order

Wonhoonchoi_Lab_Dev has accepted the order.

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Order

Wonhoonchoi_Lab_Dev has accepted the order.

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Work

Wonhoonchoi_Clinice_Dev_Member registered a new case test's Case.

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Work

Wonhoonchoi_Clinice_Dev_Member registered a new case test's Case.

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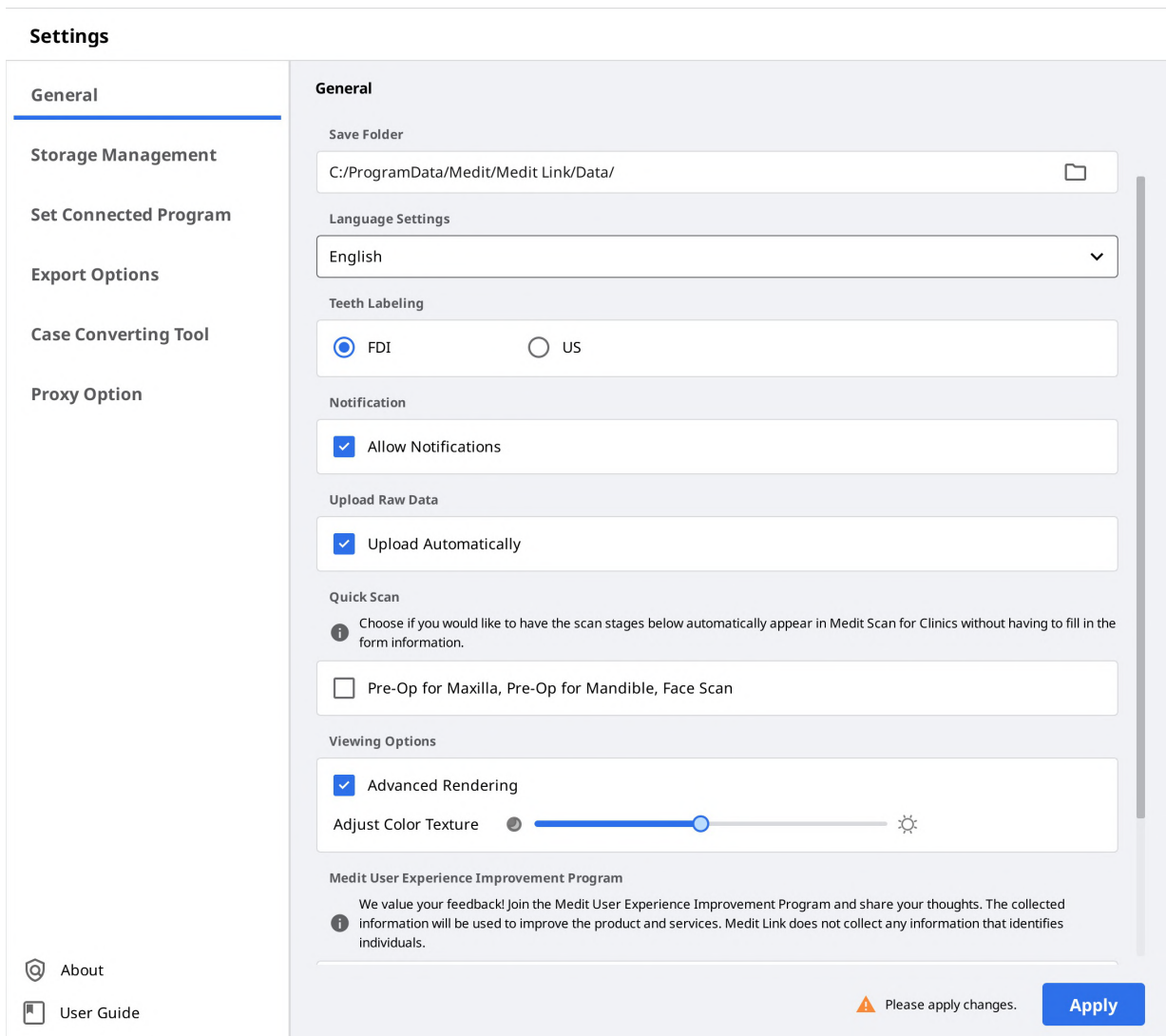
Network speed is good.

Network speed is fast enough.

Network speed is poor. Can be switched to offline mode automatically.

Either network is not connected or Medit Link is in offline mode.

Synchronization is in progress. Some case lists may not match.

Synchronization failed. The case lists may not match.

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Delete

Raw Data

Remove raw data created when using scan programs.

26.89 GB

Result Files & Attachments

Remove result files and attachments such as mesh, image, or video files.

2.54 GB

Others

It shows space used by other programs.

666.76 GB

Cloud Storage

Total Selected 0.00 KB

Delete

Raw Data

Remove raw data created when using scan programs. Raw data are not counted in cloud storage.

837.95 GB

Result Files & Attachments

Remove result files and attachments such as mesh, image, or video files.

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☒Compress Files

☒Combine Individual Mesh

Options when same folder or file already exists.

☒Always ask when exporting

☐Export a new file and overwrite the existing file. The file currently stored will be lost.

☐Create a folder and export the files

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Status	Case Name	Patient Name	Lab Name	Deletion Date	Deleted By	Size
Scan Completed	Jira_외출판트 46, 47's Case - Clone - Clone	외출판트 46, 47		2023-04-26 오전 9:12	wonhoon choi	66,654 KB
Form Only	G11598492304180810429	test test		2023-04-18 오후 5:15	wonhoon choi	0 KB
Form Only	G11598492304180806313	test test		2023-04-18 오후 5:15	wonhoon choi	0 KB
Form Only	G11598492304180811349	tes test		2023-04-18 오후 5:15	wonhoon choi	0 KB
Form Only	G11598492304180754083	Glidewell_fastdesign		2023-04-18 오후 5:00	wonhoon choi	0 KB
Form Only	G11598492304180733332	Glidewell_fastdesign		2023-04-18 오후 5:00	wonhoon choi	0 KB
Form Only	G11598492304180648146	test test		2023-04-18 오후 5:00	wonhoon choi	0 KB
Case Completed	Multi Occlusion Normal Case's Case - Clone	Multi Occlusion Normal Case		2023-04-18 오전 10:50	wonhoon choi	23,652 KB
CAD	Multi Occlusion Normal Case's Case - Clone - Clone	Multi Occlusion Normal Case		2023-04-18 오전 10:50	wonhoon choi	23,652 KB
Form Only	Occlusion Test's Case	Occlusion Test		2023-04-17 오후 6:31	wonhoon choi	0 KB
Form Only	Test Occlusion 's Case	Test Occlusion		2023-04-17 오후 6:31	wonhoon choi	0 KB
Form Only	Occlusion Test's Case	Occlusion Test		2023-04-17 오후 6:31	wonhoon choi	0 KB
Form Only	Occlusion Test's Case	Occlusion Test		2023-04-17 오후 6:26	wonhoon choi	0 KB

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No Date

Data TypeAll

☐	Case Name	Patient Name	Scanned Date↓	Size	Status
☐	BLATRIX Charline's Case	BLATRIX Charline	2025-01-06 PM 6:28	729.30 MB	Upload Required
☐	BLATRIX Charline's Case	BLATRIX Charline	2025-01-06 PM 6:28	729.30 MB	Upload Required
☐	BLATRIX Charline's Case - 복제	BLATRIX Charline	2025-01-06 PM 6:28	733.35 MB	Upload Required
☐	BLATRIX Charline's Case - Clone	BLATRIX Charline	2025-01-06 PM 6:28	734.01 MB	Upload Required
☐	test의 케이스	Patient TEST	2023-04-26 PM 7:58	1.25 GB	Upload Required
☐	i700 HD's Case	i700 HD	2023-04-24 PM 5:23	911.24 MB	Upload Required
☐	Mandibular Movement	MEDIT	2023-04-24 PM 1:59	1.88 GB	Upload Required
☐	Mandibular Movement Case	MEDIT	2023-04-24 PM 1:59	1.75 GB	Upload Required
☐	Mandibular Movement Case - Clo...	MEDIT	2023-04-24 PM 1:59	1.75 GB	Upload Required
☐	서라_Mandibular Movement_TEST	Mandibular Movement_export_test	2023-04-10 PM 3:39	34.73 MB	Upload Required
☐	서라_Mandibular Movement_TEST...	Mandibular Movement_export_test	2023-04-10 PM 3:39	731.43 MB	Upload Required

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
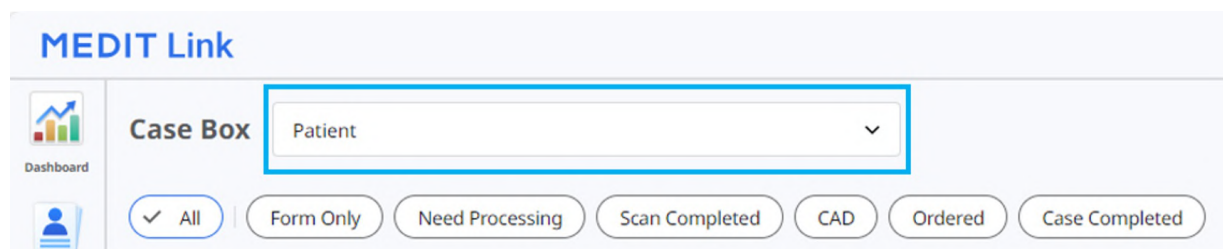
The screenshot shows the MEDIT Link interface. On the left is a sidebar with a 'Dashboard' icon and a user icon. The main area has a 'Case Box' dropdown menu currently set to 'Patient'. Below this are several filter buttons: 'All' (selected), 'Form Only', 'Need Processing', 'Scan Completed', 'CAD', 'Ordered', and 'Case Completed'.

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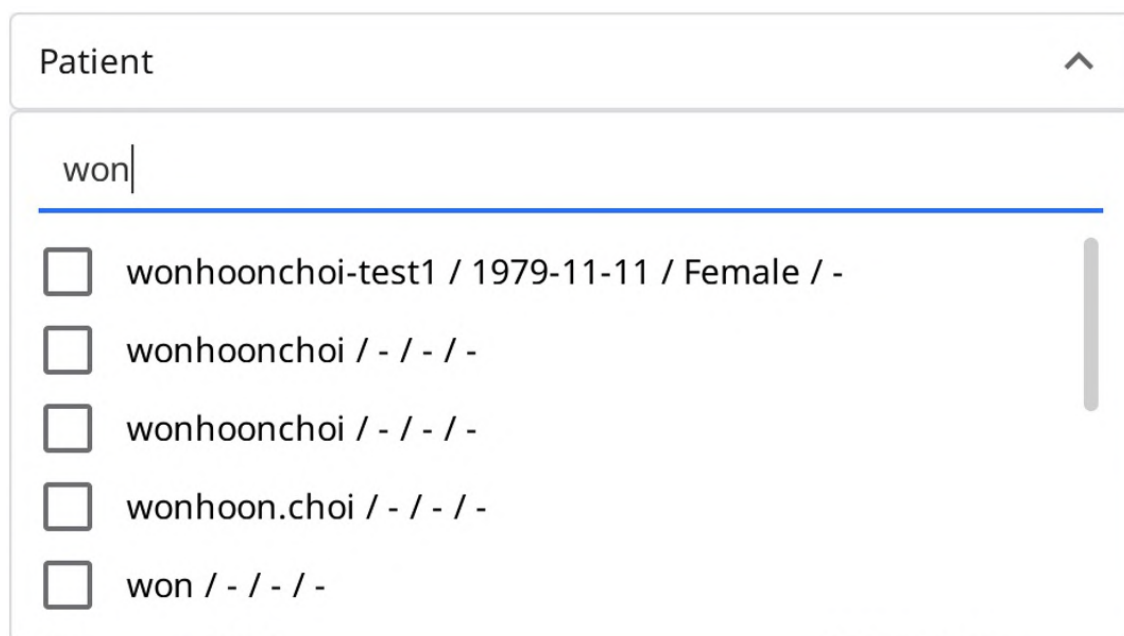This screenshot shows the expanded search dropdown for the 'Patient' case box. The search input field contains the text 'won'. Below the input field is a list of five search results, each preceded by a checkbox:

- ☐ wonhoonchoi-test1 / 1979-11-11 / Female / -
- ☐ wonhoonchoi / - / - / -
- ☐ wonhoonchoi / - / - / -
- ☐ wonhoon.choi / - / - / -
- ☐ won / - / - / -

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wonhoonchoi / - / - / -

+ 2

☒ All
 ☐ Form Only
 ☐ Need Processing
 ☐ Scan Completed
 ☐ CAD
 ☐ Ordered
 ☐ Case Completed

Total Cases: 20

2018-01-01~2023-10-24

All / Contains

Status	Patient Name	Case Name	Form Information	Last Modified Date	Scan Date	Lab Name	Order Date
Form Only	wonhoonchoi	originalCase	-	2023-10-05 오전 8:10			
Form Only	wonhoonchoi	wonhoonchoi's Case	-	2023-10-05 오전 8:09			
Scan Completed	wonhoonchoi	FileViwerTest_1	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-06-15 오전 8:14	2023-02-09 오전 8:35		
Scan Completed	wonhoonchoi	FileViwerTest_1 - Clone	2-Inlay / 3-Inlay / 4-Inlay / 5-Inlay / 6-Inlay / 8...	2023-05-25 오후 2:05	2023-02-09 오전 8:35		
Form Only	wonhoonchoi	FileViwerTest_1 - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-03-15 오후 3:06			
Scan Completed	wonhoonchoi	FileViwerTest_1 - Clone - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-03-09 오후 4:41	2023-02-09 오전 8:35		
Need Processing	wonhoonchoi	FileViwerTest_1 - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-02-14 오전 8:13	2023-02-09 오전 8:35		
Scan Completed	wonhoonchoi	FileViwerTest_1 - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-02-10 오전 11:27	2023-02-09 오전 8:35		
Form Only	wonhoonchoi	testsamplecase	-	2023-02-10 오전 7:39			
Scan Completed	wonhoonchoi	Multi Occlusion	-	2023-02-09 오전 9:26	2023-02-09 오전 9:01		
Scan Completed	wonhoonchoi	FileViwerTest_1 - Clone - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-02-09 오전 8:50	2023-02-09 오전 8:35		
Case Completed	wonhoonchoi	wonhoonchoi's Case	12-Inlay / 13-Inlay / 14-Inlay	2022-10-13 오후 10:30	2022-10-11 오후 5:16		
Case Completed	wonhoonchoi	wonhoonchoi's Case	-	2022-10-13 오후 10:28			

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
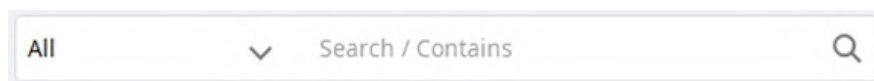
A search bar with a dropdown menu currently showing 'All'. To the right of the dropdown is a search icon (magnifying glass). The text 'Search / Contains' is visible inside the bar.

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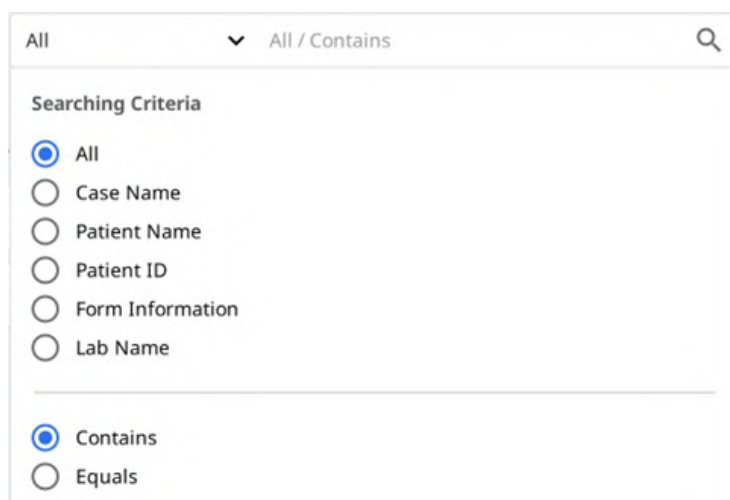
An expanded search criteria dropdown menu. It has a title 'Searching Criteria'. Under this title, there are two groups of radio buttons. The first group lists search criteria: 'All' (selected), 'Case Name', 'Patient Name', 'Patient ID', 'Form Information', and 'Lab Name'. The second group lists logical operators: 'Contains' (selected) and 'Equals'.

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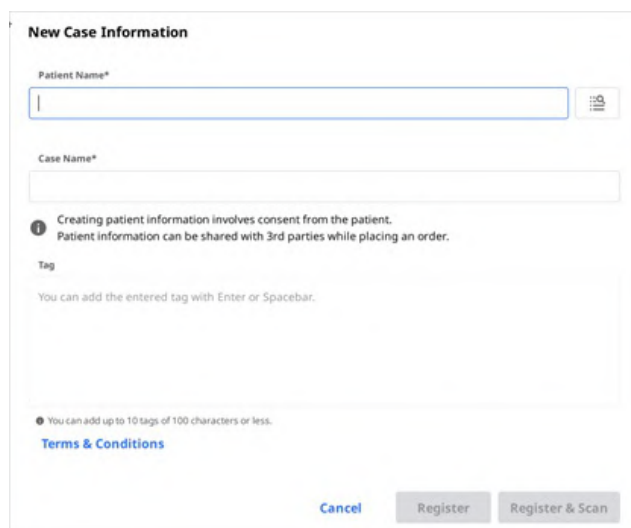The form titled "New Case Information" contains the following fields and elements:


- Patient Name***: A text input field with a search icon on the right.
- Case Name***: A text input field.
- Tag**: A text input field with a note: "You can add the entered tag with Enter or Spacebar."
- Terms & Conditions**: A link below the tag field.
- Buttons**: "Cancel", "Register", and "Register & Scan" at the bottom.

Below the Case Name field, there is a note: "Creating patient information involves consent from the patient. Patient information can be shared with 3rd parties while placing an order."

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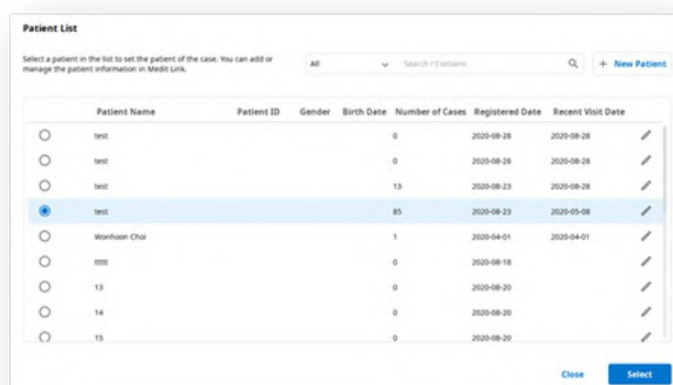
The "Patient List" dialog box shows a table of existing patients. The table has columns: Patient Name, Patient ID, Gender, Birth Date, Number of Cases, Registered Date, and Recent Visit Date. A search bar and a "+ New Patient" button are at the top right. The "Select" button is at the bottom right.

	Patient Name	Patient ID	Gender	Birth Date	Number of Cases	Registered Date	Recent Visit Date	
☐	test				0	2020-09-28	2020-09-28	
☐	test				0	2020-09-28	2020-09-28	
☐	test				13	2020-09-23	2020-09-28	
☒	test				85	2020-09-23	2020-05-08	
☐	Wanhuan Choi				1	2020-04-01	2020-04-01	
☐	test				0	2020-09-18		
☐	13				0	2020-09-20		
☐	14				0	2020-09-20		
☐	15				0	2020-09-20		

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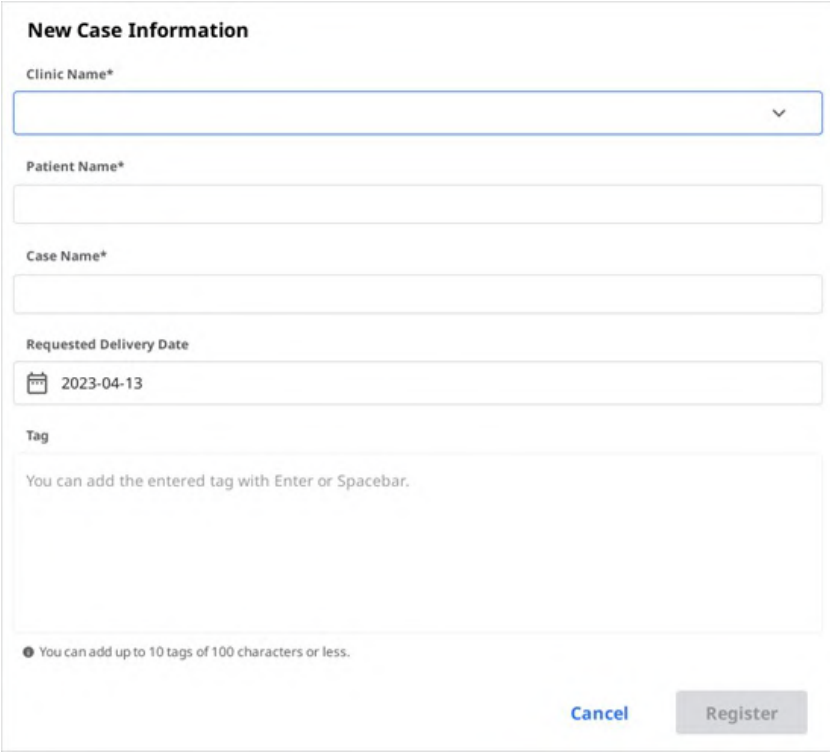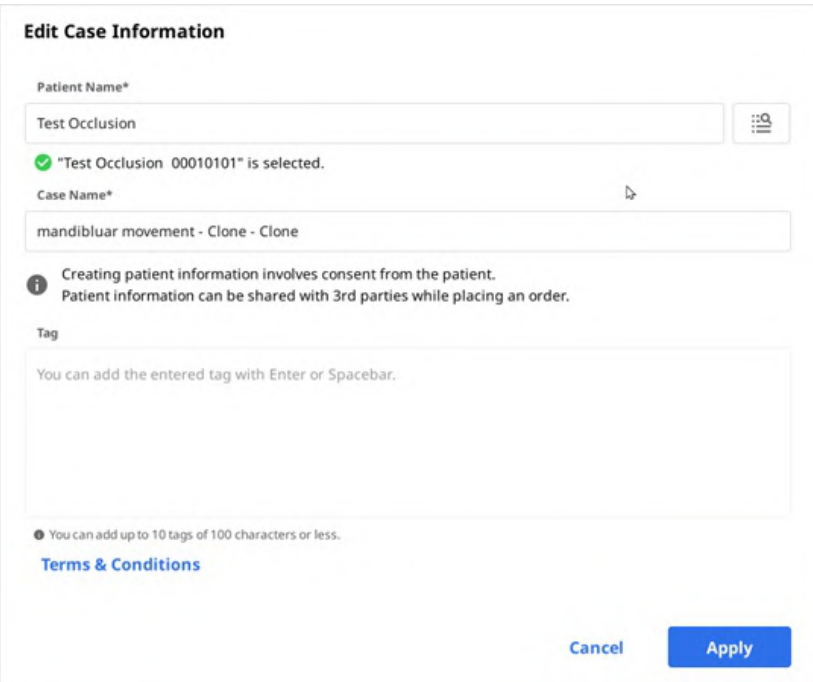
Case Name*

mandibluar movement - Clone - Clone

Creating patient information involves consent from the patient.
Patient information can be shared with 3rd parties while placing an order.

Tag

You can add the entered tag with Enter or Spacebar.

You can add up to 10 tags of 100 characters or less.

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Patient List

Select a patient in the list to set the patient of the case. You can add or manage the patient information in Medit Link.

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MEDIT Link 3.2.0

Case Box wonhoonchoi-test1 / 1979-11-11 / Female / - + 2 ▾

☒ All |
 ☐ Form Only |
 ☐ Need Processing |
 ☐ Scan Completed |
 ☐ CAD |
 ☐ Ordered |
 ☐ Case Completed

☒ **Export** ▾ |
 Status ▾ |
 |
 Selected Cases: 4

	Status	Patient Name	Case Name	Form Information
☒	Scan Completed	wonhoonchoi-test1	wonhoonchoi-test1's Case	-
☒	Form Only	wonhoonchoi-test1	wonhoonchoi-test1's Case	-
☒	Form Only	wonhoonchoi	originalcase	-
☒	Scan Completed	wonhoonchoi-test1	wonhoonchoi-test1's Case	-

Dashboard
 Patient
Case Box
 Order Box
 Case Talk
 App Box

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Selected Cases: 4

Case Completed

	Status		Case Name
☐	Ordered	test	test's Case
☒	Scan Completed	wonhoonchoi-test1	wonhoonchoi-test1's Case
☒	Scan Completed	test test(1234)	G11598492310180503424
☒	Scan Completed	test	test's Case
☒	Scan Completed	test	test's Case

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Patient Name*

Case Name*

 Creating patient information involves consent from the patient.
Patient information can be shared with 3rd parties while placing an order.

Tag

Doctorname × Specialmaterial × SpecialComment ×

 You can add up to 10 tags of 100 characters or less.

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Patient Name*

Case Name*

i Creating patient information involves consent from the patient.
Patient information can be shared with 3rd parties while placing an order.

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Tags

Doctorname x Specialmaterial x SpecialComment x

i You can add up to 10 tags of 100 characters or less.

☒ Include all occlusion relationships in the order

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