USER GUIDE

Medit Link (Revision 10) 2025.03.18

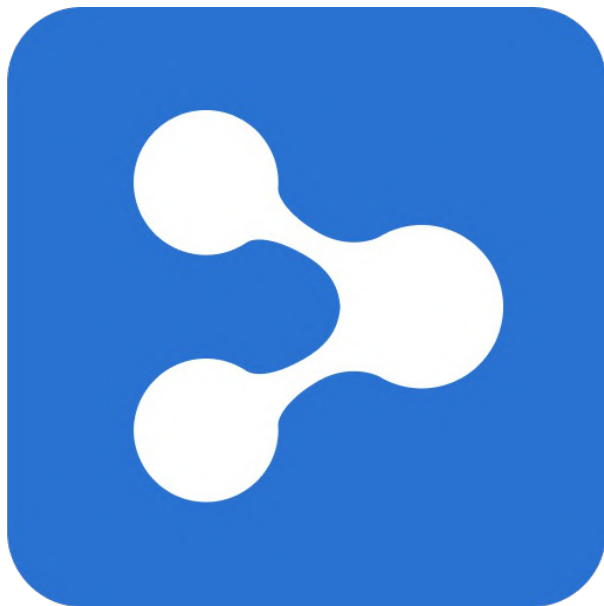
Table of Contents

Medit Link

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
Administrator

Dentist, Manager

Have the authority to create a new dental clinic account and have full access to Medit Link.

- ✓ Scan/Order/Manage orders
- ✓ Manage members

Register

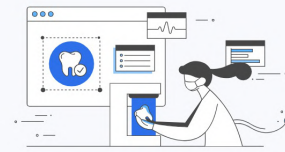Lab

Administrator

Laboratory Director, Manager

Have the authority to create a new dental lab account and have full access to Medit Link.

- ✓ Scan/Manage ordered cases
- ✓ Manage members

Register

Member

Dental Hygienist, Dental Staff

Belong to the dental clinic their admin created and have limited access to Medit Link after the admin's approval.

- ✓ Scan/Order/Manage orders

Register

Member

Dental Technician, Laboratory Staff

Belong to the dental lab their admin created and have limited access to Medit Link after the admin's approval.

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Email Address

Check

Please ensure your email address is active. Medit Link will send a verification email.

Registering with a disposable and/or fake e-mail address can make you vulnerable to security issues and make it difficult to receive information from Medit Link.

Password

Password must be 8 to 16 characters long, including a combination of at least three letters, numbers, and special characters. (Passwords are accepted in English language only.)

Confirm Password

Please type in the same password as above.

Select Country/Region

☐

I agree to all of the Terms and Conditions.

☐

Personal Information Collection Agreement (Mandatory)

☐

Terms & Conditions (Mandatory)

☐

Agree to the Marketing and Promotion Agreement (optional)

Check this option to receive marketing emails, which include product updates, software news, and promotional events.

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
Welcome, stage-test-wh!

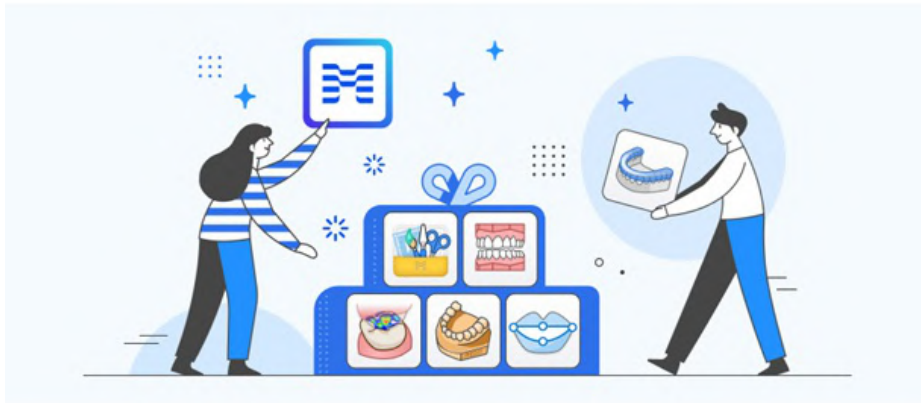Utilize Various Apps

Download Medit Link to your desktop to utilize scanning programs and various apps.
All data generated is synced to the cloud and is also available on the web.

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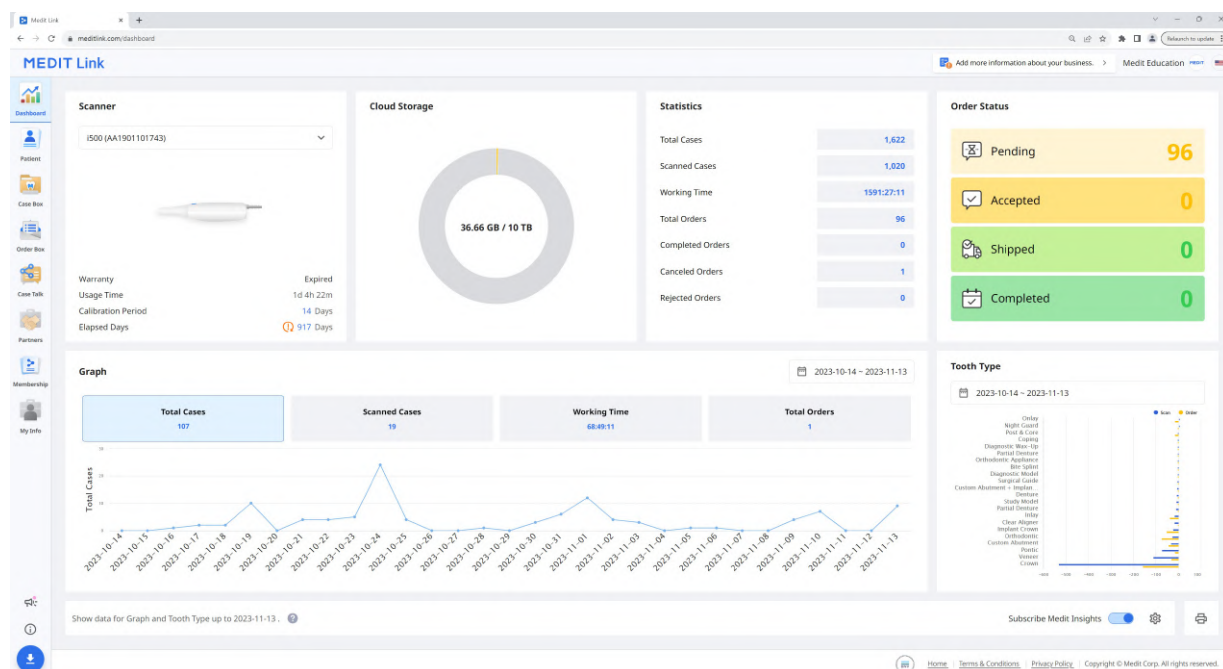
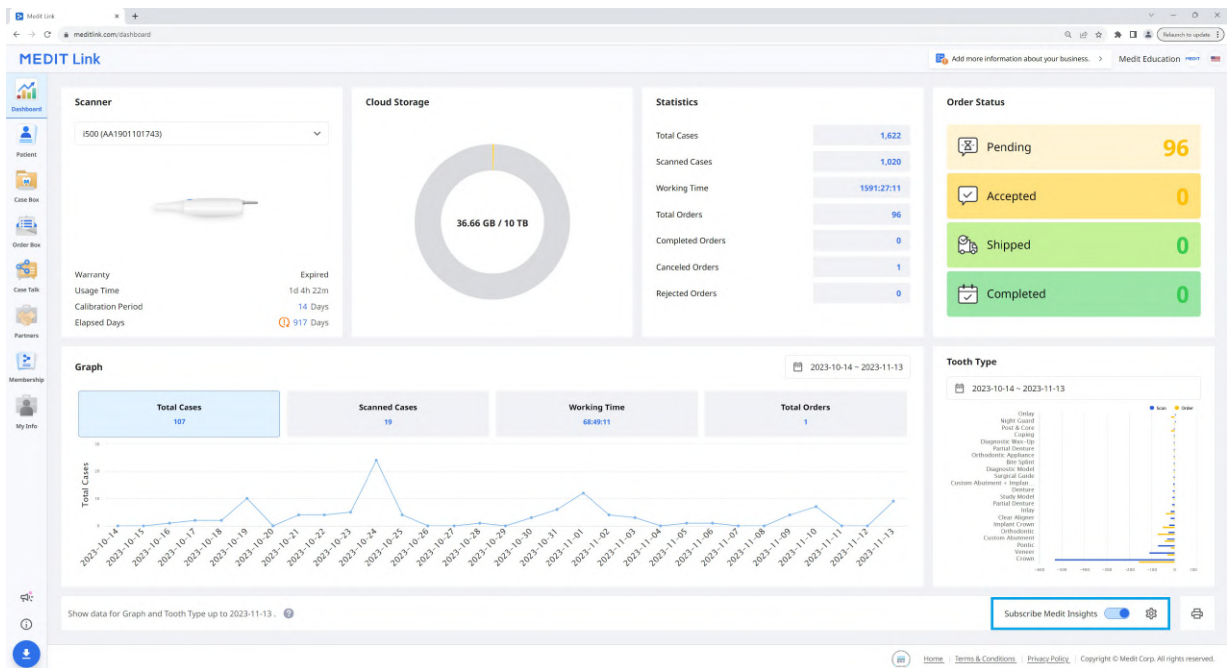
☒ Medit Link account

☐ Another email address


Enter another email address to receive the report email.

The report will be sent on the 2nd day of every month.

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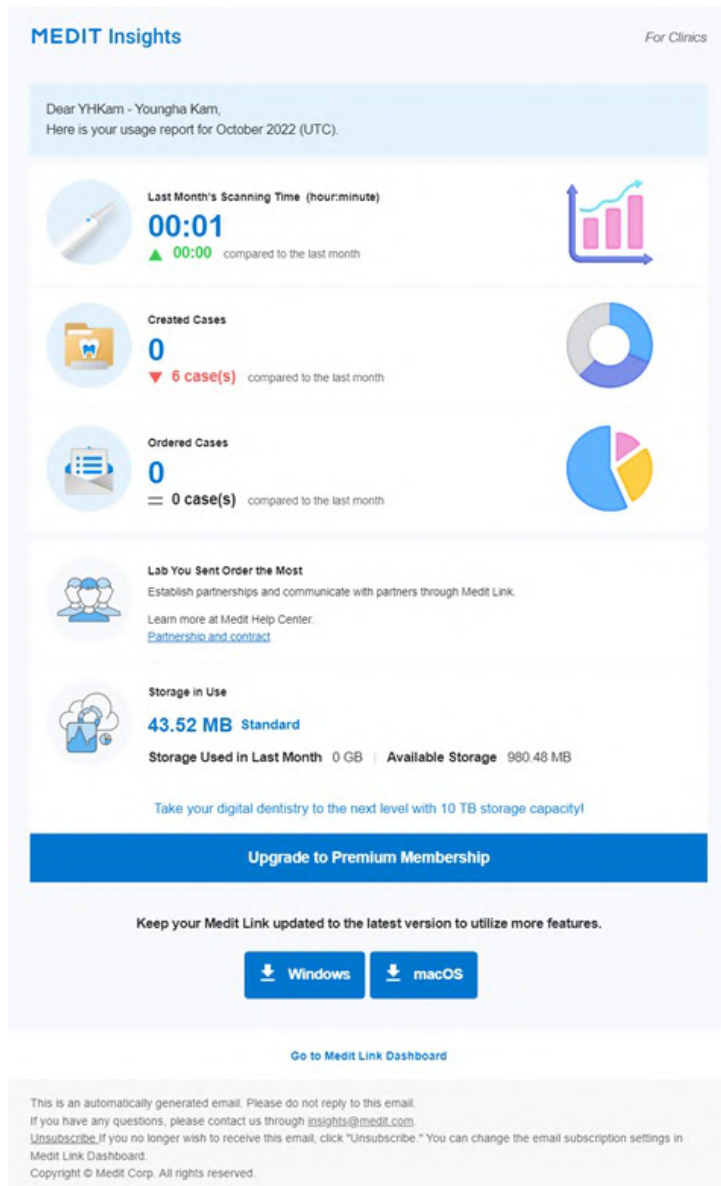
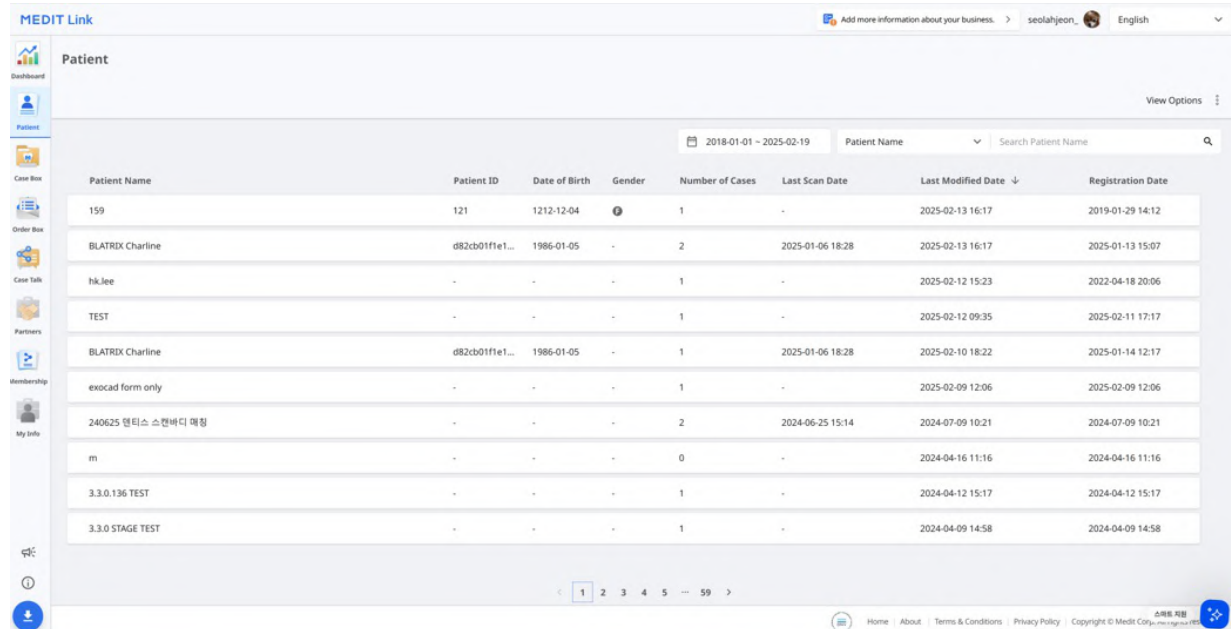
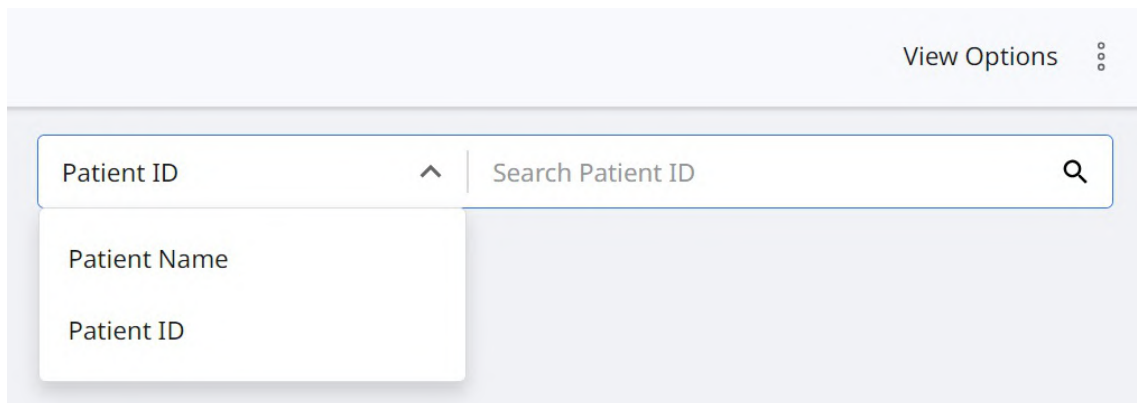
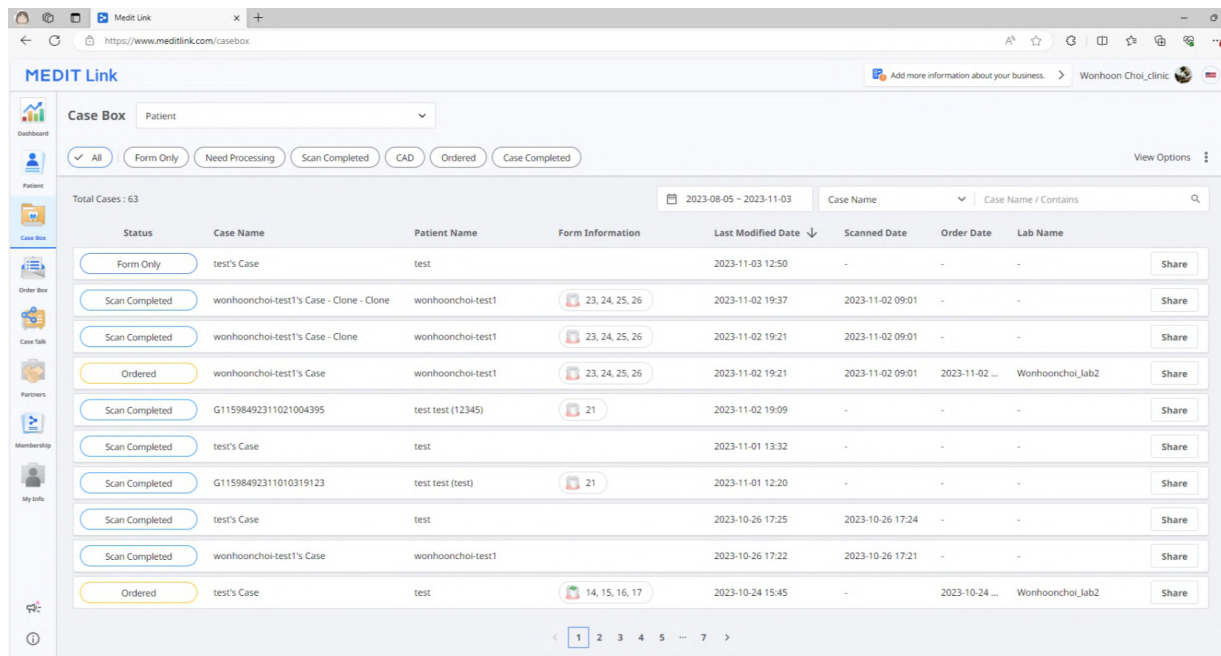
The screenshot shows the Medit Link Case Box interface. At the top, there's a header with the Medit Link logo and a user profile. Below the header, there's a sidebar with navigation icons for Dashboard, Patient, Case Box, Order Box, Case Talk, Partners, Membership, and My Info. The main content area is titled 'Case Box' and features a dropdown menu set to 'Patient'. Below this, there are filter buttons: 'All', 'Form Only', 'Need Processing', 'Scan Completed', 'CAD', 'Ordered', and 'Case Completed'. A 'View Options' link is also present. The table below shows a list of cases with columns for Status, Case Name, Patient Name, Form Information, Last Modified Date, Scanned Date, Order Date, and Lab Name. The table is paginated, showing cases 1 through 7.

Status	Case Name	Patient Name	Form Information	Last Modified Date	Scanned Date	Order Date	Lab Name
Form Only	test's Case	test		2023-11-03 12:50	-	-	-
Scan Completed	wonhoonchoi-test1's Case - Clone - Clone	wonhoonchoi-test1	23, 24, 25, 26	2023-11-02 19:37	2023-11-02 09:01	-	-
Scan Completed	wonhoonchoi-test1's Case - Clone	wonhoonchoi-test1	23, 24, 25, 26	2023-11-02 19:21	2023-11-02 09:01	-	-
Ordered	wonhoonchoi-test1's Case	wonhoonchoi-test1	23, 24, 25, 26	2023-11-02 19:21	2023-11-02 09:01	2023-11-02 ...	Wonhoonchoi_Lab2
Scan Completed	G11598492311021004395	test test (12345)	21	2023-11-02 19:09	-	-	-
Scan Completed	test's Case	test		2023-11-01 13:32	-	-	-
Scan Completed	G11598492311010319123	test test (test)	21	2023-11-01 12:20	-	-	-
Scan Completed	test's Case	test		2023-10-26 17:25	2023-10-26 17:24	-	-
Scan Completed	wonhoonchoi-test1's Case	wonhoonchoi-test1		2023-10-26 17:22	2023-10-26 17:21	-	-
Ordered	test's Case	test	14, 15, 16, 17	2023-10-24 15:45	-	2023-10-24 ...	Wonhoonchoi_Lab2

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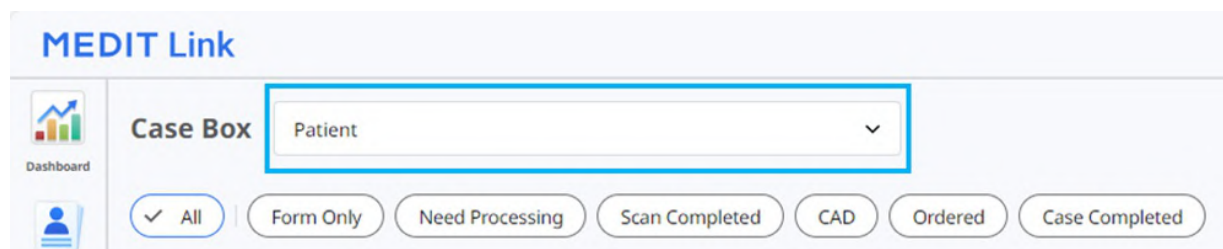The screenshot shows the Medit Link Case Box interface. The 'Case Box' dropdown menu is highlighted with a red box, and the 'Patient' option is selected. Below the dropdown, there are filter buttons: 'All', 'Form Only', 'Need Processing', 'Scan Completed', 'CAD', 'Ordered', and 'Case Completed'.

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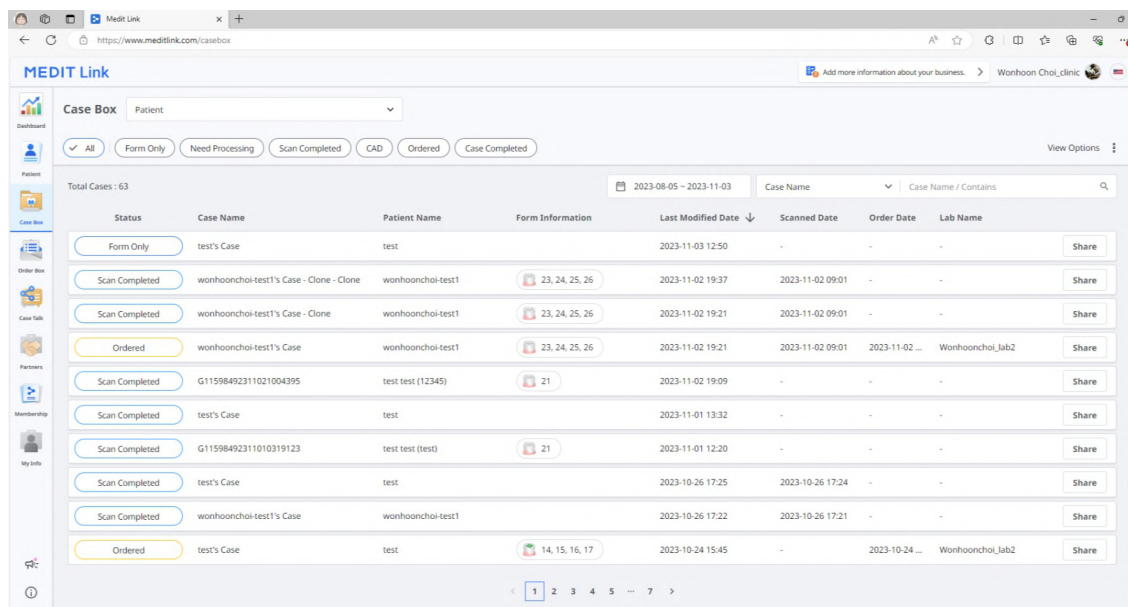
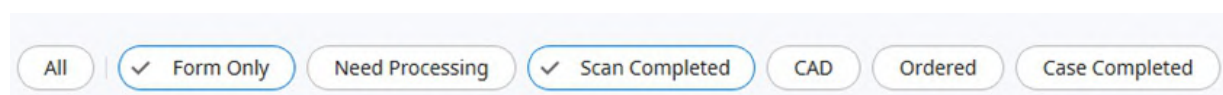
Case Box(診所)	All / Form Only / Need Processing / Scan Completed / CAD / Ordered / Case Completed
Work Box(技工所)	All / Form Only / Need Processing / Scan Completed / CAD / Done

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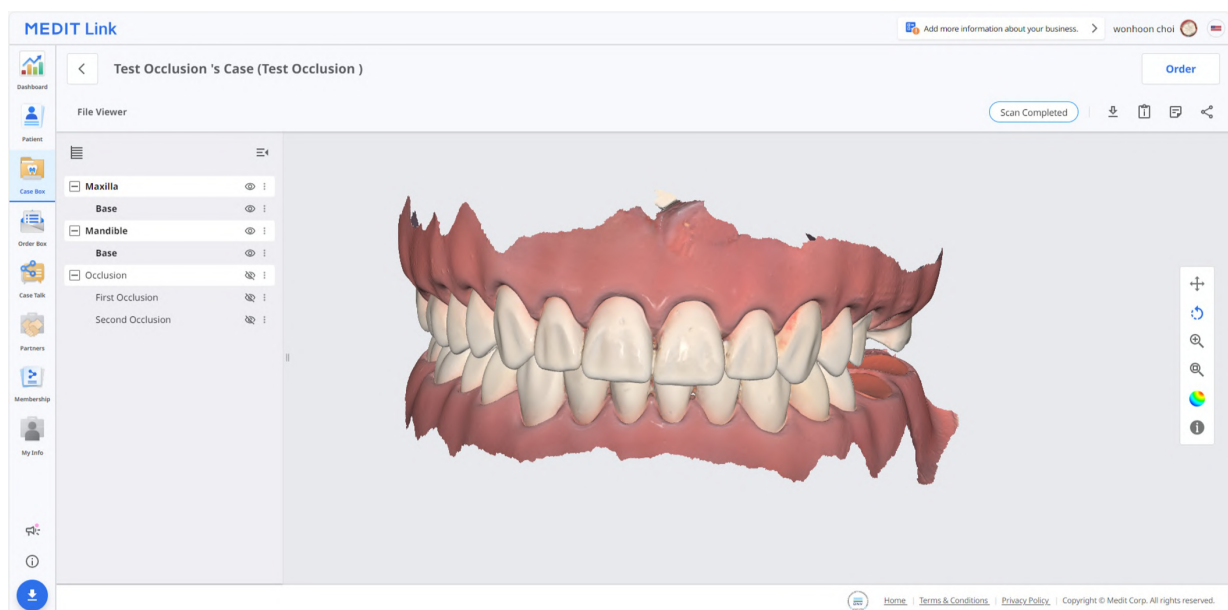
The screenshot displays the 'MEDIT Link' web application interface. The main content area is titled 'Occlusion Test's Case (Occlusion Test)'. It features a sidebar on the left with navigation icons for Dashboard, Patient, Case Box, Order Box, Case Talk, Partners, and Membership. The main area is divided into two sections: 'Order' and 'Order List'. The 'Order' section contains fields for Order ID (5996137), Patient Name (Occlusion Test), Partner Name (WonhoonChoi_Lab_Stage), Scanned Date (2023-04-20 PM 02:08), Requested Delivery Date (2023-04-27 AM 08:15), and Expected Delivery Date (2023-04-27 AM 08:15). The 'Order List' section displays a table with columns: No., Type, Method, Material, and Shade. The table contains four rows of data, all with 'Inlay' as the Type and 'Multi-Layered Zi...' as the Material. The interface also includes a 'Cancel Order' button and a 'Pending' status indicator.

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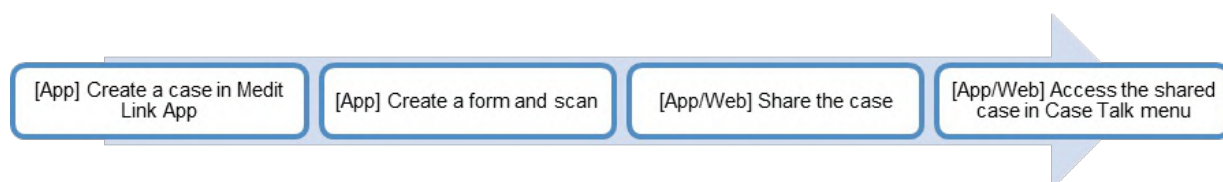
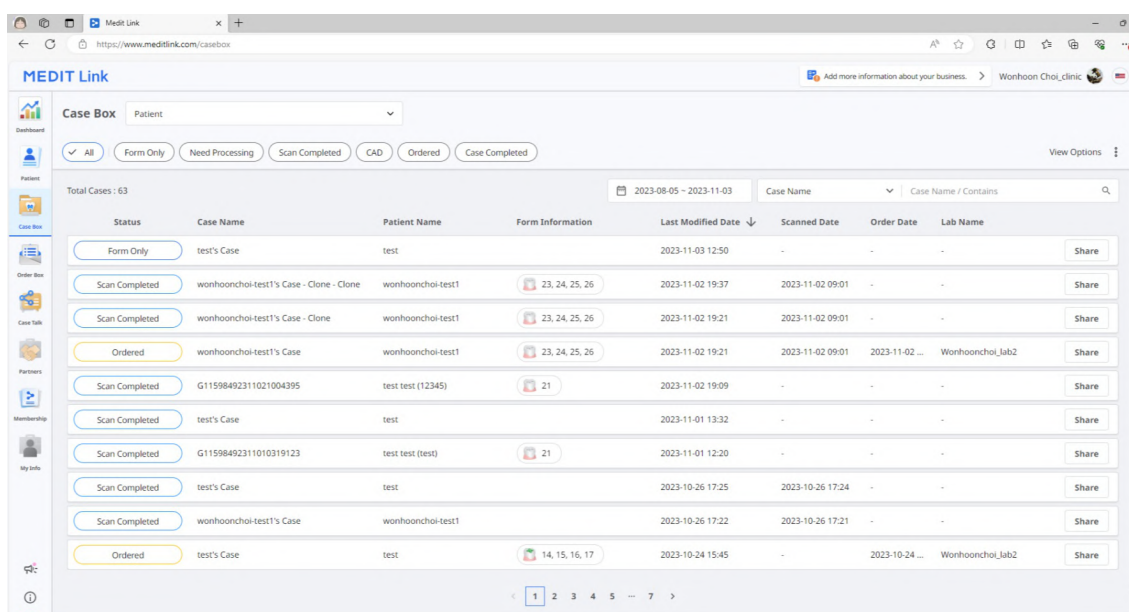
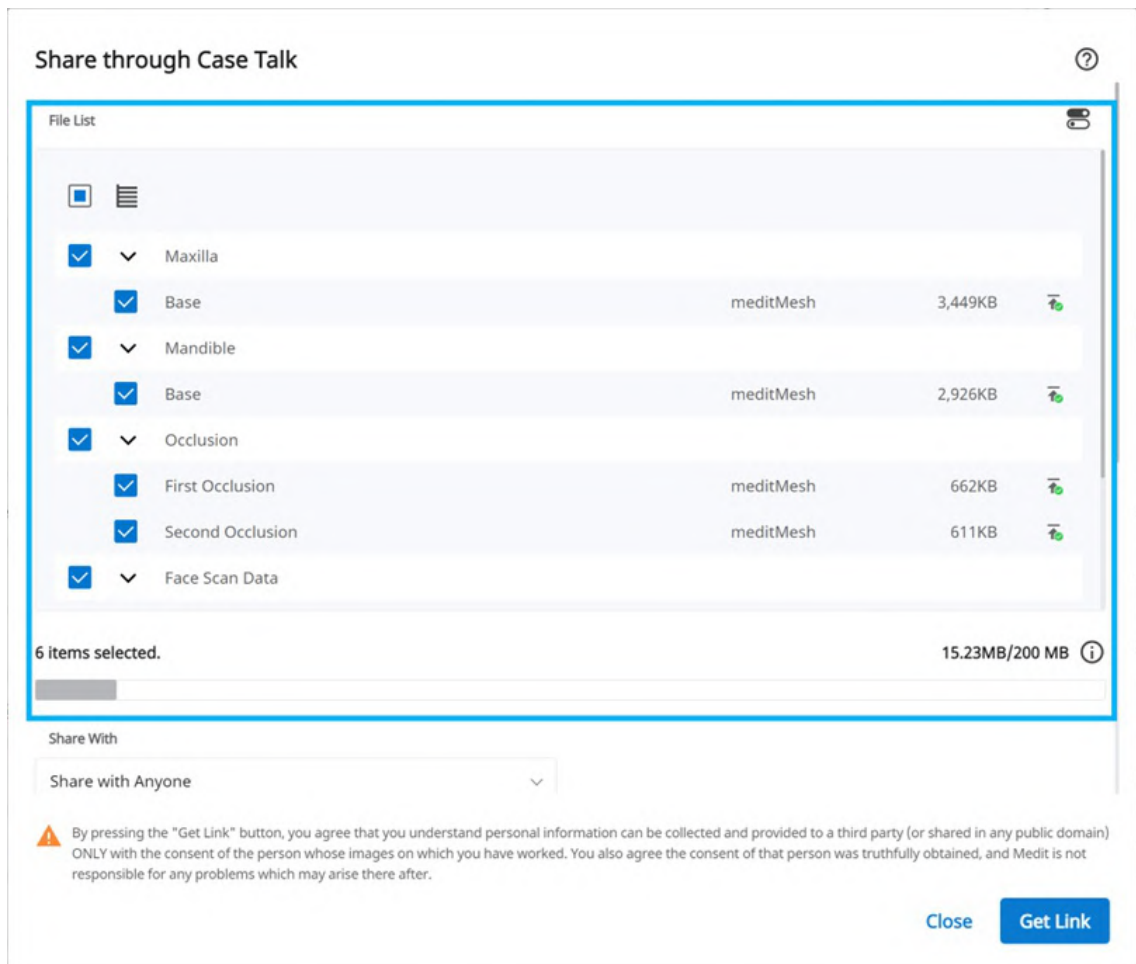
+

☰

❗ Create a shared link and share it with specific people.

Expiration Date

📅 2022-06-28

❗ You can extend the period up to 90 days from today.

PIN Code

☐ PIN Code

Enter PIN Code

❗ Set the PIN code using six digits for privacy protection.

☒ Share the patient's name

☒ Grant Download Permission [Premium Feature](#)

❗ If you check this option, all users who have access to the shared link can download the files.

❗ The premium feature is provided to all users for a limited time.
Feature availability is subject to change depending on the membership plan.

☒ Allow comments

❗ If you check this option, you will be able to communicate through the shared link.

Comment

Please enter your comment.

0/1024

⚠️ By pressing the "Get Link" button, you agree that you understand personal information can be collected and provided to a third party (or shared in any public domain) ONLY with the consent of the person whose images on which you have worked. You also agree the consent of that person was truthfully obtained, and Medit is not responsible for any problems which may arise there after.

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Invite Partner

Invite your desired partner to join Medit Link. When the lab accepts your invitation, a partnership is automatically formed.

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Membership

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The screenshot displays the MEDIT Link web interface. At the top, the 'Membership' tab is selected in the navigation bar. The main content area is divided into two sections. The left section, titled 'Medit Cloud Storage Plan', shows the '10 TB Plan' as the active subscription. It lists benefits: 'A paid plan', '10 TB of cloud storage', 'Efficient for users with many cases and various data', and 'One-month free trial included'. A 'Subscribe' button is visible. Below this, a '10 TB Plan Benefits' box shows a '10 TB Plan' with a code 'qqqq' and an expiration date of '2025-06-26'. The right section, titled 'Cloud Storage', shows a storage usage gauge at '3.46 GB / 10 TB'. It includes a 'Set Period' filter set to '2018-01-01 ~ 2025-02-17'. Below the gauge, there are checkboxes for 'Total Selected' (0.00 KB), 'Raw Data' (52.36 GB), and 'Result Files & Attachments' (3.46 GB). The footer contains links for Home, About, Terms & Conditions, Privacy Policy, and Copyright © Medit Co., along with a Korean text '스캐노 지평'.

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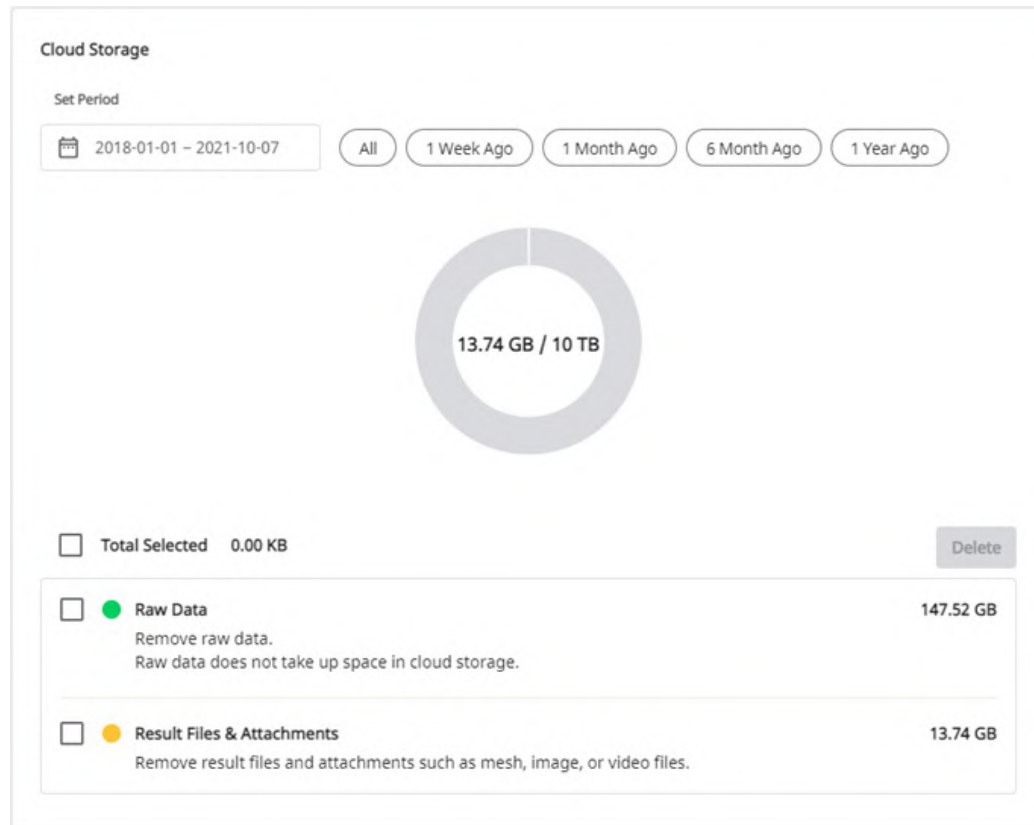
2 GB Plan **10 TB Plan**


Free (Promotion)

- ✓ A paid plan
- ✓ 10 TB of cloud storage
- ✓ Efficient for users with many cases and various data
- ✓ One-month free trial included

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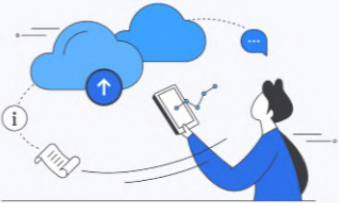10 TB Plan

9.99\$ USD / Monthly

- ✓ Subscribe now to keep your 10 TB Plan after the promotion expires.
- ✓ You will pay a subscription fee after the current promotion expires. If you don't wish to pay, you can unsubscribe before payment.

* VAT is not included.
* Learn more about the Medit Cloud Storage Plan in our Q&A [Here](#).

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←

Subscribe to Medit Cloud Storage 10TB Plan

\$0.00
Then \$9.99 per day after coupon expires

Medit Cloud Storage 10TB Plan for Clinics Medit Cloud Storage Plan Billed daily	\$9.99
Subtotal	\$9.99
1 Month Free for New Subscribers 100% off for a day	-\$9.99
Tax ⓘ	\$0.00
Total due today	\$0.00

Pay with card

Email yhgam0829@gmail.com

Card information

1234 1234 1234 1234    

MM / YY CVC 

Cardholder name

Full name on card

Country or region

South Korea ▼

Subscribe

By confirming your subscription, you allow Medit USA Inc to charge you for future payments in accordance with their terms. You can always cancel your subscription.

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Cancel Subscription

Are you sure you want to cancel your **10 TB Plan**?

✓ By confirming, you will automatically be switched to the free 2 GB Plan.
Enjoy full access to the features and storage of your current plan through to 2025-02-20.

ⓘ Please note that all uploaded files will be retained regardless of the plan you choose.
If your storage exceeds the 2 GB limit, you won't be able to upload new files.

*Click [here](#) to learn more about our membership plans.

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Confirm Subscription Cancellation

Your subscription will be canceled on 2025-02-20.

Stay Subscribed

Cancel Subscription

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Confirm Subscription Cancellation


Your subscription has been canceled.
You can use your current plan until 2025-02-20. After that, the cloud capacity will be limited to 2 GB.
If the used space exceeds the limit, data will not be backed up, and functions such as Order and Share will not be available.

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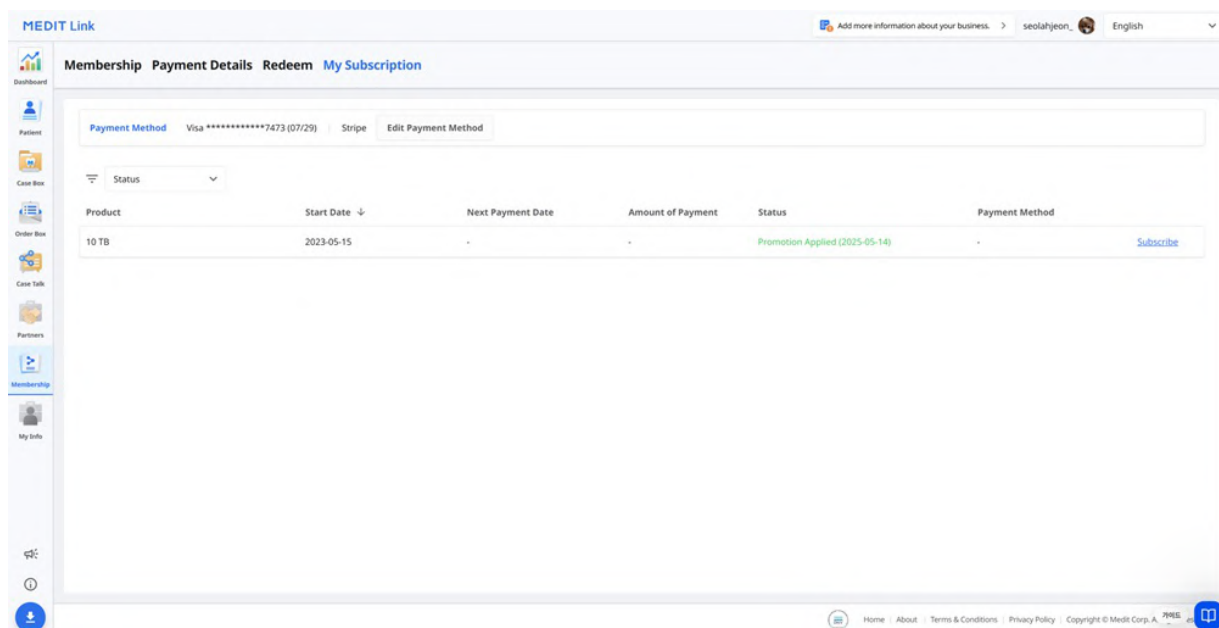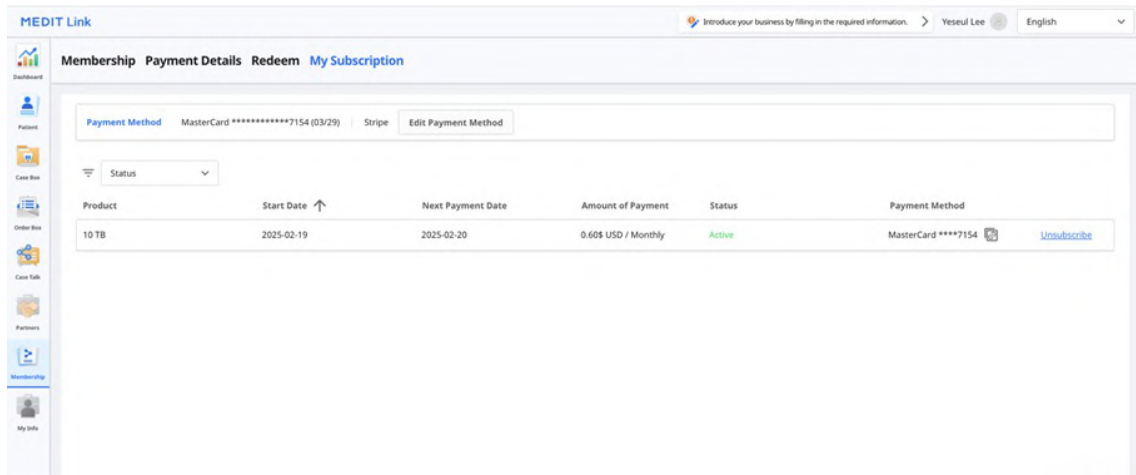
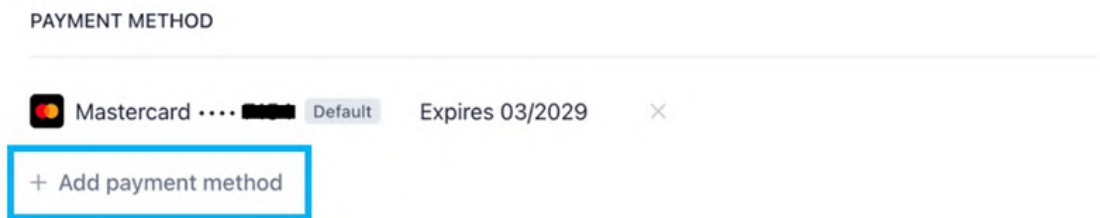
Country

South Korea

By providing your card information, you allow MEDIT USA Inc to charge your card for future payments in accordance with their terms.

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CURRENT SUBSCRIPTION

[DEV] Medit Cloud Storage 10TB Plan for Clinics
\$9.99 per day

[View details](#) ▾

Your subscription renews on February 20, 2025.

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MEDIT Link Add more information about your business. > seolah.jeon_dev_admin English

Membership Payment Details Redeem My Subscription

Register Redeem Code

Register

Redeem Code

All

	Product	Promotion	Code Number	Registration Date	From	To	Status	
4	Premium	4 Clinic Premium Promotion	2A6F-0125-6747-4A87	2022-03-16	2021-10-01	2021-10-05	Waiting	Cancel

1

Home About Terms & Conditions Privacy Policy Copyright © Medit Corp. GUIDES

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Register Redeem Code

You can register for the promotion or apply immediately.
 Register allows you to apply and activate when you need it, If applied immediately, it will start after the expiration of the active promotion.
 Periodic payments are stopped automatically.

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MEDIT Link

Introduce your business by filling in the required information.
Yeseul Lee
English

Dashboard
Patient
Case Box
Order Box
Case Talk
Partners
Membership
My Info

Membership
Payment Details
Redeem
My Subscription

All
Membership
Medit Apps

2024-11-21 ~ 2025-02-19

Payment Date	Product	Payment Type	Duration	Details	Amount of Payment	Status
2025-02-19	10 TB	Subscription	2025-02-19 ~ 2025-02-20	Monthly payment	0.00 \$ (USD)	Payment Completed

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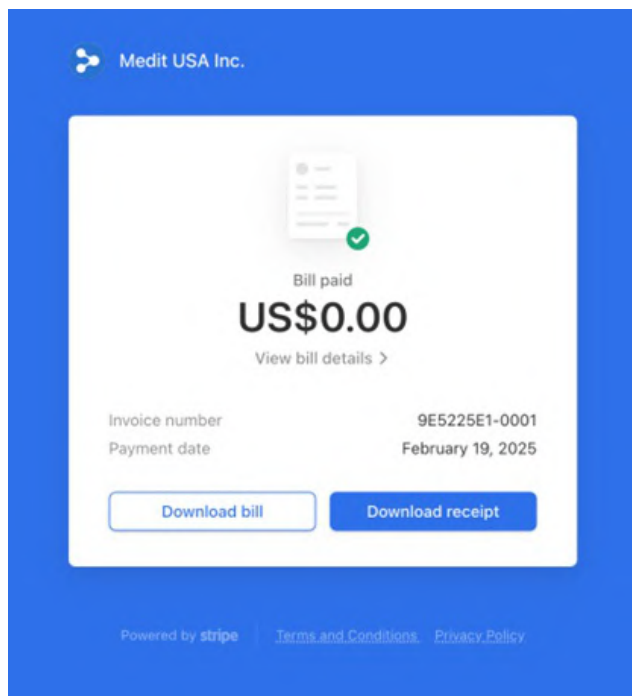
X Close Bill Details

Paid on February 19, 2025

SUMMATION				
Target	yeseul.lee.medit.clinic@gmail.com			
source	Medit USA Inc.			
bill	#9E5225E1-0001			
ITEM				
FEBRUARY 19 - FEBRUARY 20, 2025				
Medit Cloud Storage 10TB Plan for Clinics		US\$0.60		
Quantity 1 piece				
subtotal		US\$0.60		
1 Month Free for New Subscribers(100% 할인)				
		-US\$0.60		
Payment amount		US\$0.00		
Amount paid		US\$0.00		
Remaining amount		US\$0.00		

Do you have any questions?Contact Medit USA Inc.

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Country/Region

Korea, Republic of ▼

📘 Your chosen language is applied to the email sent by Medit Link and the notification.

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
Edit

Cover and Logo

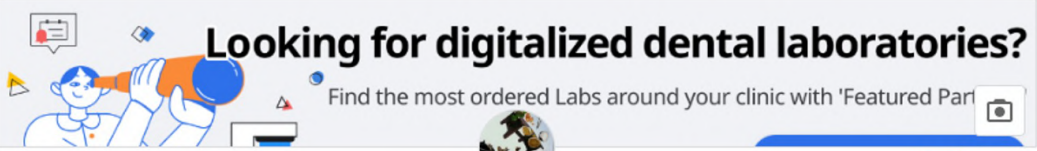Looking for digitalized dental laboratories?
Find the most ordered Labs around your clinic with 'Featured Partners'

1

2

3

Services

Select

▼

Specialties (choose up to 3)

Select

▼

 Select the services your business specializes in from the list.

About Us

e.g., Welcome to OOO Dental Clinic.
We are a licensed dental clinic situated in California specializing mainly in aesthetic restorations and implantology.

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Notice
NFT Service Termination Notice
2023-04-19

To our users,

We want to thank all Medit Link users who have supported and used our NFT service.

Unfortunately, we must inform you that this service will be terminated following the upcoming program update.

Please see the details below for more information:

- **Feature in Question**
 - Share via NFT (in Medit Link Case Talk)
- **Effective Date**
 - Starting from the update in the first half of 2023 (scheduled for April 27)

※ The update schedule is subject to change.

We'll continuously work to provide better service.

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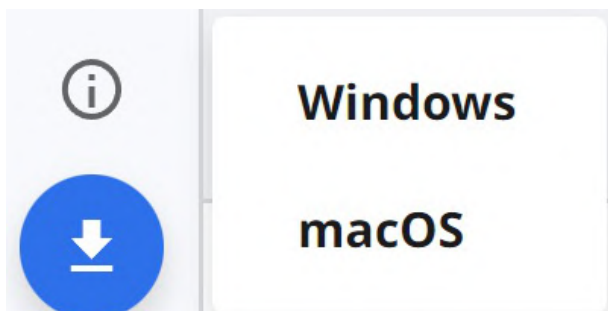
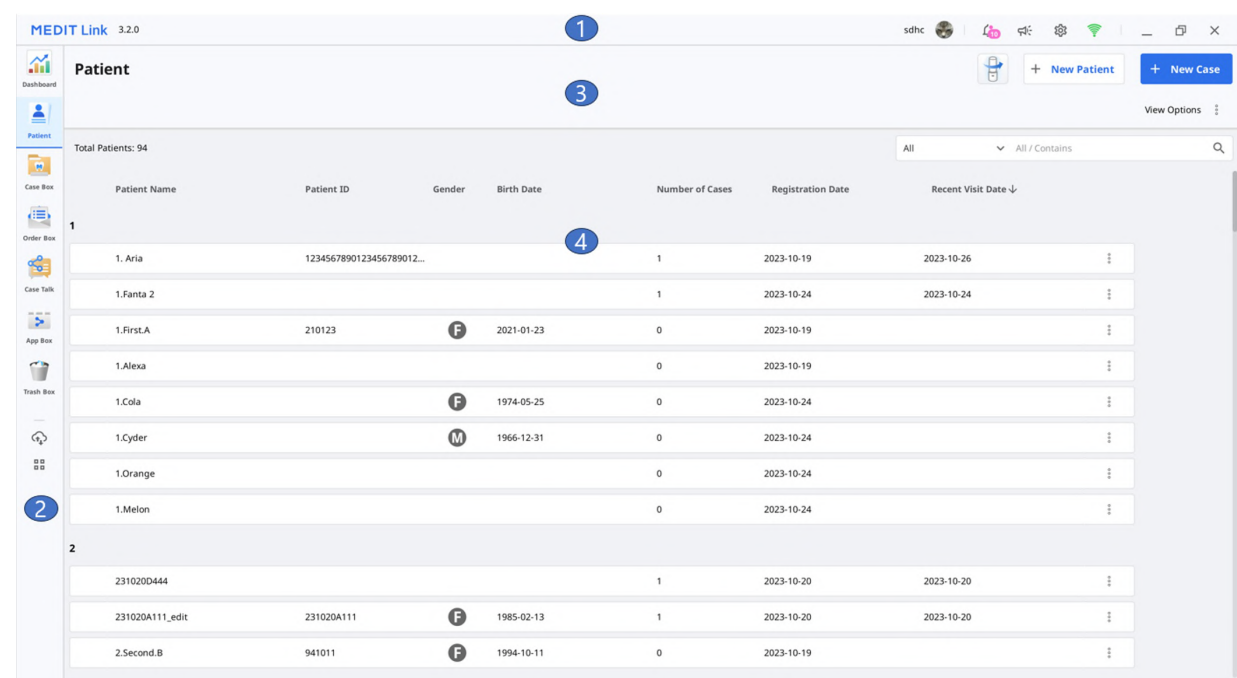
2

Delivery

Wonhoonchoi_Lab_Dev has dispatched the ordered package.

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Order

Wonhoonchoi_Lab_Dev has accepted the order.

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Order

Wonhoonchoi_Lab_Dev has accepted the order.

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Order

Wonhoonchoi_Lab_Dev has accepted the order.

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Order

Wonhoonchoi_Lab_Dev has accepted the order.

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Order

Wonhoonchoi_Lab_Dev has accepted the order.

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Work

Wonhoonchoi_Clinice_Dev_Member registered a new case test's Case.

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Work

Wonhoonchoi_Clinice_Dev_Member registered a new case test's Case.

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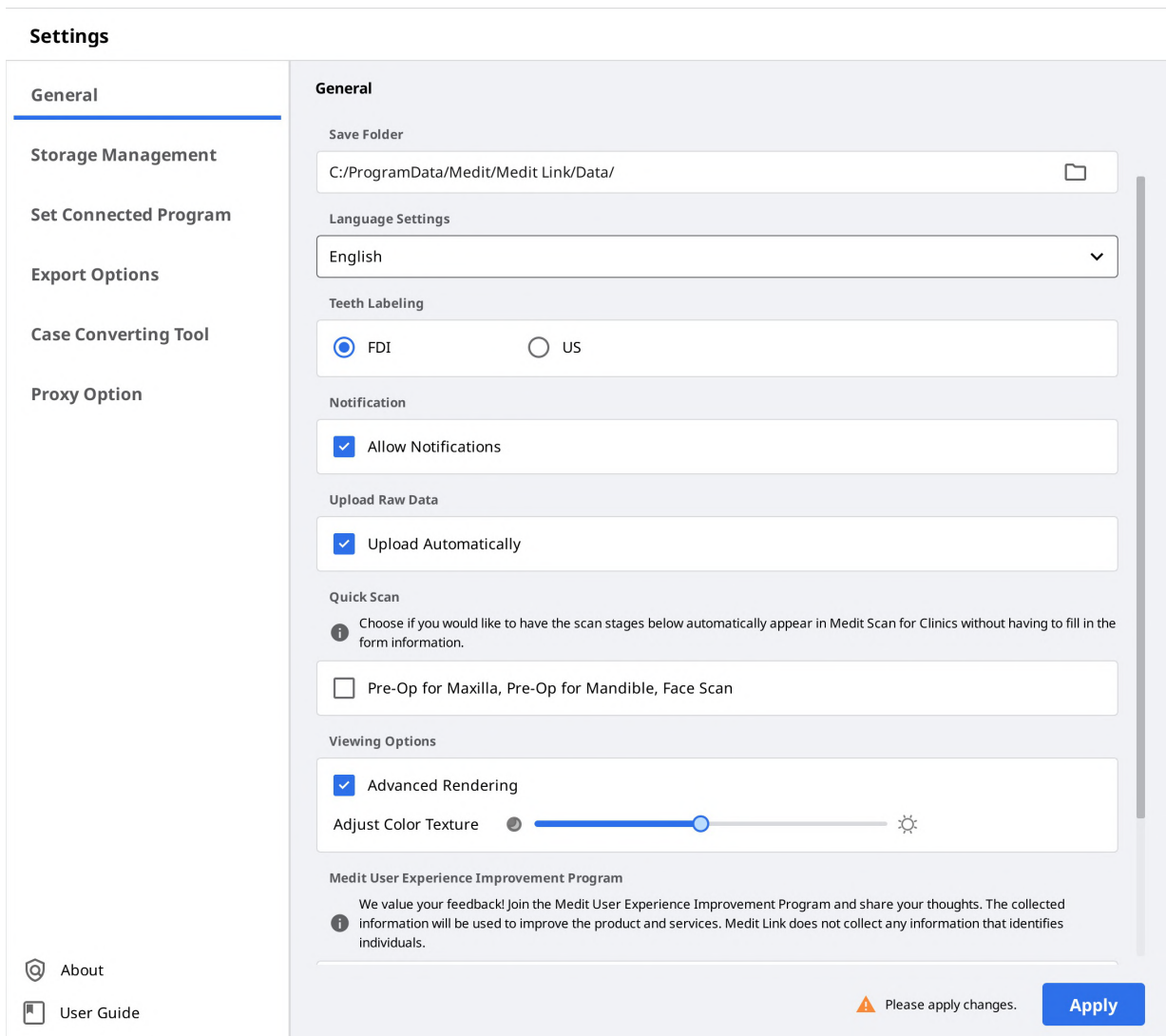
Network speed is good.

Network speed is fast enough.

Network speed is poor. Can be switched to offline mode automatically.

Either network is not connected or Medit Link is in offline mode.

Synchronization is in progress. Some case lists may not match.

Synchronization failed. The case lists may not match.

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Delete

Raw Data

Remove raw data created when using scan programs.

26.89 GB

Result Files & Attachments

Remove result files and attachments such as mesh, image, or video files.

2.54 GB

Others

It shows space used by other programs.

666.76 GB

Cloud Storage

Total Selected 0.00 KB

Delete

Raw Data

Remove raw data created when using scan programs. Raw data are not counted in cloud storage.

837.95 GB

Result Files & Attachments

Remove result files and attachments such as mesh, image, or video files.

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Compress & Combine

☒ Compress Files
☒ Combine Individual Mesh

Options when same folder or file already exists.

☒ Always ask when exporting
☐ Export a new file and overwrite the existing file. The file currently stored will be lost.
☐ Create a folder and export the files

Rule for Folder Name

☒ Case Creation Date
Scanned Date
Order Date

☒ Patient Name
Case Name
Partner Name

☐ Patient ID
☐ Order ID

Folder Name

Case Creation Date - Patient Name

Rule for File Name

☒ Case Creation Date
Scanned Date
Order Date

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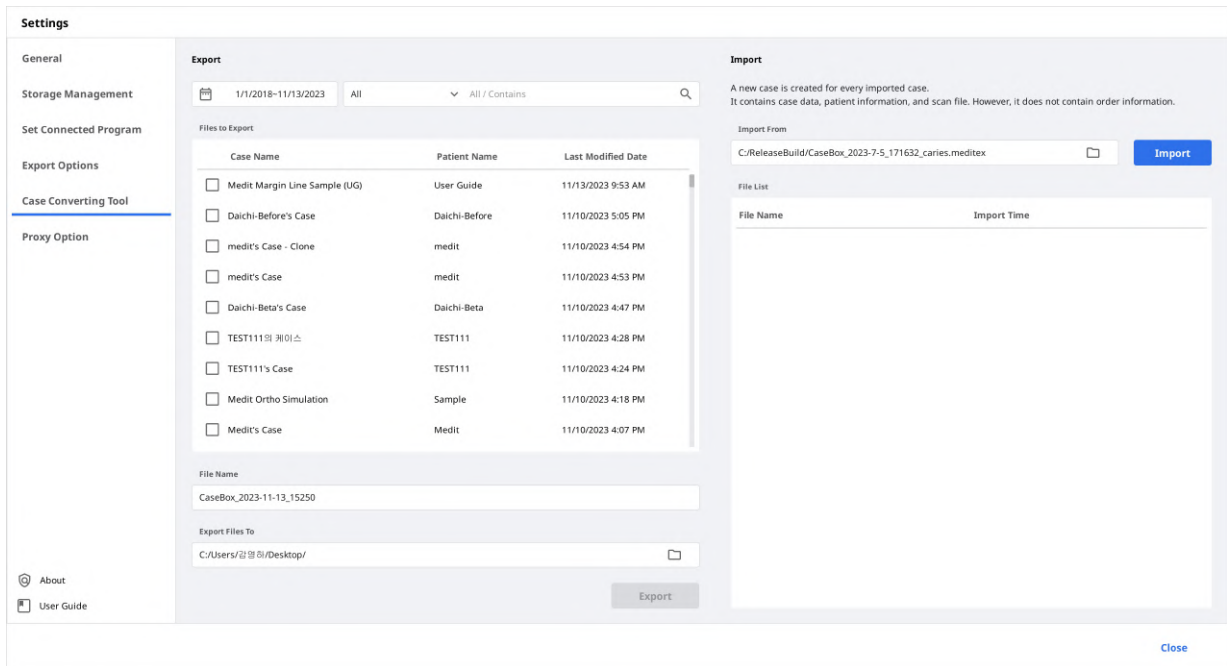
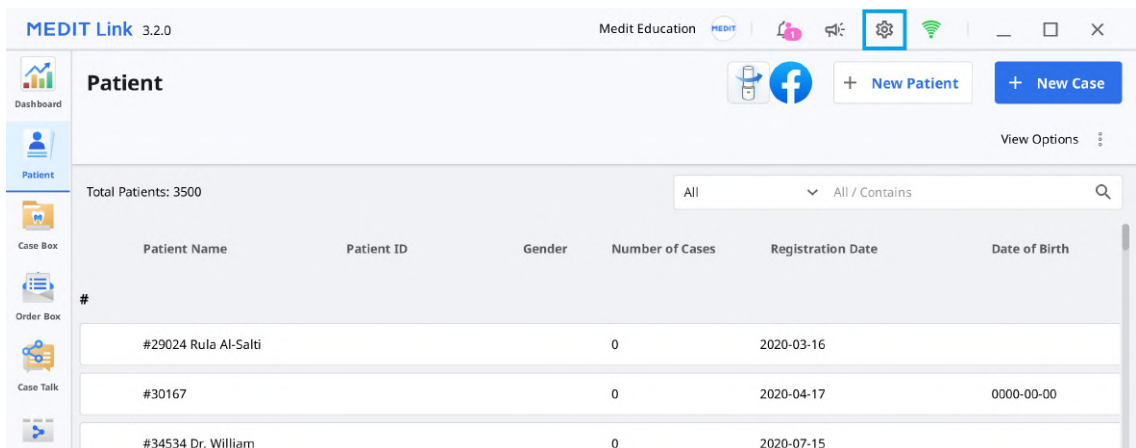
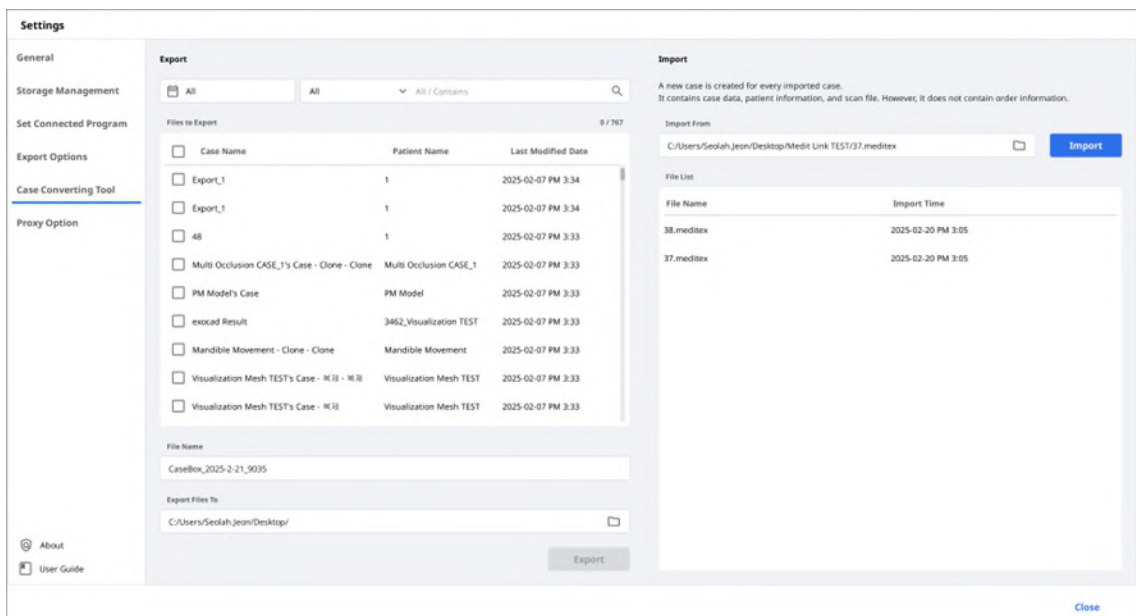
Question

Exporting successfully completed.
Do you want to open the folder containing the exported file?

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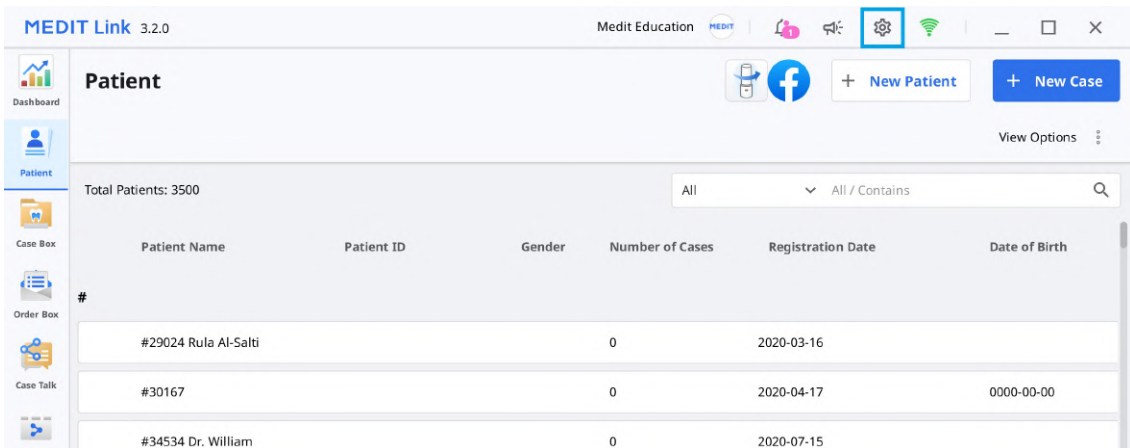
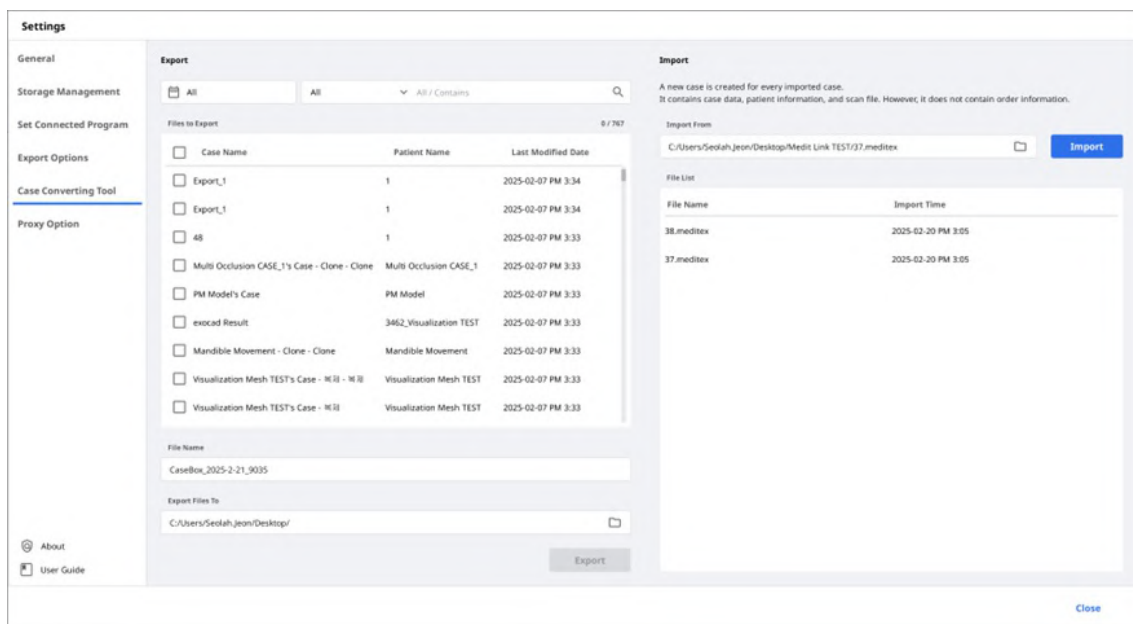
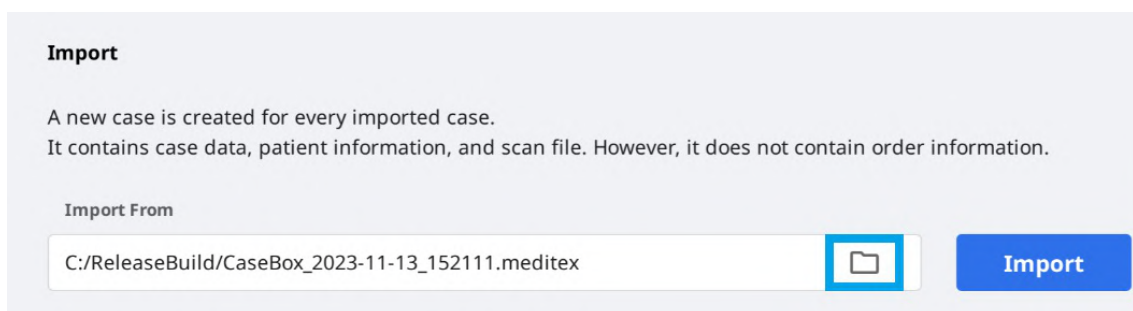
Import

A new case is created for every imported case.
It contains case data, patient information, and scan file. However, it does not contain order information.

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Case Name Case Name / Contains

Total Cases: 210 ☒ No Date Data Type All

☐	Case Name	Patient Name	Scanned Date↓	Size	Status
☐	BLATRIX Charline's Case	BLATRIX Charline	2025-01-06 PM 6:28	729.30 MB	Upload Required
☐	BLATRIX Charline's Case	BLATRIX Charline	2025-01-06 PM 6:28	729.30 MB	Upload Required
☐	BLATRIX Charline's Case - 복제	BLATRIX Charline	2025-01-06 PM 6:28	733.35 MB	Upload Required
☐	BLATRIX Charline's Case - Clone	BLATRIX Charline	2025-01-06 PM 6:28	734.01 MB	Upload Required
☐	test의 케이스	Patient TEST	2023-04-26 PM 7:58	1.25 GB	Upload Required
☐	i700 HD's Case	i700 HD	2023-04-24 PM 5:23	911.24 MB	Upload Required
☐	Mandibular Movement	MEDIT	2023-04-24 PM 1:59	1.88 GB	Upload Required
☐	Mandibular Movement Case	MEDIT	2023-04-24 PM 1:59	1.75 GB	Upload Required
☐	Mandibular Movement Case - Clo...	MEDIT	2023-04-24 PM 1:59	1.75 GB	Upload Required
☐	서라_Mandibular Movement_TEST	Mandibular Movement_export_test	2023-04-10 PM 3:39	34.73 MB	Upload Required
☐	서라_Mandibular Movement_TEST...	Mandibular Movement_export_test	2023-04-10 PM 3:39	731.43 MB	Upload Required

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
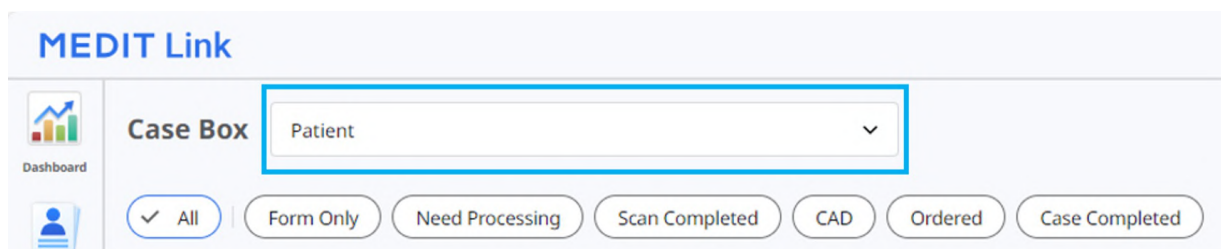
The screenshot shows the MEDIT Link interface. On the left is a sidebar with a 'Dashboard' icon and a user profile icon. The main area has a 'Case Box' dropdown menu currently set to 'Patient'. Below this are several filter buttons: 'All' (selected), 'Form Only', 'Need Processing', 'Scan Completed', 'CAD', 'Ordered', and 'Case Completed'.

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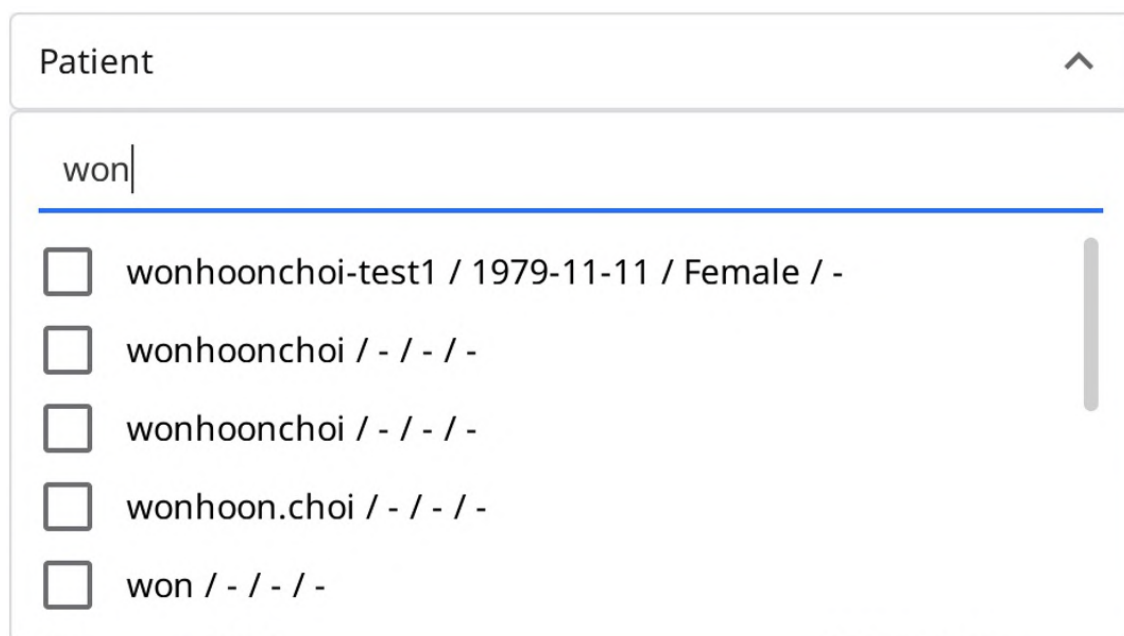This screenshot shows the search results dropdown menu. The search input field contains 'won'. The results list shows five entries, each with a checkbox and a patient identifier followed by slashes and dashes: 'wonhoonchoi-test1 / 1979-11-11 / Female / -', 'wonhoonchoi / - / - / -', 'wonhoonchoi / - / - / -', 'wonhoon.choi / - / - / -', and 'won / - / - / -'.

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MEDIT Link 3.2.0 Wonhoon Choi Clinic

Case Box wonhoonchoi / - / - / - + 2

☒ All
 ☐ Form Only
 ☐ Need Processing
 ☐ Scan Completed
 ☐ CAD
 ☐ Ordered
 ☐ Case Completed

Total Cases: 20 2018-01-01-2023-10-24 All / Contains

Status	Patient Name	Case Name	Form Information	Last Modified Date	Scan Date	Lab Name	Order Date
Form Only	wonhoonchoi	originalCase	-	2023-10-05 오전 8:10			
Form Only	wonhoonchoi	wonhoonchoi's Case	-	2023-10-05 오전 8:09			
Scan Completed	wonhoonchoi	FileViwerTest_1	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-06-15 오전 8:14	2023-02-09 오전 8:35		
Scan Completed	wonhoonchoi	FileViwerTest_1 - Clone	2-Inlay / 3-Inlay / 4-Inlay / 5-Inlay / 6-Inlay / 8...	2023-05-25 오후 2:05	2023-02-09 오전 8:35		
Form Only	wonhoonchoi	FileViwerTest_1 - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-03-15 오후 3:06			
Scan Completed	wonhoonchoi	FileViwerTest_1 - Clone - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-03-09 오후 4:41	2023-02-09 오전 8:35		
Need Processing	wonhoonchoi	FileViwerTest_1 - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-02-14 오전 8:13	2023-02-09 오전 8:35		
Scan Completed	wonhoonchoi	FileViwerTest_1 - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-02-10 오전 11:27	2023-02-09 오전 8:35		
Form Only	wonhoonchoi	testsamplecase	-	2023-02-10 오전 7:39			
Scan Completed	wonhoonchoi	Multi Occlusion	-	2023-02-09 오전 9:26	2023-02-09 오전 9:01		
Scan Completed	wonhoonchoi	FileViwerTest_1 - Clone - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-02-09 오전 8:50	2023-02-09 오전 8:35		
Case Completed	wonhoonchoi	wonhoonchoi's Case	12-Inlay / 13-Inlay / 14-Inlay	2022-10-13 오후 10:30	2022-10-11 오후 5:16		
Case Completed	wonhoonchoi	wonhoonchoi's Case	-	2022-10-13 오후 10:28			

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
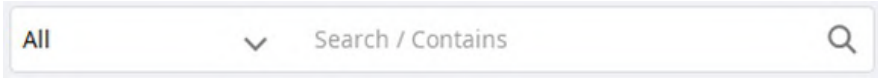
A search bar with a dropdown menu currently showing 'All'. To the right of the dropdown is the text 'Search / Contains' and a magnifying glass icon.

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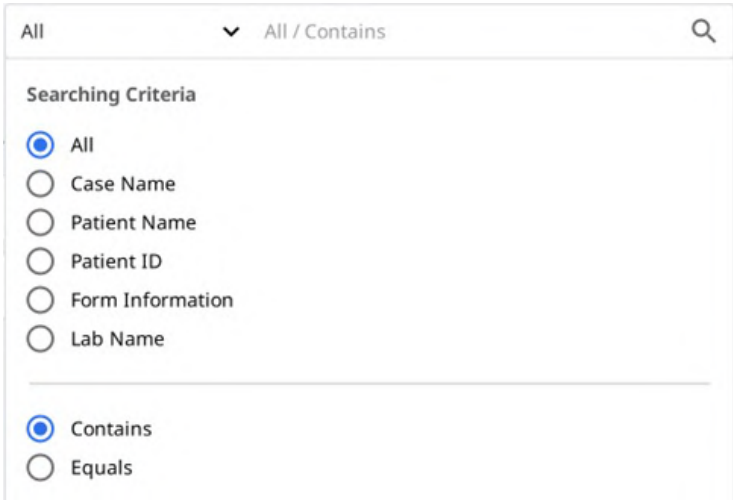
An expanded search criteria dropdown menu. It has a header 'Searching Criteria'. Under this header, there are two groups of radio buttons. The first group lists search criteria: 'All' (selected), 'Case Name', 'Patient Name', 'Patient ID', 'Form Information', and 'Lab Name'. The second group lists logical operators: 'Contains' (selected) and 'Equals'.

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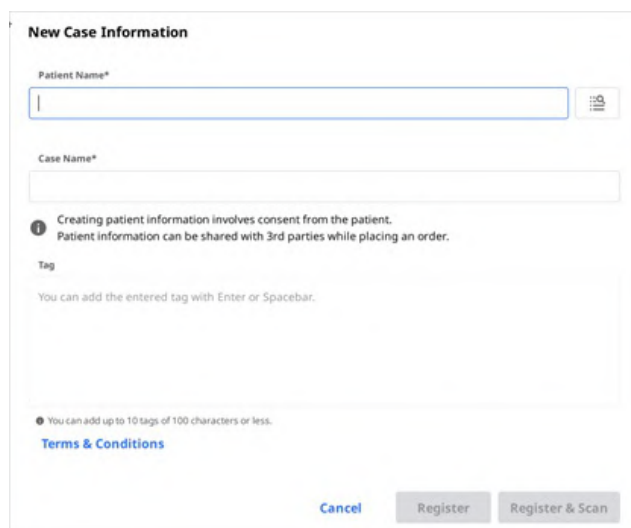
The 'New Case Information' form contains the following fields and elements:

- Patient Name***: A text input field with a search icon on the right.
- Case Name***: A text input field.
- Information**: A note stating 'Creating patient information involves consent from the patient. Patient information can be shared with 3rd parties while placing an order.'
- Tag**: A text area with the instruction 'You can add the entered tag with Enter or Spacebar.'
- Footer**: A note 'You can add up to 10 tags of 100 characters or less.', a link to 'Terms & Conditions', and three buttons: 'Cancel', 'Register', and 'Register & Scan'.

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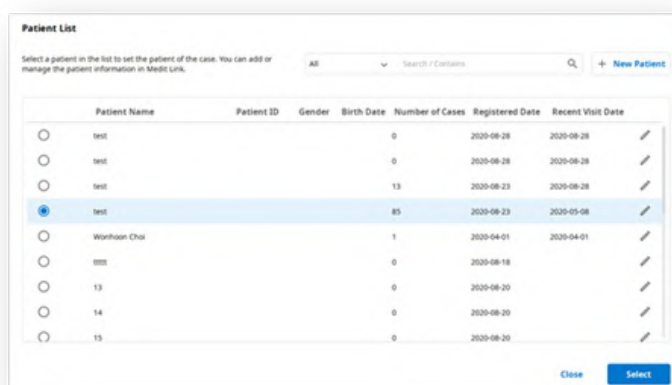
The 'Patient List' dialog box displays a table of existing patients. The table has the following columns: Patient Name, Patient ID, Gender, Birth Date, Number of Cases, Registered Date, and Recent Visit Date. The first four rows have 'test' as the patient name. The fifth row has 'Wanhsun Choi'. The last three rows have IDs 13, 14, and 15. The fourth row (with ID 85) is selected. There are radio buttons to the left of each row and an edit icon to the right. At the bottom right are 'Close' and 'Select' buttons.

	Patient Name	Patient ID	Gender	Birth Date	Number of Cases	Registered Date	Recent Visit Date	
☐	test				0	2020-09-28	2020-09-28	
☐	test				0	2020-09-28	2020-09-28	
☐	test				13	2020-09-23	2020-09-28	
☒	test				85	2020-09-23	2020-05-08	
☐	Wanhsun Choi				1	2020-04-01	2020-04-01	
☐	test				0	2020-09-18		
☐	13				0	2020-09-20		
☐	14				0	2020-09-20		
☐	15				0	2020-09-20		

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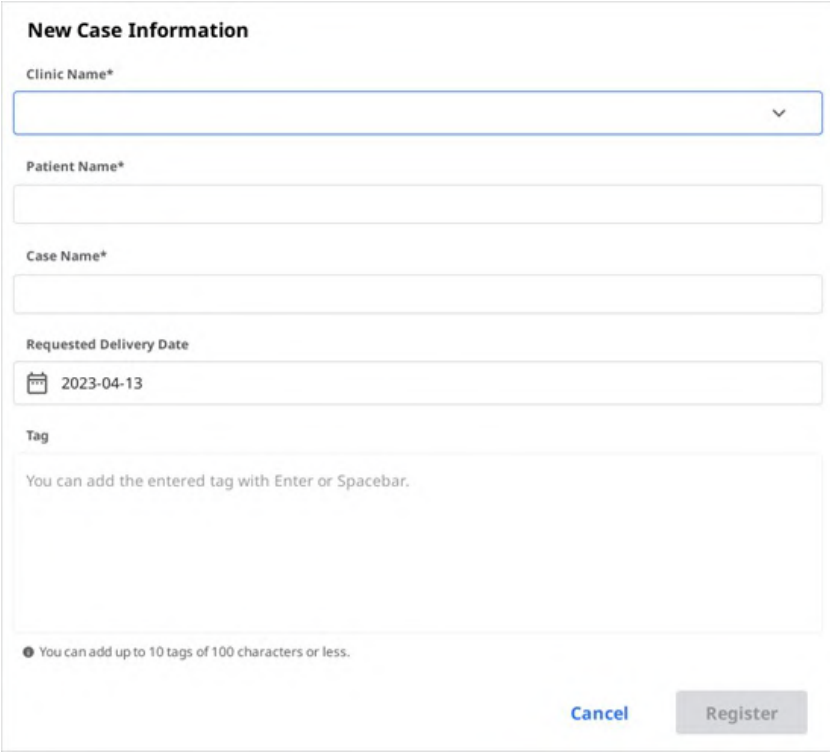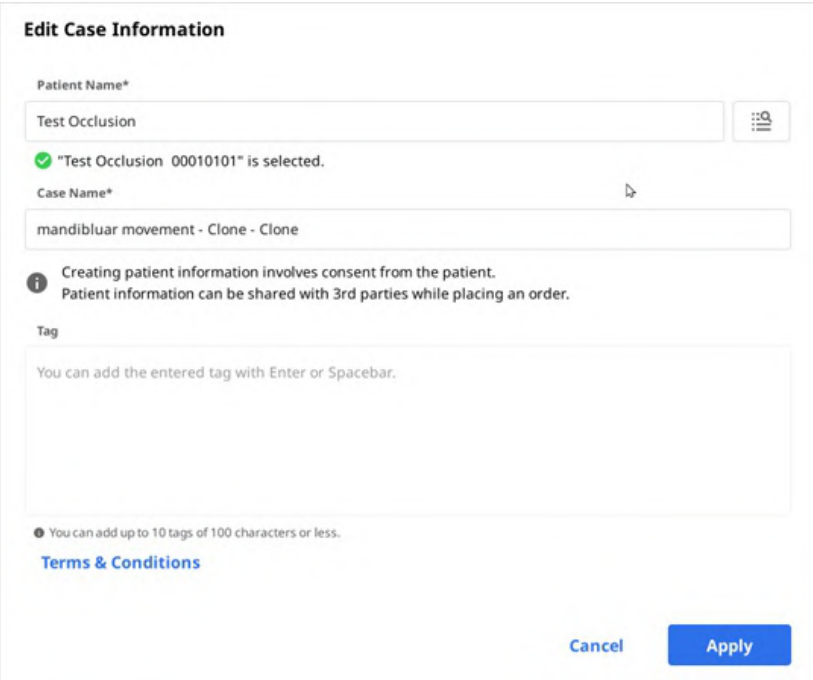
Case Name*

mandibluar movement - Clone - Clone

Creating patient information involves consent from the patient.
Patient information can be shared with 3rd parties while placing an order.

Tag

You can add the entered tag with Enter or Spacebar.

You can add up to 10 tags of 100 characters or less.

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Patient List

Select a patient in the list to set the patient of the case. You can add or manage the patient information in Medit Link.

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2022-04-18

Tag

You can add the entered tag with Enter or Spacebar.

● You can add up to 10 tags of 100 characters or less.

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Options **File List**

Export Files To
C:/Users/Medit-Wonhoonchoi/Desktop/test/

Folder Name
Multi Occlusion - Scan - Before Scanning

File Name
Multi Occlusion - Scan - Before Scanning

CAD Type
Maestro

Mesh Format
Convert mesh files to the below formats
☐ meditMesh ☐ OBJ ☒ **PLY** ☐ STL
☒ Do not convert the attached files

Image Format
☐ Convert images to the below format
☐ JPG ☐ JPEG ☒ **PNG** ☐ BMP

☐ Compress Files
☐ Combine Individual Mesh
☐ Export All Occlusal Relationship
☐ Save Case Information as PDF [Preview](#)

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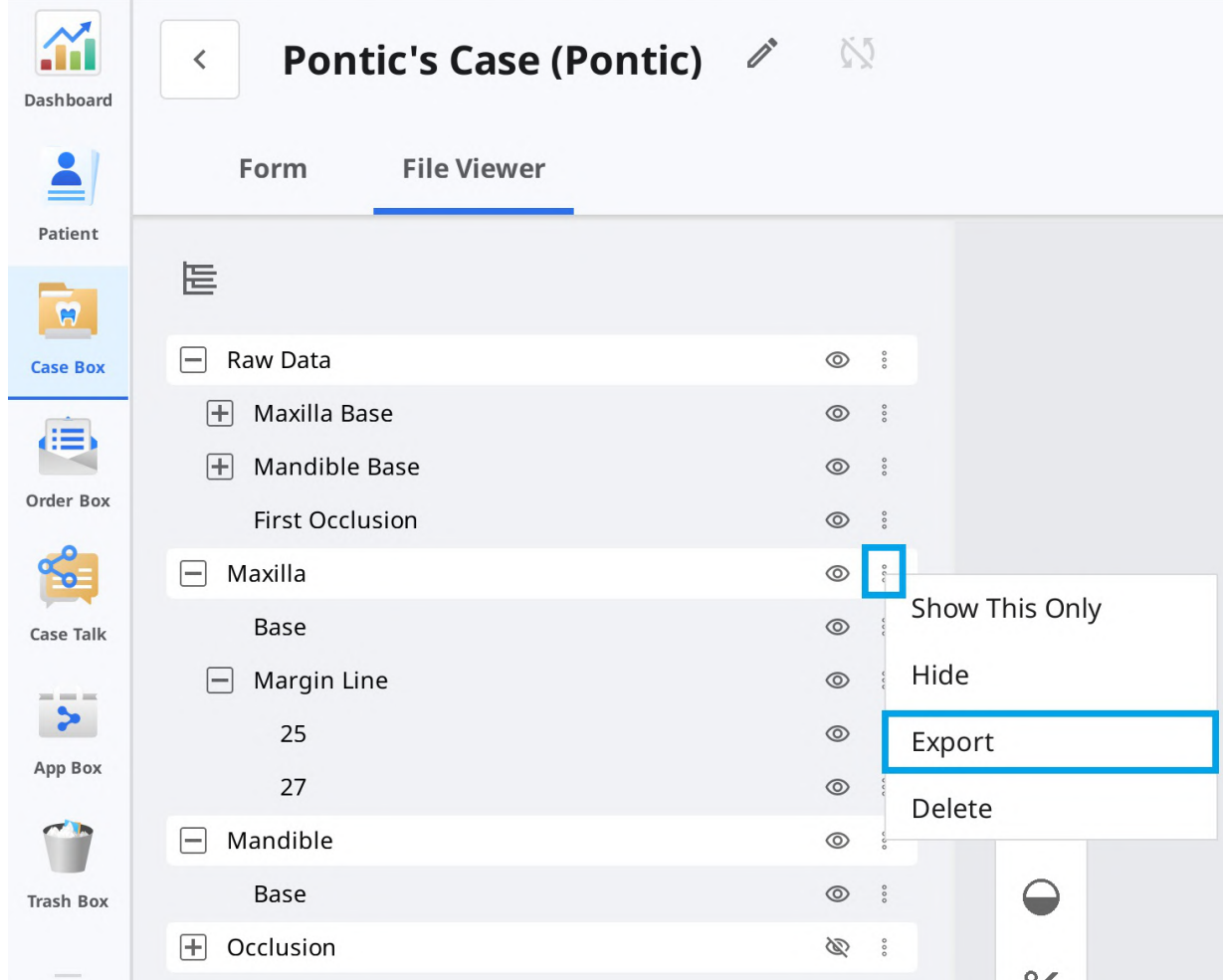
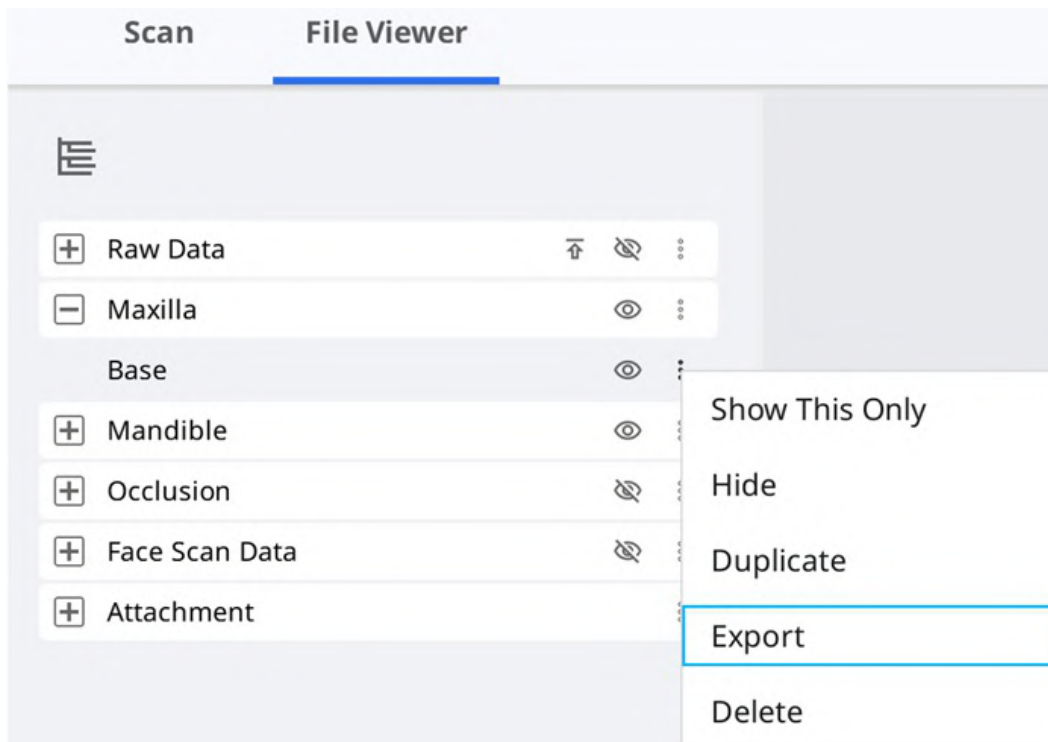
C:/Users/Medit-Wonhoonchoi/Desktop/test/

Folder Name

Multi Occlusion - Scan - Before Scanning

File Name

Multi Occlusion - Scan - Before Scanning

CAD Type

Maestro

Mesh Format

Convert mesh files to the below formats

☐ meditMesh
☐ OBJ
☒ **PLY**
☐ STL

☒ Do not convert the attached files

Image Format

☐ Convert images to the below format

☐ JPG
☐ JPEG
☒ **PNG**
☐ BMP

☐ Compress Files
☐ Combine Individual Mesh
☐ Export All Occlusal Relationship
☐ Save Case Information as PDF [Preview](#)

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Share through Case Talk

Share with Anyone

Create a shared link and share it with anyone.

Expiration Date: 2022-04-27

PIN Code: ☐ PIN Code

You can extend the period up to 90 days from today.

Set the PIN code using six digits for privacy protection.

☒ Share the patient's name

☒ **Grant Download Permission** Premium Feature

If you check this option, all users who have access to the shared link can download the files.

The premium feature is provided to all users for a limited time. Feature availability is subject to change depending on the membership plan.

☒ Allow comments

If you check this option, you will be able to communicate through the shared link.


Comment: Please enter your comment.

0/1024

By pressing the "Get Link" button, you agree that you understand personal information can be collected and provided to a third party (or shared in any public domain) ONLY with the consent of the person whose images on which you have worked. You also agree the consent of that person was truthfully obtained, and Medit is not responsible for any problems which may arise there after.

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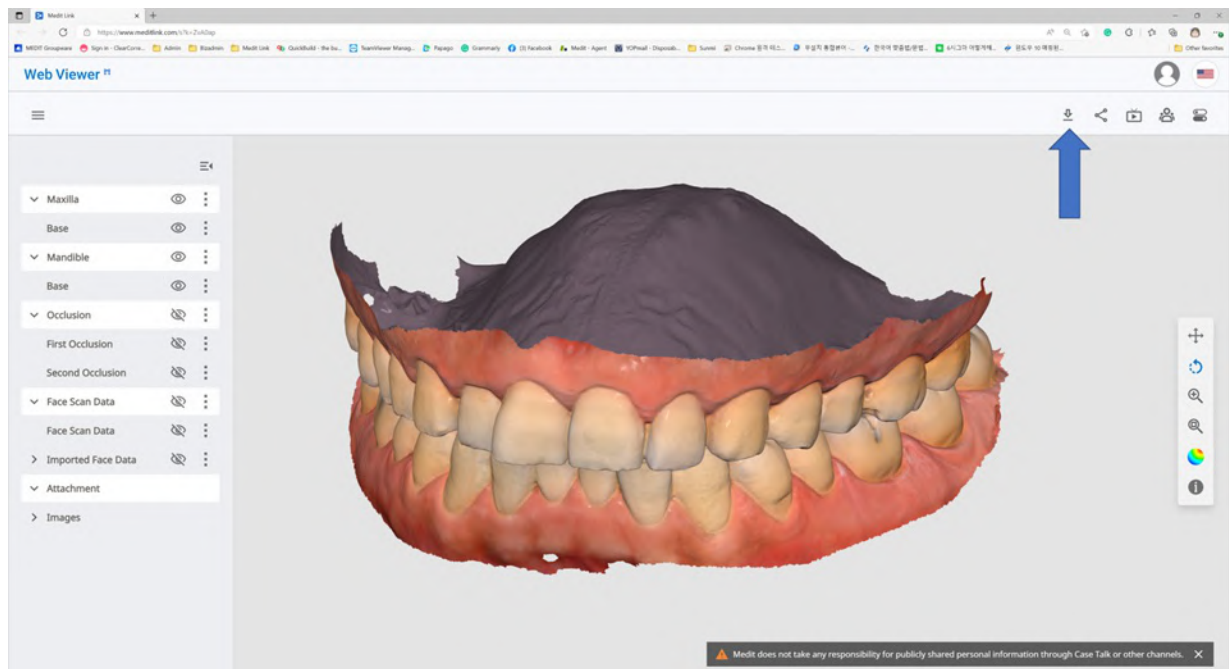Download Premium Feature

The premium feature is provided to all users for a limited time.
Feature availability is subject to change depending on the membership plan.

Mesh Format

Convert mesh files to the below formats

☐ meditMesh ☐ OBJ ☒ PLY ☐ STL

i The selected format will be applied to the result files created in the Medit Link application.

i The attachment files will not be converted but will keep their original format.

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The screenshot shows the 'Order' form with the following fields and options:

- Order** (Section Header)
- General** (Selected Tab), **Order List**, **File List**
- Partner**: WonhoonChoi_Lab_Stage
- Partner Information**:
 - Icon: A person at a computer.
 - WonhoonChoi_Lab_Stage**
 - Seoul, KR
 - [02835] 23 Incheon-ro 22-gil, Seongbuk-gu, Seoul, Korea
 - 01026088546
- Case Name**: mandibular movement - Clone - Clone
- Patient Name**: Test Occlusion
- ☒ **Share patient information**
- Requested Delivery Date**: 2023-04-13
- Delivery Time**: 오전 8:15 (Range from AM 8:00 to PM 8:00)
- Memo**: (Empty text area)
- Tags**:
 - You can add the entered tag with Enter or Spacebar.
 - You can add up to 10 tags of 100 characters or less.
- ☒ **Include all occlusion relationships in the order**
- ☒ **Include Tags In The Order**
- Buttons**: Cancel, OK

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The screenshot shows a web interface for a dental clinic. At the top, there's a header with a back arrow, the text 'test's Case / Wonhoonchoi-Clinic', and two buttons: 'Reject' and 'Confirm Order'. The 'Confirm Order' button is highlighted with a red rectangle. Below the header, there are tabs for 'Order' and 'File Viewer'. At the bottom right, there's a 'Pending' status indicator.

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The screenshot shows a web interface for a dental clinic. At the top, there's a header with a back arrow, the text 'test's Case (test)', and a 'Cancel' button highlighted with a red rectangle. Below the header, there are tabs for 'Order' and 'File Viewer'. At the bottom right, there's a 'Pending' status indicator.

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The screenshot shows a web interface for a dental clinic. At the top, there's a header with a back arrow, the text 'test's Case / Wonhoonchoi-Clinic', and two buttons: 'Reject' and 'Confirm Order'. The 'Reject' button is highlighted with a red rectangle. Below the header, there are tabs for 'Order' and 'File Viewer'. At the bottom right, there's a 'Pending' status indicator.

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The screenshot shows the MEDIT Link 3.2.0 interface. The main window displays the details of an order titled 'test's Case - Clone / wh-clinic-ds...'. The order details include: Order ID (13455074), Patient Name (***), Partner Name (wh-clinic-dso-3), Scanned Date (2023-08-04 09:15:58), Order Date (2023-10-26 09:43:38), Requested Delivery Date (2023-11-02), and Expected Delivery Date (2023-11-02). There is also a Memo field and a Tags field. On the right side, there is a table titled 'Order List' showing a list of orders with columns: No., Type, Method, Material, and Shade. The table contains four rows of data.

No.	Type	Method	Material	Shade
3	Inlay	-	Multi-Layered Zirconia	None Shade
4	Inlay	-	Multi-Layered Zirconia	None Shade
5	Inlay	-	Multi-Layered Zirconia	None Shade
6	Inlay	-	Multi-Layered Zirconia	None Shade

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Shipment Registration

Select a delivery method and enter the required information to register.
The information you enter is saved in the Shipment tab and can be checked at any time.

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Case Name: mandibular movement - Clone - Clone Patient Name: Test Occlusion

☒ Share patient information

Requested Delivery Date: 2023-04-13 Delivery Time: 오전 8:15 (AM 8:00 to PM 8:00)

Memo

Tags
You can add the entered tag with Enter or Spacebar.

☒ Include all occlusion relationships in the order
☒ Include Tags In The Order

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Order **File Viewer**

Order ID: 10413603

Patient Name: test Partner Name: Wonhoonchoi-Clinic

Scanning Date: 2022-09-29 오전 7:56 Order Date: 2022-09-29 오전 8:13

Requested Delivery Date: 2022-10-01 오전 10:30 **Expected Delivery Date: 2022-10-01 오전 10:30** (with edit icon)

Memo
<Scanbody Library>
16 : ODENTIS_Custom_(JNMR) ver 1.41 - Custom_Mini_Extra Long_16mm - N-HEX_T10
15 : ODENTIS_Custom_(JNMR) ver 1.41 - Custom_Mini_Extra Long_16mm - N-HEX_T10
14 : ODENTIS_Custom_(JNMR) ver 1.41 - Custom_Mini_Extra Long_16mm - N-HEX_T10

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Selected Cases: 4

Case Completed

	Status		Case Name
☐	Ordered	test	test's Case
☒	Scan Completed	wonhoonchoi-test1	wonhoonchoi-test1's Case
☒	Scan Completed	test test(1234)	G11598492310180503424
☒	Scan Completed	test	test's Case
☒	Scan Completed	test	test's Case

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Tag

Doctorname × Specialmaterial × SpecialComment ×

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Patient Name*

Case Name*

*Creating patient information involves consent from the patient.
Patient information can be shared with 3rd parties while placing an order.*

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Memo

Tags

Doctorname × Specialmaterial × SpecialComment ×

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☒ Include all occlusion relationships in the order

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