USER GUIDE

Medit Link

Table of Contents

Medit Link

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Administrator

Dentist (Manager)

All Medit Link functions are available.

- ✓ Order
- ✓ Member Management

Register

Administrator

Representative/Leader (Manager)

All Medit Link functions are available.

- ✓ Order Management
- ✓ Member Management

Register

Member

Dental Hygienist and Dental Staff

In order to use Medit Link, join your group as a member. Case management and basic functions will be available after you join the group.

- ✓ Oral Scanning

Register

Member

Laboratory Employees

In order to use Medit Link, join your lab's group as a member. Case management and basic functions will be available after you join the group.

- ✓ Receiving Data for Making Dental Prostheses
- ✓ Sharing work files and production processes.

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resistertest@medit.com

Check

Please ensure your email address is active. Medit Link will send a verification email.

Registering with a disposable and/or fake e-mail address can make you vulnerable to security issues and make it difficult to receive information from Medit Link.

Password

Confirm Password

Confirm Password

Please type in the same password as above.

Select Country

Select Country

☐ I agree to all of the Terms and Conditions.

☐ Personal Information Collection Agreement (Mandatory)

☐ Terms & Conditions (Mandatory)

☐ Agree to the Marketing and Promotion Agreement (optional)

Check this option to receive marketing emails, which include product updates, software news, and promotional events.

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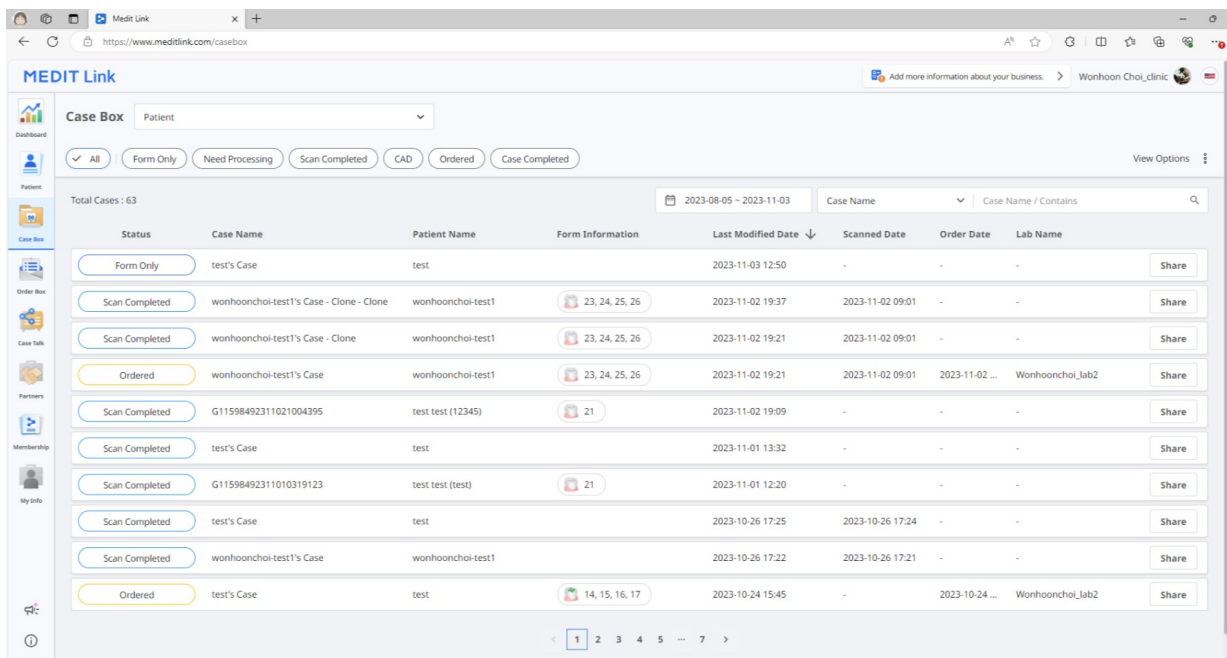
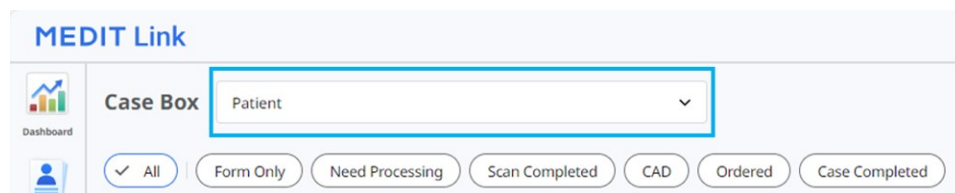
The screenshot shows the Medit Link Case Box interface. At the top, there's a 'Case Box' dropdown menu set to 'Patient'. Below it are filter buttons: 'All', 'Form Only', 'Need Processing', 'Scan Completed', 'CAD', 'Ordered', and 'Case Completed'. A date range filter is set to '2023-08-05 - 2023-11-03'. A search bar for 'Case Name' is also present. The main table lists cases with columns: Status, Case Name, Patient Name, Form Information, Last Modified Date, Scanned Date, Order Date, and Lab Name. Each row has a 'Share' button.

Status	Case Name	Patient Name	Form Information	Last Modified Date	Scanned Date	Order Date	Lab Name
Form Only	test's Case	test		2023-11-03 12:50	-	-	-
Scan Completed	wonhoonchoi-test1's Case - Clone - Clone	wonhoonchoi-test1	23, 24, 25, 26	2023-11-02 19:37	2023-11-02 09:01	-	-
Scan Completed	wonhoonchoi-test1's Case - Clone	wonhoonchoi-test1	23, 24, 25, 26	2023-11-02 19:21	2023-11-02 09:01	-	-
Ordered	wonhoonchoi-test1's Case	wonhoonchoi-test1	23, 24, 25, 26	2023-11-02 19:21	2023-11-02 09:01	2023-11-02 ...	Wonhoonchoi_lab2
Scan Completed	G11598492311021004395	test test (12345)	21	2023-11-02 19:09	-	-	-
Scan Completed	test's Case	test		2023-11-01 13:32	-	-	-
Scan Completed	G11598492311010319123	test test (test)	21	2023-11-01 12:20	-	-	-
Scan Completed	test's Case	test		2023-10-26 17:25	2023-10-26 17:24	-	-
Scan Completed	wonhoonchoi-test1's Case	wonhoonchoi-test1		2023-10-26 17:22	2023-10-26 17:21	-	-
Ordered	test's Case	test	14, 15, 16, 17	2023-10-24 15:45	-	2023-10-24 ...	Wonhoonchoi_lab2

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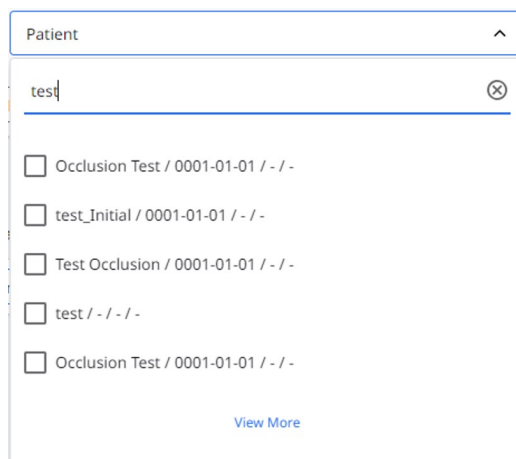
The screenshot shows the Medit Link Case Box interface. The 'Case Box' dropdown menu is open, showing 'Patient' as the selected option. Below the dropdown are filter buttons: 'All', 'Form Only', 'Need Processing', 'Scan Completed', 'CAD', 'Ordered', and 'Case Completed'.

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The screenshot shows the search results dropdown menu. The search field contains 'test'. Below the search field is a list of search results, each with a checkbox and a label: 'Occlusion Test / 0001-01-01 / - / -', 'test_initial / 0001-01-01 / - / -', 'Test Occlusion / 0001-01-01 / - / -', 'test / - / - / -', and 'Occlusion Test / 0001-01-01 / - / -'. A 'View More' link is at the bottom.

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MEDIT Link

Case Box Patient

☒ All
 ☐ Form Only
 ☐ Need Processing
 ☐ Scan Completed
 ☐ CAD
 ☐ Ordered
 ☐ Case Completed

Total Cases: 63

2023-08-05 ~ 2023-11-03 Case Name Case Name / Contains

Status	Case Name	Patient Name	Form Information	Last Modified Date	Scanned Date	Order Date	Lab Name	Share
Form Only	test's Case	test		2023-11-03 12:50	-	-	-	Share
Scan Completed	wonhoonchoi-test1's Case - Clone - Clone	wonhoonchoi-test1	23, 24, 25, 26	2023-11-02 19:37	2023-11-02 09:01	-	-	Share
Scan Completed	wonhoonchoi-test1's Case - Clone	wonhoonchoi-test1	23, 24, 25, 26	2023-11-02 19:21	2023-11-02 09:01	-	-	Share
Ordered	wonhoonchoi-test1's Case	wonhoonchoi-test1	23, 24, 25, 26	2023-11-02 19:21	2023-11-02 09:01	2023-11-02 ...	Wonhoonchoi Lab2	Share
Scan Completed	G11598492311021004395	test test (12345)	21	2023-11-02 19:09	-	-	-	Share
Scan Completed	test's Case	test		2023-11-01 13:32	-	-	-	Share
Scan Completed	G11598492311010319123	test test (test)	21	2023-11-01 12:20	-	-	-	Share
Scan Completed	test's Case	test		2023-10-26 17:25	2023-10-26 17:24	-	-	Share
Scan Completed	wonhoonchoi-test1's Case	wonhoonchoi-test1		2023-10-26 17:22	2023-10-26 17:21	-	-	Share
Ordered	test's Case	test	14, 15, 16, 17	2023-10-24 15:45	-	2023-10-24 ...	Wonhoonchoi Lab2	Share

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☐ All
 ☒ Form Only
 ☐ Need Processing
 ☒ Scan Completed
 ☐ CAD
 ☐ Ordered
 ☐ Case Completed

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Case Box(歯科医院)	All / Form Only / Need Processing / Scan Completed / CAD / Ordered / Case Completed
Work Box(技工所)	All / Form Only / Need Processing / Scan Completed / CAD / Done

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Standard

Free

- ✓ Free plan for every member
- ✓ 1 GB of cloud storage
- ✓ Efficient for new users

If you have any difficulties with payment, select "Click" on the right. [Click](#)

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[All](#) [1 Week Ago](#) [1 Month Ago](#) [6 Month Ago](#) [1 Year Ago](#)

☐ Total Selected 0.00 KB [Delete](#)

☐ Raw Data	147.52 GB
Remove raw data. Raw data does not take up space in cloud storage.	
☐ Result Files & Attachments	13.74 GB
Remove result files and attachments such as mesh, image, or video files.	

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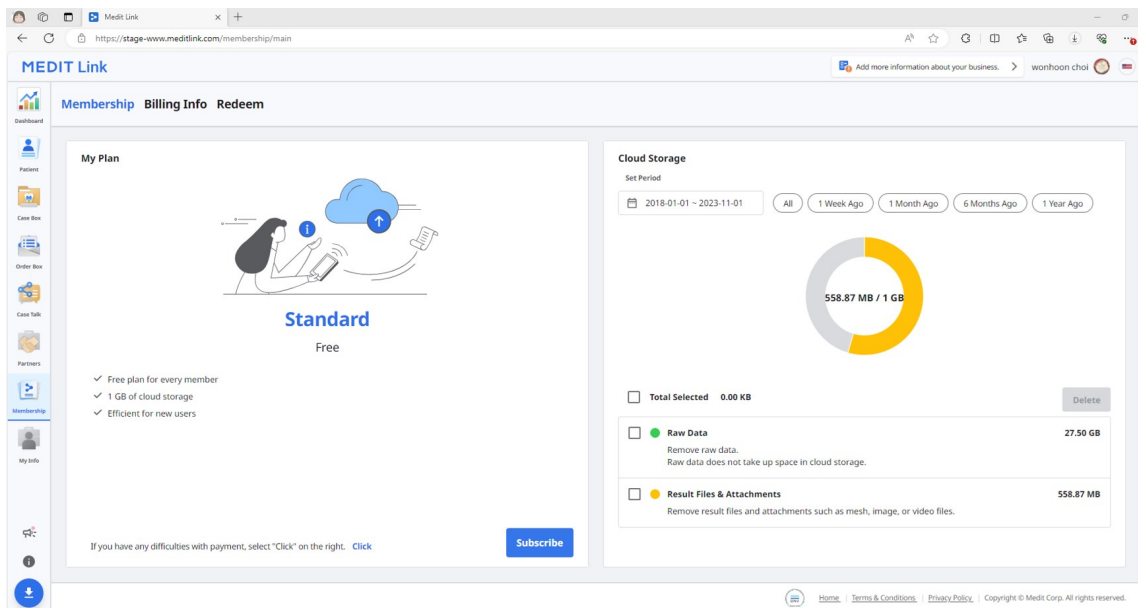
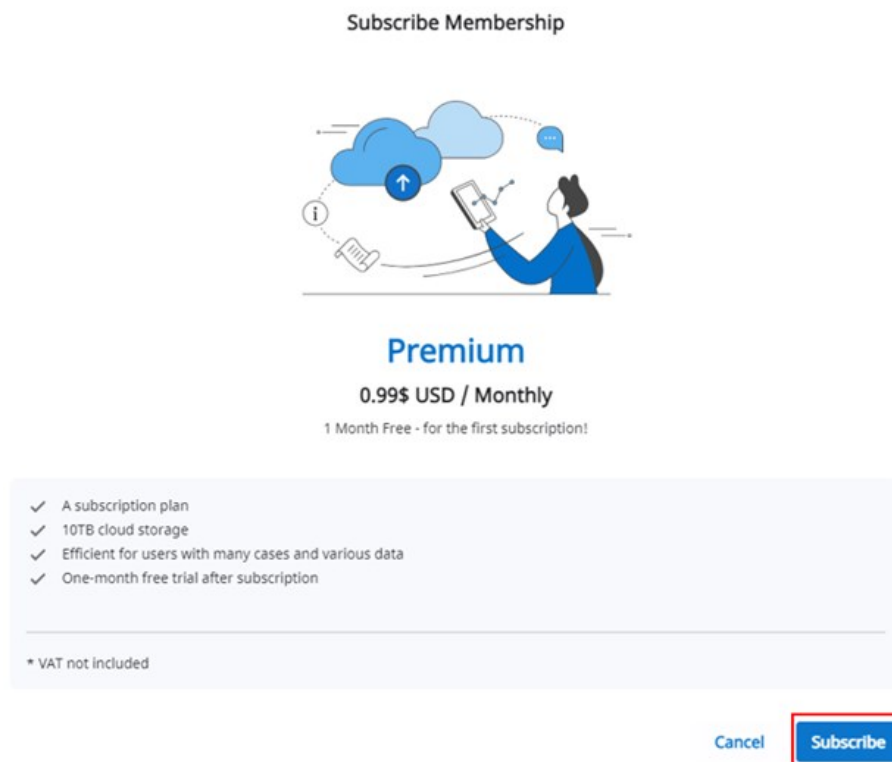
← MEDIT Link TEST MODE

Subscribe to Premium membership for Clinic

\$0.00

Then \$0.99 per day

Medit Link Membership

Premium membership for Clinic	\$0.99
Billed daily	
Subtotal	\$0.99
1 MONTH FREE – for first subscription 100% off for a day	
Total due today	\$0.00

Powered by stripe

Terms Privacy

Pay with card

Email

Card information

1234 1234 1234 1234

VISA

MM / YY

CVC

Name on card

Country or region

South Korea

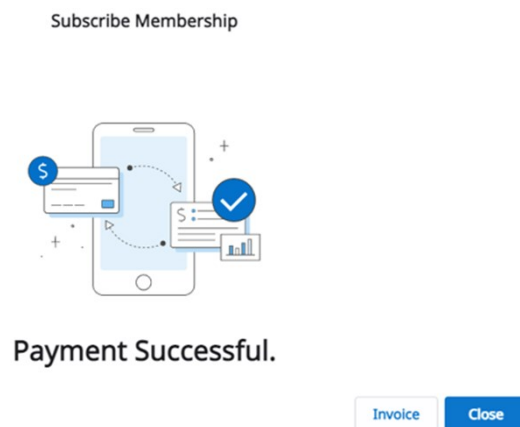
☐

Save information to pay faster next time

Subscribe

By confirming your subscription, you allow MEDIT SG PTE. LTD. to charge your card for this payment and future payments in accordance with their terms.

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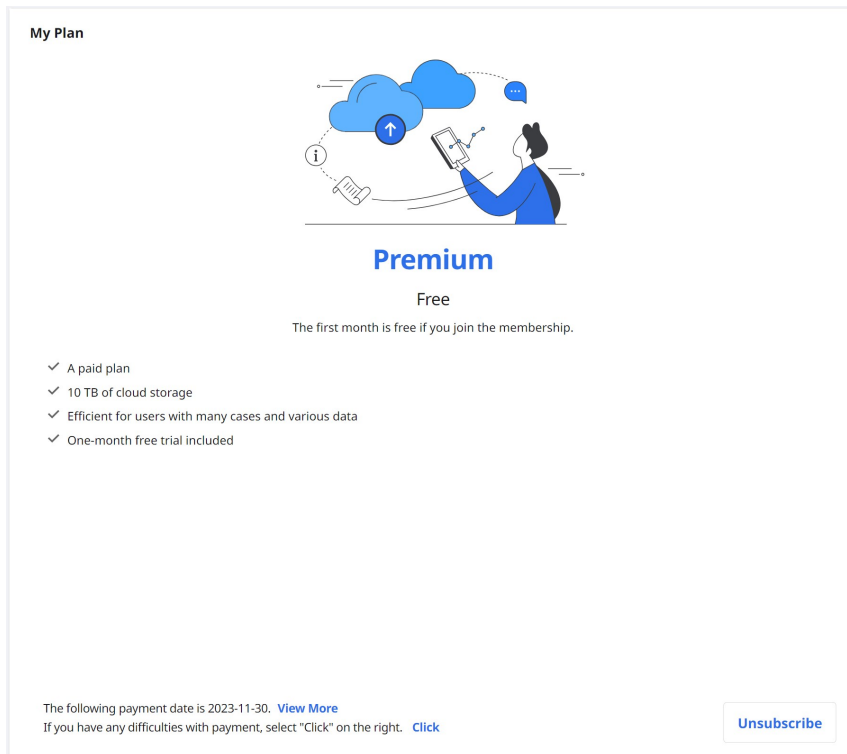
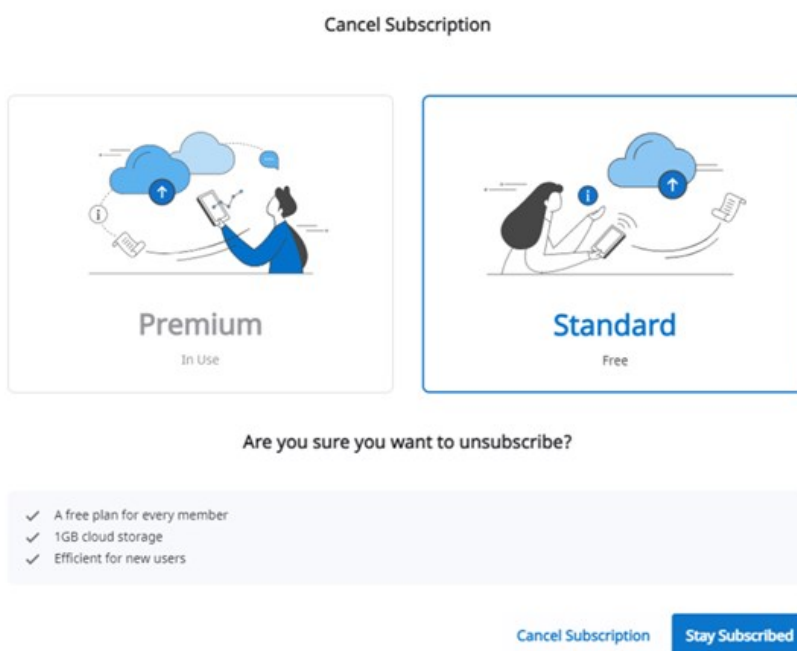
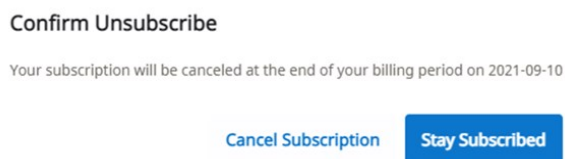
Confirm Unsubscribe

Your subscription has been canceled.
Your current plan can use up to 2021-09-10. After that, the cloud capacity is limited to 1GB.
If the cloud capacity is exceeded, the files will not be backed up, and the order and share functions will not be available.

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Membership **Billing Info** Redeem

My Plan and Billing Information
Plan Premium
Amount 0.99\$ USD / Monthly
Next Payment Date 2021-09-09
Payment Method Visa *****4242 (11/21) | Stripe [Edit Payment Method](#)
If you haven't subscribed yet?
Experience a differentiated premium.

Billing History

Payment Date	Membership	From	To	Payment Amount	Status
2021-09-09	Premium	2021-09-09	2021-09-10	0.99 \$ (USD)	Payment Completed
2021-09-09	Premium	2021-09-09	2021-09-09	0 \$ (USD)	Payment Completed

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Add

Go back

By adding a payment method, you agree to MEDIT SG PTE. LTD.'s Terms of Service and Privacy Policy.

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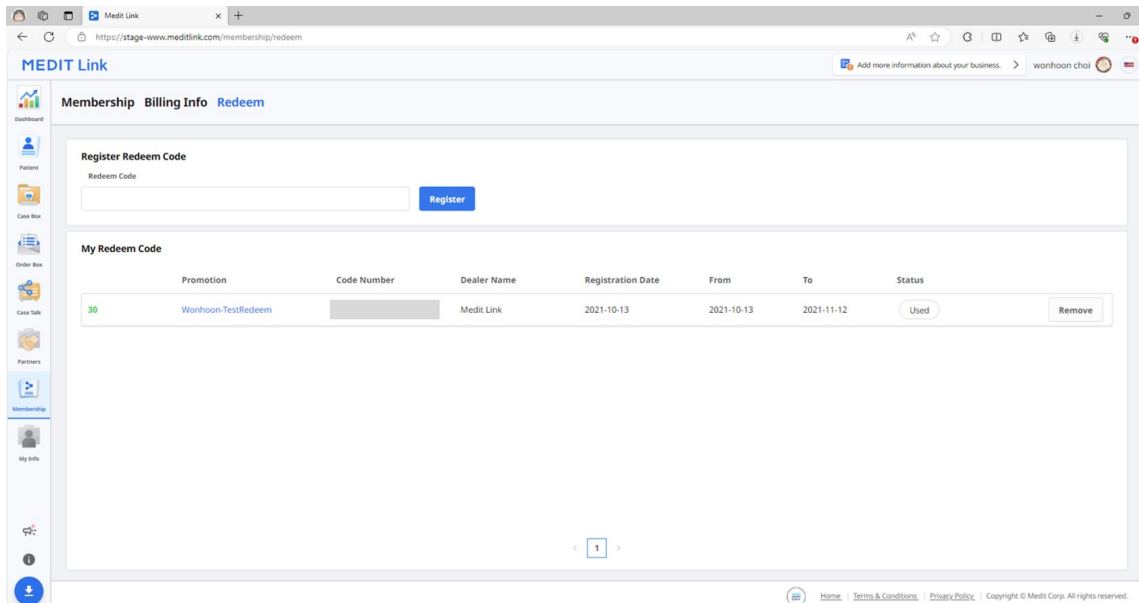
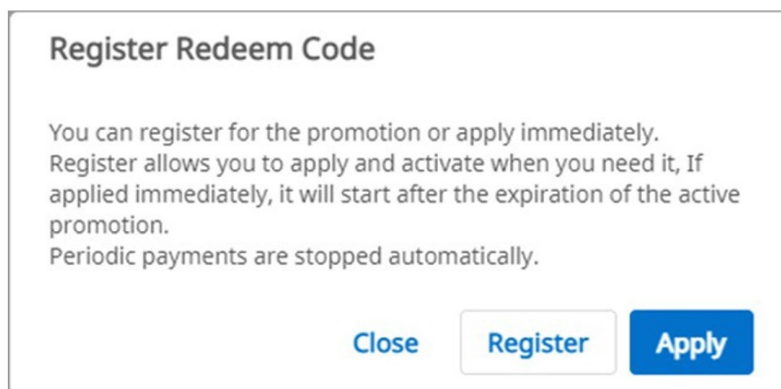
Country
Korea, Republic of ▼

Your chosen language is applied to the email sent by Medit Link and the notification.

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Edit

Cover and Logo

Looking for digitalized dental laboratories?

Find the most ordered Labs around your clinic with 'Featured Part

Services

Select

Specialties (choose up to 3)

Select

About Us

e.g., Welcome to OOO Dental Clinic.
We are a licensed dental clinic situated in California specializing mainly in aesthetic restorations and implantology.

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-

Buy

ⓘ The scanner information is automatically registered once your scanner is connected to a Medit scan program.

Home

Terms & Conditions

Privacy Policy

Copyright © Medit Corp. All rights reserved.

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Notice
NFT Service Termination Notice
2023-04-19

To our users,

We want to thank all Medit Link users who have supported and used our NFT service.

Unfortunately, we must inform you that this service will be terminated following the upcoming program update.

Please see the details below for more information:

- **Feature in Question**
 - Share via NFT (in Medit Link Case Talk)
- **Effective Date**
 - Starting from the update in the first half of 2023 (scheduled for April 27)

※ The update schedule is subject to change.

We'll continuously work to provide better service.

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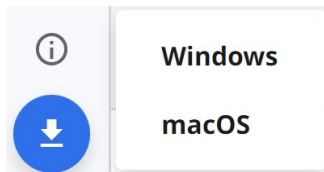
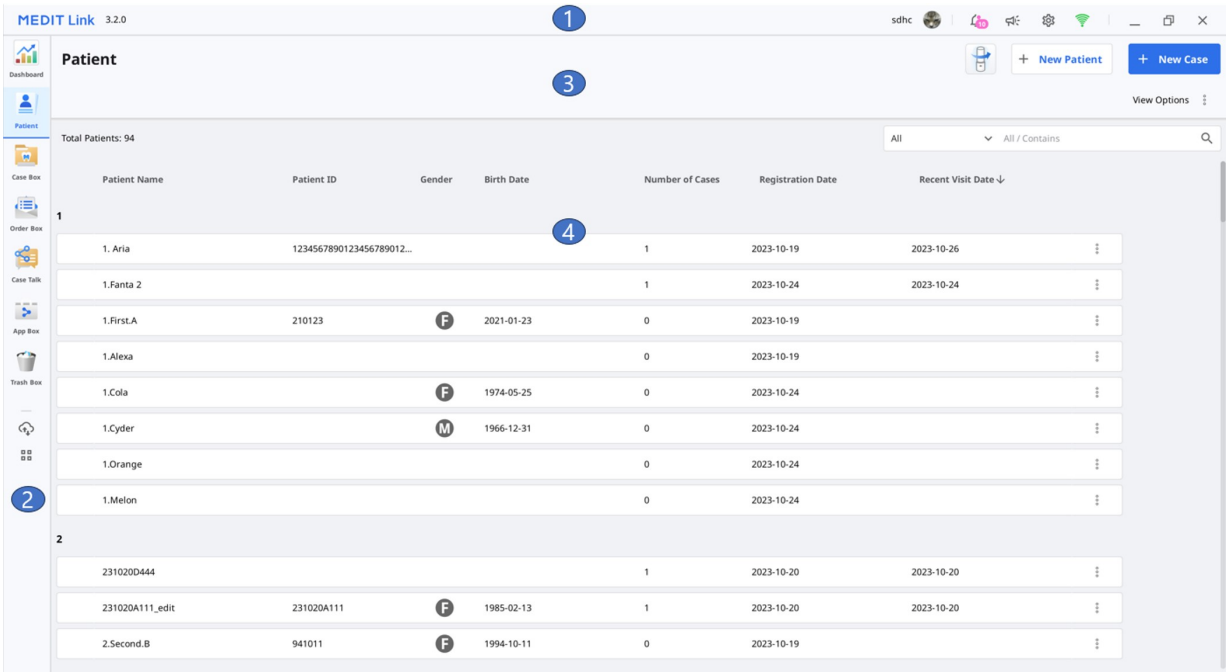
2

Delivery 2020-08-20 오후 5:20

Wonhoonchoi_Lab_Dev has dispatched the ordered package.

Order 2020-08-20 오후 5:18

Wonhoonchoi_Lab_Dev has accepted the order.

Order 2020-08-17 오후 4:05

Wonhoonchoi_Lab_Dev has accepted the order.

Order 2020-08-17 오후 4:05

Wonhoonchoi_Lab_Dev has accepted the order.

Order 2020-08-17 오후 4:05

Wonhoonchoi_Lab_Dev has accepted the order.

Order 2020-08-17 오후 4:05

Wonhoonchoi_Lab_Dev has accepted the order.

Work 2020-08-14 오후 6:33

Wonhoonchoi_Clinice_Dev_Member registered a new case test's Case.

Work 2020-08-14 오전 10:15

Wonhoonchoi_Clinice_Dev_Member registered a new case test's Case.

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Network Status

-  Network speed is good.
-  Network speed is fast enough.
-  Network speed is poor. Can be switched to offline mode automatically.
-  Either network is not connected or Medit Link is in offline mode.
-  Synchronization is in progress. Some case lists may not match.
-  Synchronization failed. The case lists may not match.

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The screenshot shows the 'Settings' window with the 'General' tab selected. The left sidebar contains links for 'General', 'Storage Management', 'Set Connected Program', 'Export Options', 'Case Converting Tool', and 'Proxy Option'. The main area displays the following settings:

- Save Folder:** C:/ProgramData/Medit/Medit Link/Data/
- Language Settings:** English
- Teeth Labeling:** FDI (selected), US
- Notification:** Allow Notifications (checked)
- Upload Raw Data:** Upload Automatically (checked)
- Quick Scan:** Choose if you would like to have the scan stages below automatically appear in Medit Scan for Clinics without having to fill in the form information. Pre-Op for Maxilla, Pre-Op for Mandible, Face Scan (unchecked)
- Viewing Options:** Advanced Rendering (checked), Adjust Color Texture (slider)
- Medit User Experience Improvement Program:** We value your feedback! Join the Medit User Experience Improvement Program and share your thoughts. The collected information will be used to improve the product and services. Medit Link does not collect any information that identifies individuals.

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Compress & Combine

☒ Compress Files
☒ Combine Individual Mesh

Options when same folder or file already exists.

☒ Always ask when exporting
☐ Export a new file and overwrite the existing file. The file currently stored will be lost.
☐ Create a folder and export the files

Rule for Folder Name

☒ Case Creation Date
☐ Scanned Date
☐ Order Date

☒ Patient Name
☐ Case Name
☐ Partner Name

☐ Patient ID
☐ Order ID

Folder Name

Case Creation Date - Patient Name

Rule for File Name

☒ Case Creation Date
☐ Scanned Date
☐ Order Date

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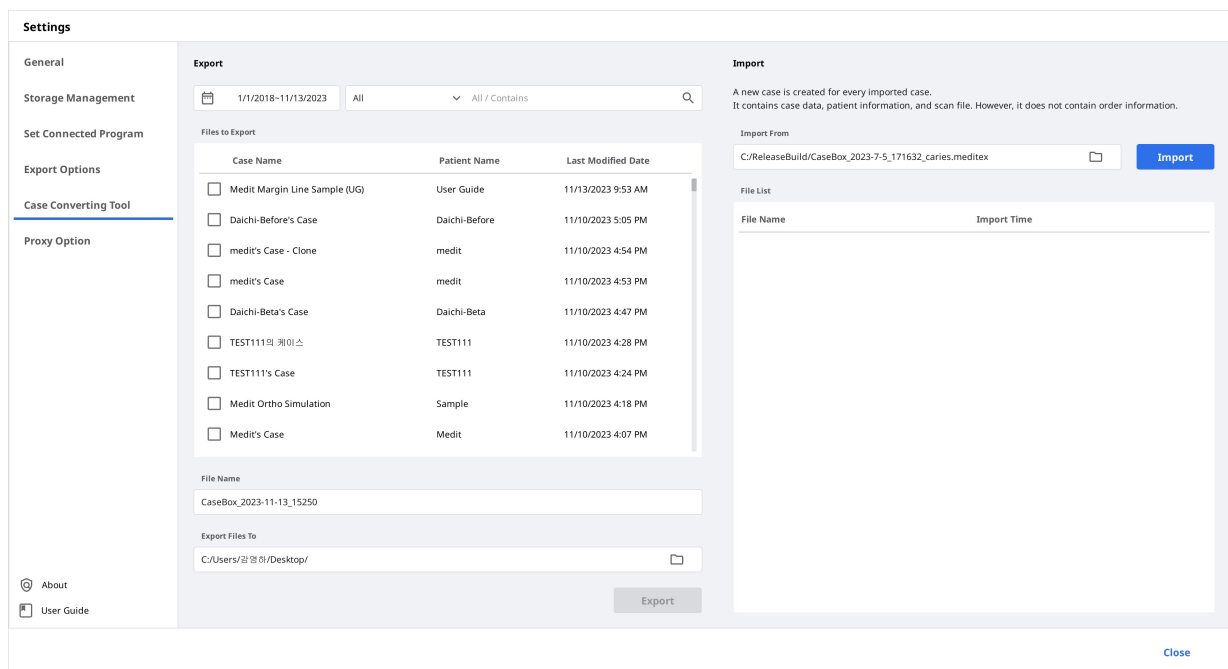
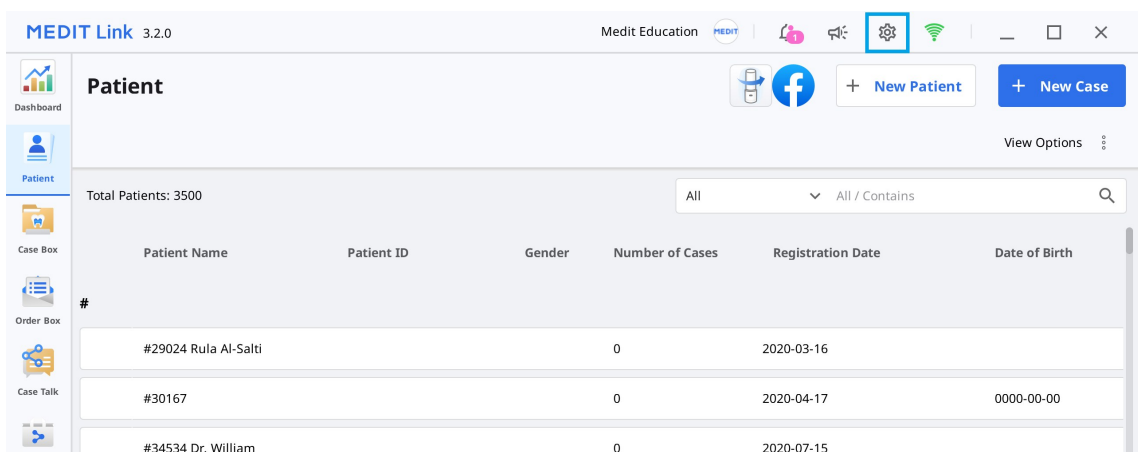
Export

Import

A new case is created for every imported case. It contains case data, patient information, and scan file. However, it does not contain order information.

Import From
C:/ReleaseBuild/CaseBox_2023-7-5_171632_caries.meditex

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1/1/2018~11/13/2023 All / Contains

Files to Export

Case Name	Patient Name	Last Modified Date
☒ Medit Margin Line Sample (UG)	User Guide	11/13/2023 9:53 AM
☒ Daichi-Before's Case	Daichi-Before	11/10/2023 5:05 PM

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Question

Exporting successfully completed.
Do you want to open the folder containing the exported file?

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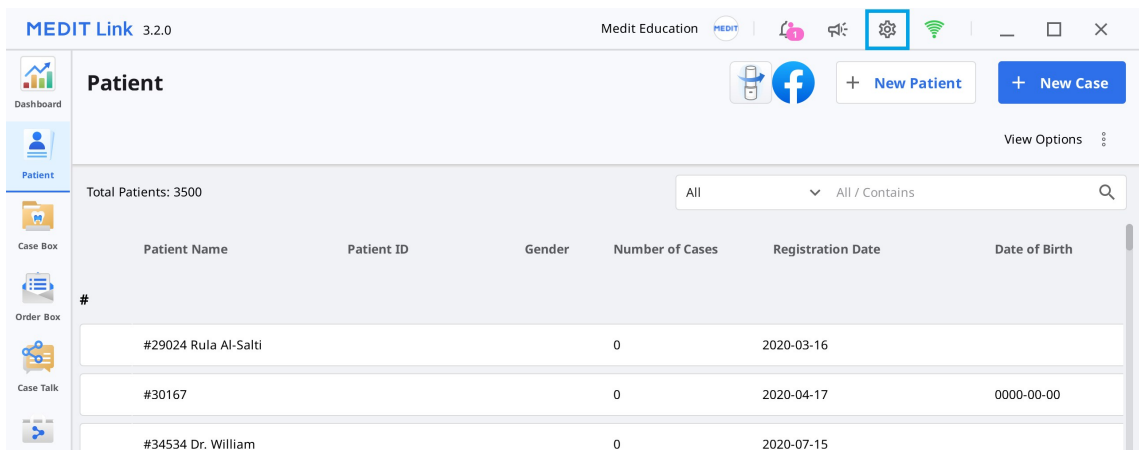
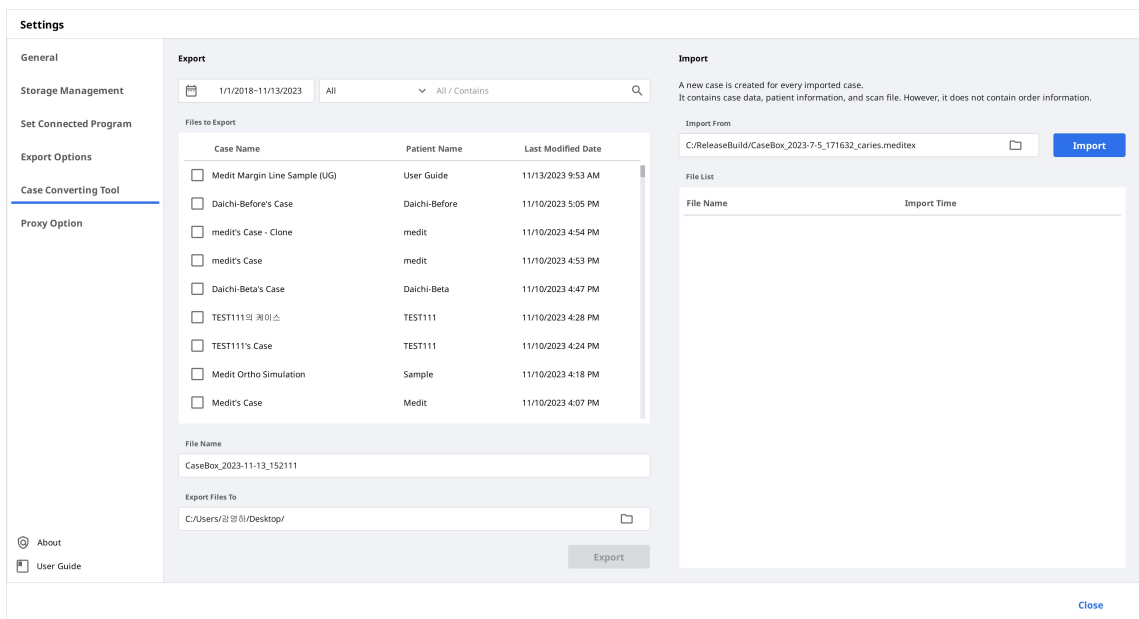
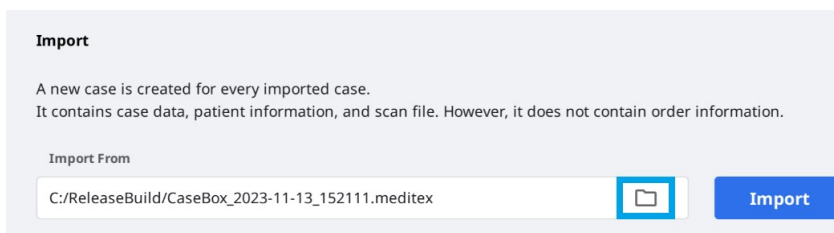
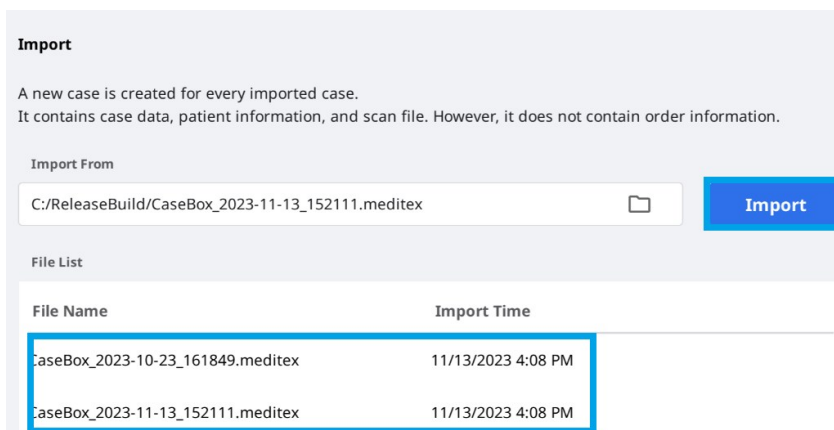
Operation Korea

<https://www.medit.com/>

Copyright Medit Corp. All rights reserved.

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Merge all patient information to

Please select the patient you'd like to merge all data to.
After merging all data together, all information in other patient files will be deleted.

	Patient Name	Patient ID	Gender	Birth Date	Number of Cases	Date Registered	Recent Visit Date
☒	5-9_jScan_CAMOrd...				1	2020-09-18	2020-09-18
☐	test				2	2020-09-17	2020-09-17
☐	test				22	2019-07-26	2020-09-17

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Case Box/Work Box

1 Case name and patient name display.

2 Display of main functions: Order/Scan (Open)/Order button and app icon.

3 Display of available tabs: Form, CAD, File Viewer.

4 Display of case and file management functions.

5 Display of the Form registration UI for CAD program selection.

6 Display of Type/Method/Shade information for the selected tooth.

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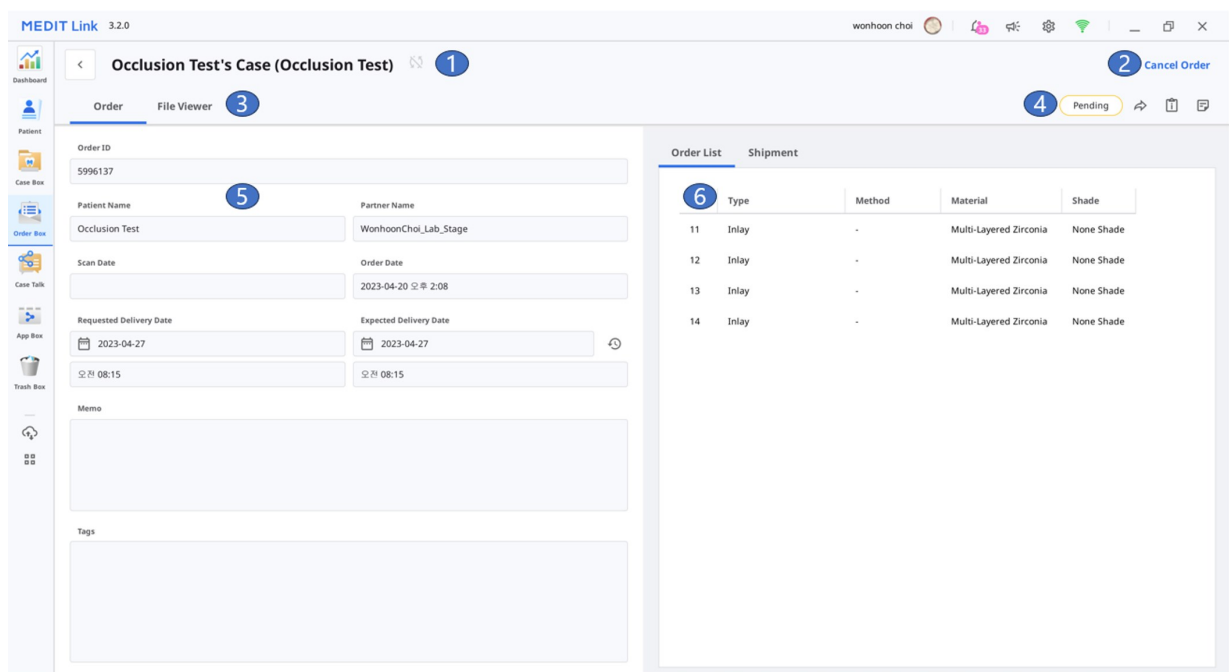
The screenshot shows the Medit Link Case Talk interface. At the top, there's a header with 'MEDIT Link 3.2.0' and a user profile 'wonhoon choi'. Below the header, there's a sidebar with icons for Dashboard, Patient, Case Box, Order Box, Case Talk (selected), App Box, and Trash Box. The main area displays a table of shared cases. The table has columns: Case Talk Name, Patient Name, Form Information, Scanned Date, Shared Date, and Expiration Date. There are three rows of data.

Case Talk Name	Patient Name	Form Information	Scanned Date	Shared Date	Expiration Date
BS_4_Scan for Clinics's Case - Clone_2023-10-31 15:04:08	BS***** f* C*****	Maxilla, Mandible	2020-09-25 14:12	2023-10-31 15:04	2024-01-29 23:59
BS_4_Scan for Clinics's Case - Clone_2023-10-20 11:59:15	BS_4_Scan for Clinics	Maxilla, Mandible	2020-09-25 14:12	2023-10-20 11:59	2023-11-20 23:59
1111	Email_Push	24, 25	2022-09-30 13:48	2022-09-30 16:35	unlimited

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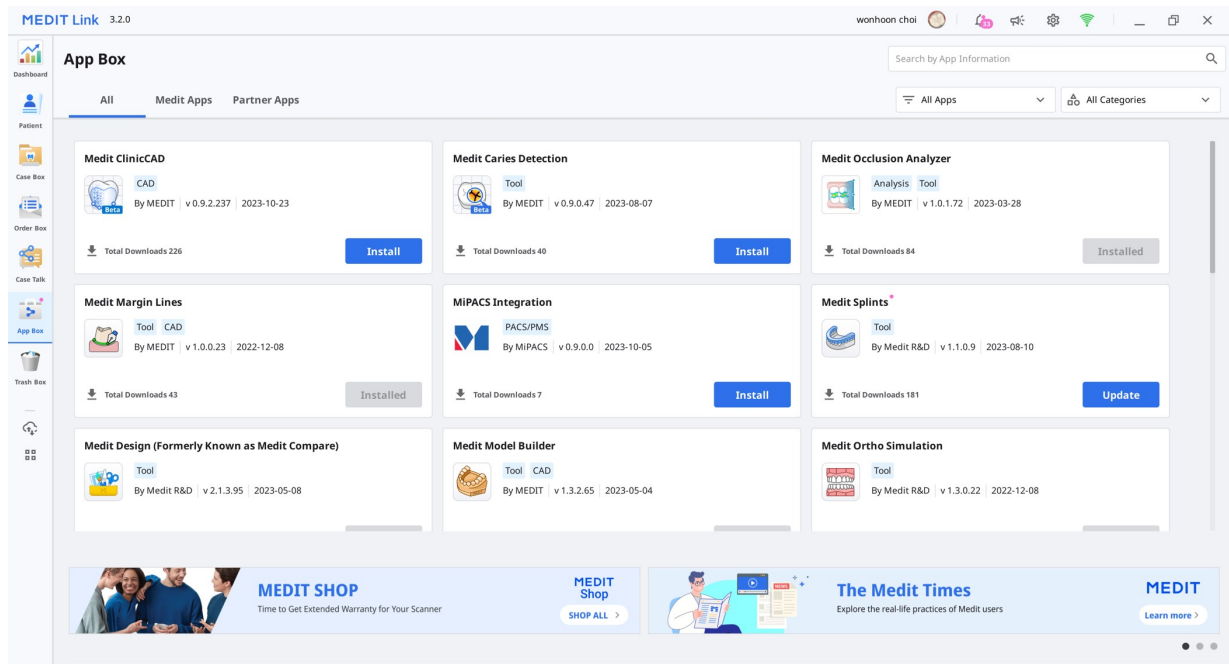
Restore | Selected Cases: 0

All / Contains

Status	Case Name	Patient Name	Lab Name	Deletion Date ↓	Deleted By	Size
☐ Scan Completed	Jira_임플란트 46, 47's Case - Clone - Clone	임플란트 46, 47		2023-04-26 오전 9:12	wonhoon choi	66,654 KB
☐ Form Only	G11598492304180810429	test test		2023-04-18 오후 5:15	wonhoon choi	0 KB
☐ Form Only	G11598492304180806313	test test		2023-04-18 오후 5:15	wonhoon choi	0 KB
☐ Form Only	G11598492304180811349	tes test		2023-04-18 오후 5:15	wonhoon choi	0 KB
☐ Form Only	G11598492304180754083	Glidewell_fastdesign		2023-04-18 오후 5:00	wonhoon choi	0 KB
☐ Form Only	G11598492304180733332	Glidewell_fastdesign		2023-04-18 오후 5:00	wonhoon choi	0 KB
☐ Form Only	G11598492304180648146	test test		2023-04-18 오후 5:00	wonhoon choi	0 KB
☐ Case Completed	Multi Occlusion Normal Case's Case - Clone	Multi Occlusion Normal Case		2023-04-18 오전 10:50	wonhoon choi	23,652 KB
☐ CAD	Multi Occlusion Normal Case's Case - Clone - Clone	Multi Occlusion Normal Case		2023-04-18 오전 10:50	wonhoon choi	23,652 KB
☐ Form Only	Occlusion Test's Case	Occlusion Test		2023-04-17 오후 6:31	wonhoon choi	0 KB
☐ Form Only	Test Occlusion 's Case	Test Occlusion		2023-04-17 오후 6:31	wonhoon choi	0 KB
☐ Form Only	Occlusion Test's Case	Occlusion Test		2023-04-17 오후 6:31	wonhoon choi	0 KB
☐ Form Only	Occlusion Test's Case	Occlusion Test		2023-04-17 오후 6:26	wonhoon choi	0 KB

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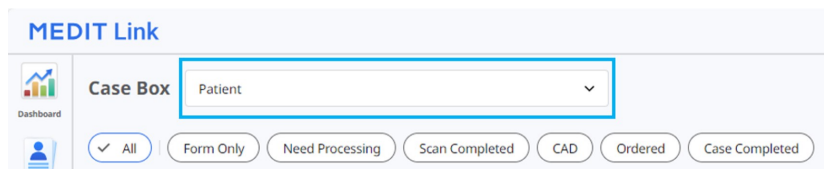
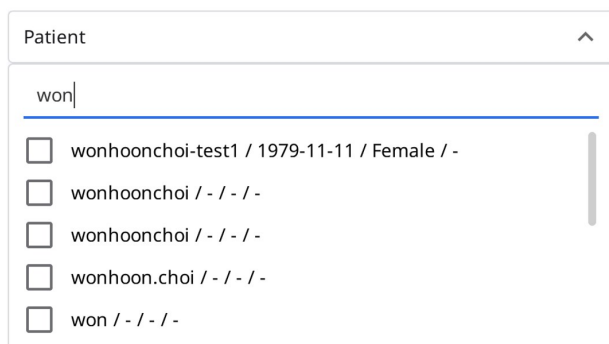
MEDIT Link

Case Box: Patient

Filters: All (selected), Form Only, Need Processing, Scan Completed, CAD, Ordered, Case Completed

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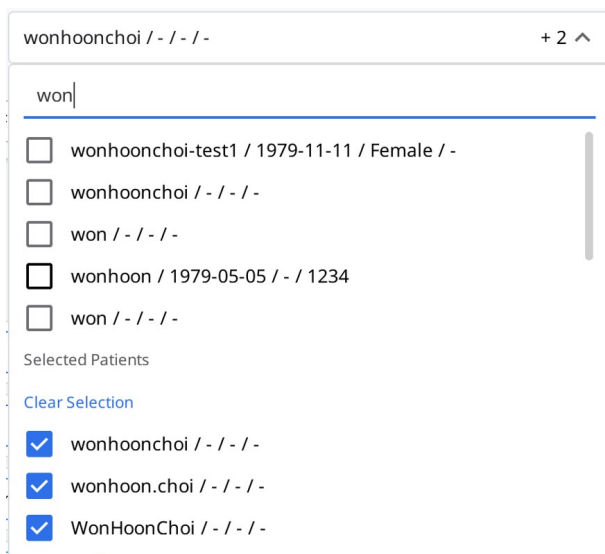
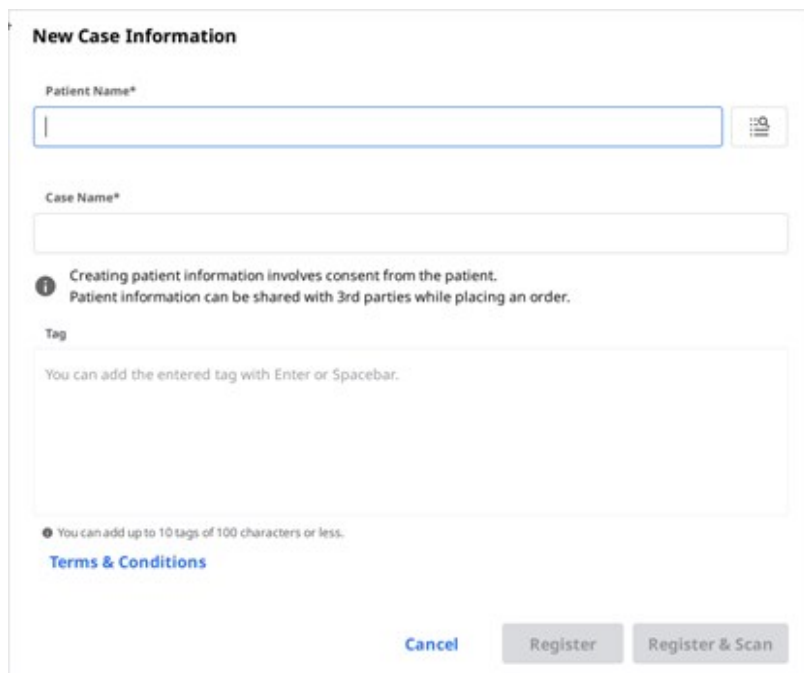
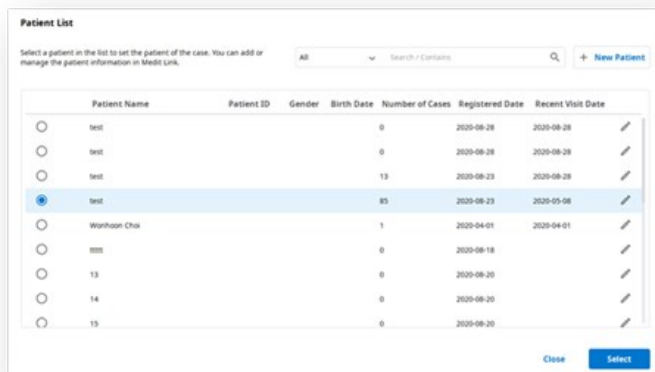
The form titled "New Case Information" contains the following fields and elements:

- Patient Name***: A text input field with a dropdown icon on the right.
- Case Name***: A text input field.
- Information**: A note stating "Creating patient information involves consent from the patient. Patient information can be shared with 3rd parties while placing an order."
- Tag**: A text area with the instruction "You can add the entered tag with Enter or Spacebar."
- Footer**: Includes a note "You can add up to 10 tags of 100 characters or less.", a link for "Terms & Conditions", and three buttons: "Cancel", "Register", and "Register & Scan".

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The "Patient List" dialog shows a table of existing patients. The table has columns for Patient Name, Patient ID, Gender, Birth Date, Number of Cases, Registered Date, and Recent Visit Date. A search bar and a "New Patient" button are at the top. A "Select" button is at the bottom right.

	Patient Name	Patient ID	Gender	Birth Date	Number of Cases	Registered Date	Recent Visit Date	
☐	test				0	2020-08-28	2020-08-28	
☐	test				0	2020-08-28	2020-08-28	
☐	test				13	2020-08-23	2020-08-28	
☒	test				85	2020-08-23	2020-05-08	
☐	Wanhuan Choi				1	2020-04-01	2020-04-01	
☐	test				0	2020-08-18		
☐	13				0	2020-08-20		
☐	14				0	2020-08-20		
☐	15				0	2020-08-20		

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The screenshot shows a form titled "New Case Information". It contains the following fields and elements:

- Clinic Name***: A dropdown menu.
- Patient Name***: A text input field.
- Case Name***: A text input field.
- Requested Delivery Date**: A date picker showing "2023-04-13".
- Tag**: A text area with the instruction "You can add the entered tag with Enter or Spacebar." and a note "You can add up to 10 tags of 100 characters or less."
- At the bottom right, there are two buttons: "Cancel" and "Register".

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The screenshot shows a form titled "Edit Case Information". It contains the following fields and elements:

- Patient Name***: A text input field containing "Test Occlusion". To its right is a search icon.
- A green checkmark and text: "Test Occlusion 00010101" is selected.
- Case Name***: A text input field containing "mandibular movement - Clone - Clone".
- An information icon and text: "Creating patient information involves consent from the patient. Patient information can be shared with 3rd parties while placing an order."
- Tag**: A text area with the instruction "You can add the entered tag with Enter or Spacebar." and a note "You can add up to 10 tags of 100 characters or less."
- A link: "Terms & Conditions".
- At the bottom right, there are two buttons: "Cancel" and "Apply".

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Patient List

Select a patient in the list to set the patient of the case. You can add or manage the patient information in Medit Link.

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2022-04-18

Tag

You can add the entered tag with Enter or Spacebar.

● You can add up to 10 tags of 100 characters or less.

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C:/Users/Medit-Wonhoonchoi/Desktop/test/

Folder Name

Multi Occlusion - Scan - Before Scanning

File Name

Multi Occlusion - Scan - Before Scanning

CAD Type

Maestro

Mesh Format

Convert mesh files to the below formats

☐ meditMesh

☐ OBJ

☒ PLY

☐ STL

☒ Do not convert the attached files

Image Format

☐ Convert images to the below format

☐ JPG

☐ JPEG

☒ PNG

☐ BMP

☐ Compress Files

☐ Combine Individual Mesh

☐ Export All Occlusal Relationship

☐ Save Case Information as PDF [Preview](#)

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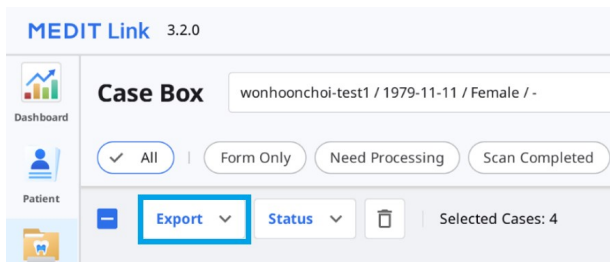
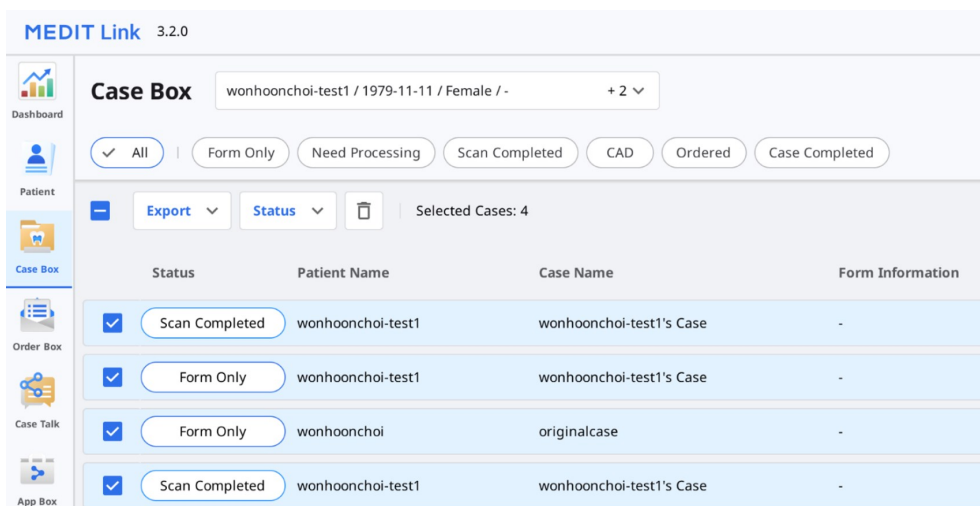
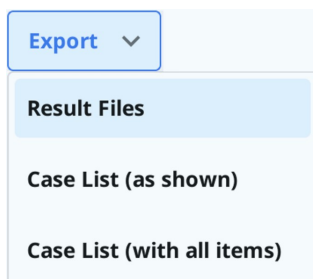
Options

Export Files To
C:/Users/Wonhoon Choi Laptop/Desktop/ExportTest/3/

Folder Name
Case Creation Date - Case Name

File Name
Case Creation Date - Case Name

Type
Medit

Mesh Format
Convert mesh files to the below formats
☐ meditMesh ☐ OBJ ☒ PLY ☐ STL
☐ Do not convert the attached files

Image Format
Convert images to the below format
☐ JPG ☐ JPEG ☒ PNG ☐ BMP

☐ Compress Files
☐ Combine Individual Mesh
☐ Export All Occlusal Relationship
☐ Save Case Information as PDF

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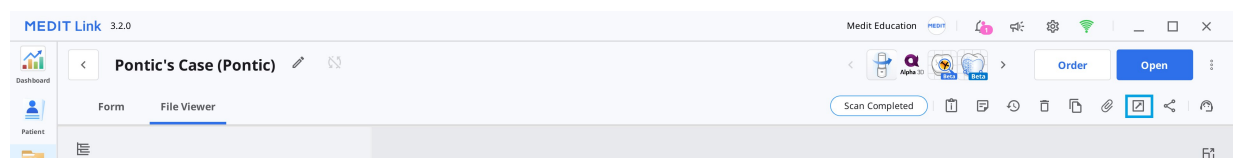
Options File List

Export Files To
C:/Users/Medit-Wonhoonchoi/Desktop/test/

Folder Name
Multi Occlusion - Scan - Before Scanning

File Name
Multi Occlusion - Scan - Before Scanning

CAD Type
Maestro

Mesh Format
Convert mesh files to the below formats
☐ meditMesh ☐ OBJ ☒ PLY ☐ STL
☒ Do not convert the attached files

Image Format
Convert images to the below format
☐ JPG ☐ JPEG ☒ PNG ☐ BMP

☐ Compress Files
☐ Combine Individual Mesh
☐ Export All Occlusal Relationship
☐ Save Case Information as PDF [Preview](#)

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Export Files in This Case

Options File List

☒ [Icon]
☒ Maxilla
☒ Base
☒ Mandible
☒ Base
☒ Occlusion
☒ First Occlusion
☒ Second Occlusion
☒ Dynamic Occlusion
☒ Free
☒ Maxilla
☒ Mandible
☒ Face Scan Data
☒ Face Scan Data

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Question

 Exporting successfully completed.
 Do you want to open the folder containing the exported file?

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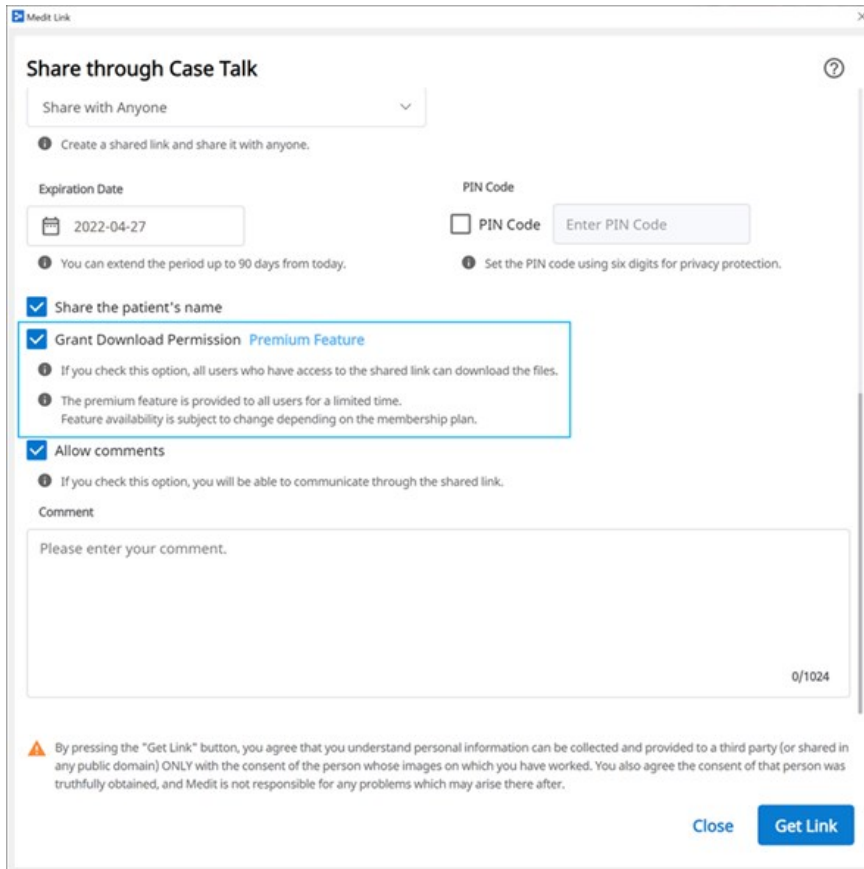
Share through Case Talk

Share with Anyone

Create a shared link and share it with anyone.

Expiration Date: 2022-04-27

PIN Code: ☐ PIN Code Enter PIN Code

You can extend the period up to 90 days from today.

Set the PIN code using six digits for privacy protection.

☒ Share the patient's name

☒ Grant Download Permission Premium Feature

If you check this option, all users who have access to the shared link can download the files.

The premium feature is provided to all users for a limited time. Feature availability is subject to change depending on the membership plan.

☒ Allow comments

If you check this option, you will be able to communicate through the shared link.

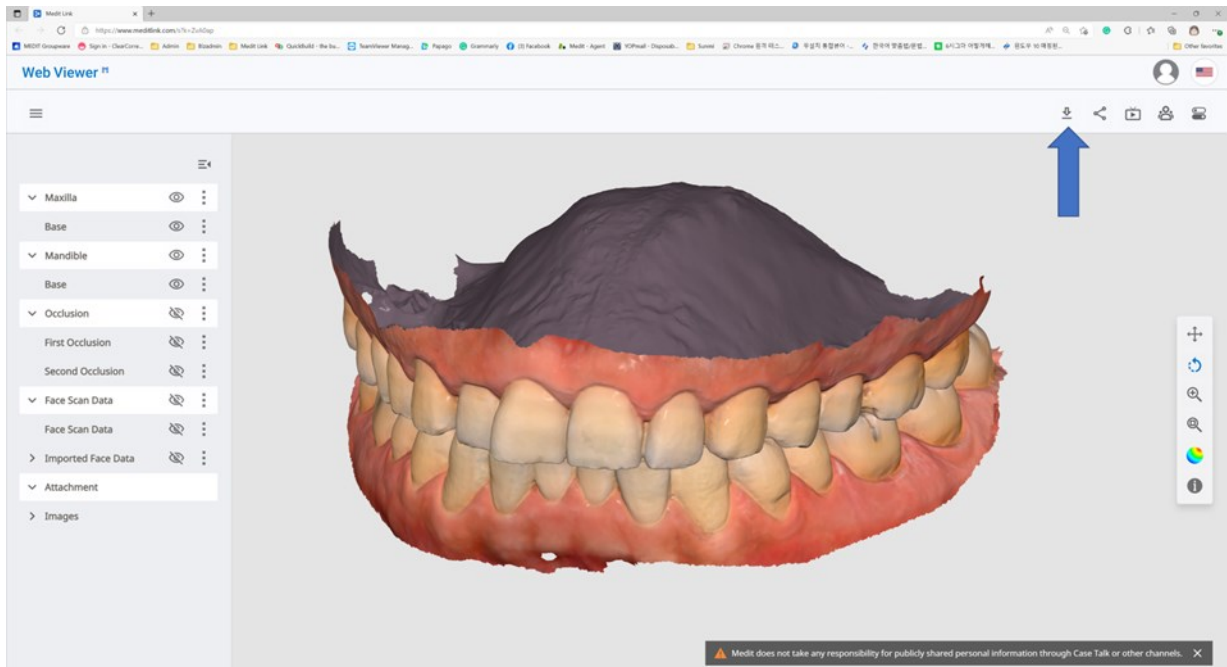
Comment: Please enter your comment.

0/1024

By pressing the "Get Link" button, you agree that you understand personal information can be collected and provided to a third party (or shared in any public domain) ONLY with the consent of the person whose images on which you have worked. You also agree the consent of that person was truthfully obtained, and Medit is not responsible for any problems which may arise there after.

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Download Premium Feature

The premium feature is provided to all users for a limited time.
Feature availability is subject to change depending on the membership plan.

Mesh Format

Convert mesh files to the below formats

☐ meditMesh ☐ OBJ ☒ PLY ☐ STL

i The selected format will be applied to the result files created in the Medit Link application.

i The attachment files will not be converted but will keep their original format.

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Memo

Tags
You can add the entered tag with Enter or Spacebar.

☒ Include all occlusion relationships in the order
☒ Include Tags In The Order

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The screenshot shows the 'Case Box' interface. On the left is a sidebar with icons for Dashboard, Patient, Case Box (highlighted), Order Box, Case Talk, and App Box. The main area is titled 'Case Box' and has a 'Patient' dropdown menu. Below this are filter buttons: 'All' (selected), 'Form Only', 'Need Processing', 'Scan Completed', 'CAD', 'Ordered', and 'Cas'. There are also buttons for 'Export', 'Status' (with a dropdown arrow), and a trash icon. A 'Selected Cases: 4' indicator is present. The 'Status' dropdown menu is open, showing 'Case Completed' as the selected option. Below the menu is a table with columns 'Status', 'Case Name', and 'Case Name'. The table contains four rows, each with a checkbox, a status label, a case ID, and a case name.

	Status	Case Name	Case Name
☐	Ordered	test	test's Case
☒	Scan Completed	wonhoonchoi-test1	wonhoonchoi-test1's Case
☒	Scan Completed	test test(1234)	G11598492310180503424
☒	Scan Completed	test	test's Case
☒	Scan Completed	test	test's Case

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The screenshot shows the 'New Case Information' form. It includes fields for 'Patient Name*' and 'Case Name*'. Below these is an information icon and text: 'Creating patient information involves consent from the patient. Patient information can be shared with 3rd parties while placing an order.' The 'Tag' field is highlighted with a blue border and contains three existing tags: 'Doctorname', 'Specialmaterial', and 'SpecialComment', each with a close button. Below the tag field is a note: 'You can add up to 10 tags of 100 characters or less.' and a link for 'Terms & Conditions'. At the bottom are buttons for 'Cancel', 'Register', and 'Register & Scan'.

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This screenshot shows the same 'New Case Information' form, but the 'Tag' field now contains the letter 't'. Below the field, a search bar shows 't (new)' and a list of suggestions is displayed: 'tet', 'tagtest', 'tag', and 'tat2'. The suggestions list is highlighted with a blue border.

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Case Name
mandibular movement - Clone - Clone

Patient Name
Test Occlusion

☒ Share patient information

Requested Delivery Date
2023-04-13

Delivery Time
오전 8:15 AM 8:00 PM 8:00

Memo

Tags
Doctorname × Specialmaterial × SpecialComment ×

You can add up to 10 tags of 100 characters or less.

☒ Include all occlusion relationships in the order

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