USER GUIDE

Medit Link

Table of Contents

Medit Link

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Lab

Administrator
Dentist (Manager)

All Medit Link functions are available.

- ✓ Order
- ✓ Member Management

Register

Administrator
Representative/Leader (Manager)

All Medit Link functions are available.

- ✓ Order Management
- ✓ Member Management

Register

Member
Dental Hygienist and Dental Staff

In order to use Medit Link, join your group as a member. Case management and basic functions will be available after you join the group.

- ✓ Oral Scanning

Register

Member
Laboratory Employees

In order to use Medit Link, join your lab's group as a member. Case management and basic functions will be available after you join the group.

- ✓ Receiving Data for Making Dental Prostheses
- ✓ Sharing work files and production processes.

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Email Address

Check

ⓘ Please ensure your email address is active. Medit Link will send a verification email.
 ⓘ Registering with a disposable and/or fake e-mail address can make you vulnerable to security issues and make it difficult to receive information from Medit Link.

Password

Confirm Password

ⓘ Please type in the same password as above.

Select Country

Select Country ▼

☐ I agree to all of the Terms and Conditions.

☐ Personal Information Collection Agreement (Mandatory)

View

☐ Terms & Conditions (Mandatory)

View

☐ Agree to the Marketing and Promotion Agreement (optional)

View

ⓘ Check this option to receive marketing emails, which include product updates, software news, and promotional events.

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
☒ Medit Link account

☐ Another email address

The report will be sent on the 2nd day of every month.

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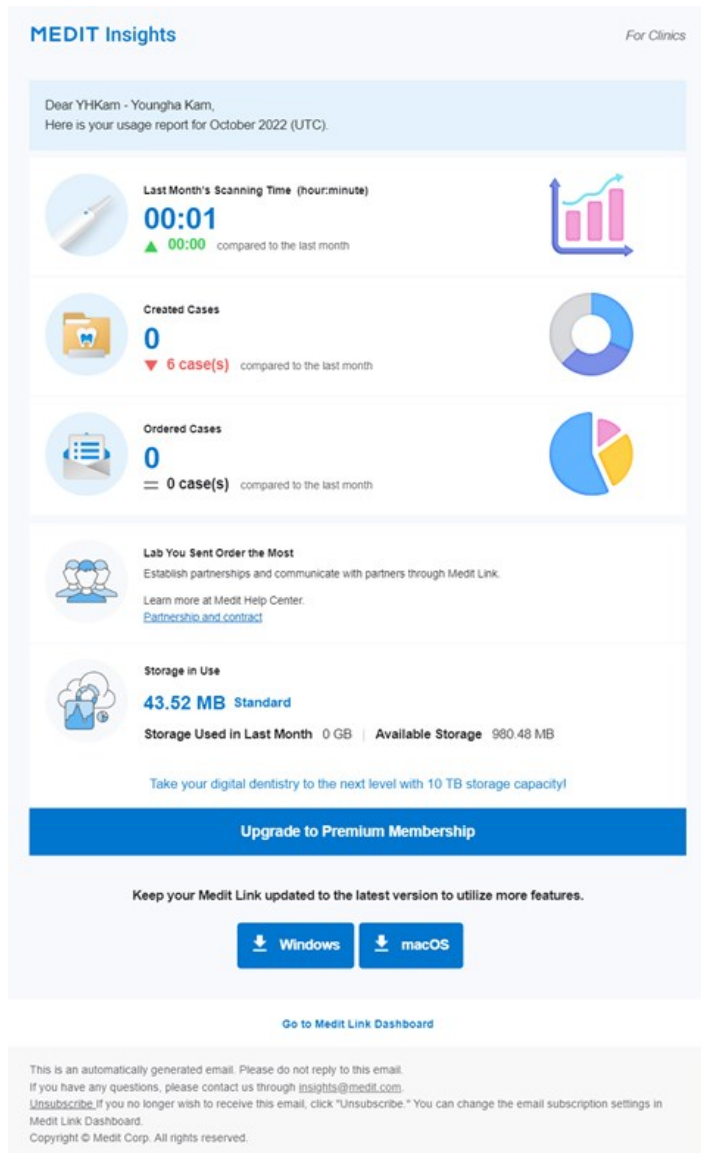
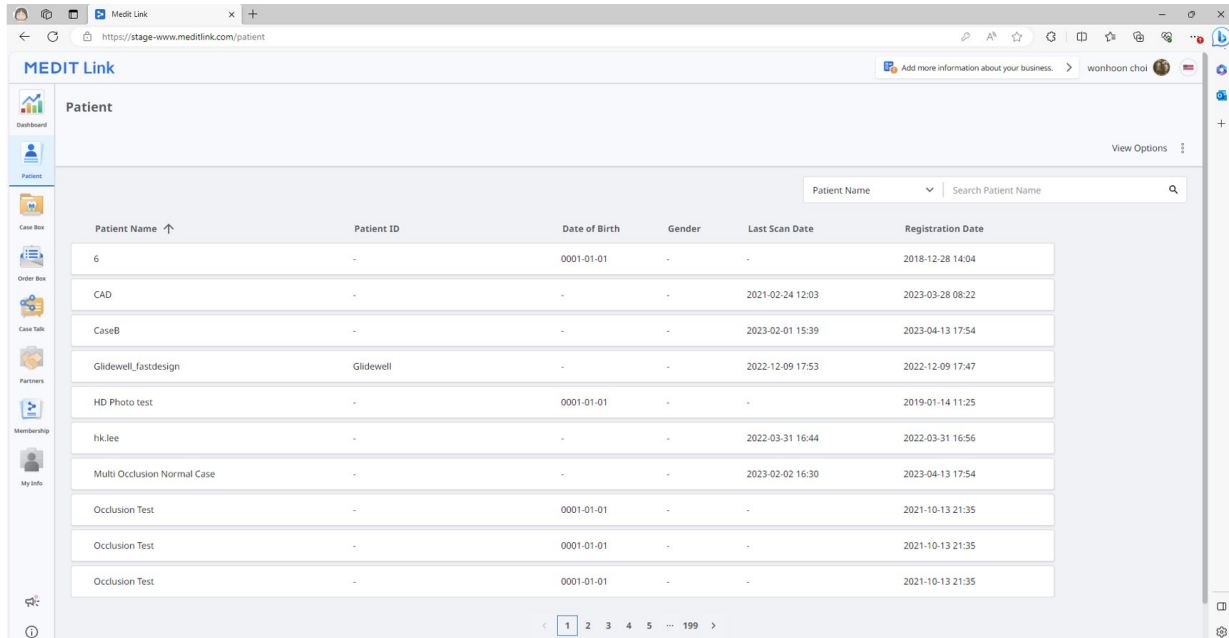
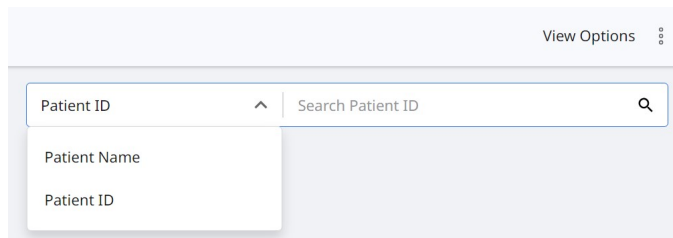
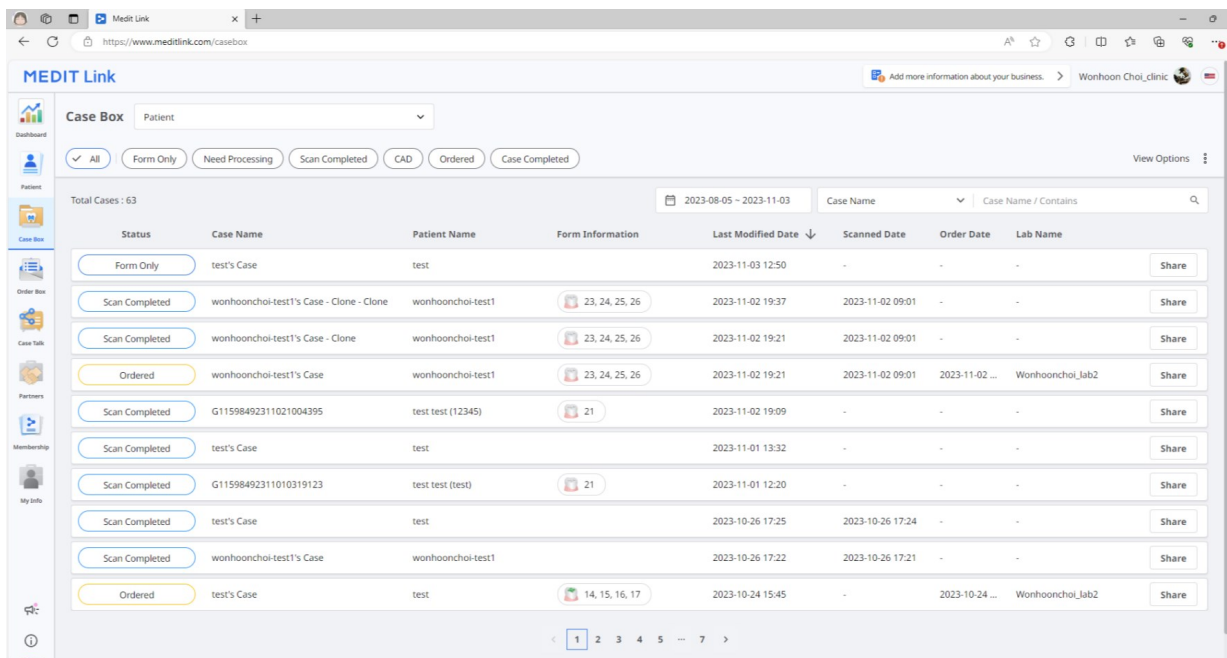
The screenshot shows the Medit Link Case Box interface. At the top, there's a 'Case Box' dropdown menu set to 'Patient'. Below it are filter buttons: 'All', 'Form Only', 'Need Processing', 'Scan Completed', 'CAD', 'Ordered', and 'Case Completed'. A date range '2023-08-05 - 2023-11-03' and a search bar for 'Case Name' are also present. The main table lists cases with columns: Status, Case Name, Patient Name, Form Information, Last Modified Date, Scanned Date, Order Date, and Lab Name. The table contains 10 rows of case data, including details like 'test's Case', 'wonhoonchoi-test1's Case - Clone - Clone', and 'G11598492311021004395'.

Status	Case Name	Patient Name	Form Information	Last Modified Date	Scanned Date	Order Date	Lab Name
Form Only	test's Case	test		2023-11-03 12:50	-	-	-
Scan Completed	wonhoonchoi-test1's Case - Clone - Clone	wonhoonchoi-test1	23, 24, 25, 26	2023-11-02 19:37	2023-11-02 09:01	-	-
Scan Completed	wonhoonchoi-test1's Case - Clone	wonhoonchoi-test1	23, 24, 25, 26	2023-11-02 19:21	2023-11-02 09:01	-	-
Ordered	wonhoonchoi-test1's Case	wonhoonchoi-test1	23, 24, 25, 26	2023-11-02 19:21	2023-11-02 09:01	2023-11-02 ...	Wonhoonchoi_lab2
Scan Completed	G11598492311021004395	test test (12345)	21	2023-11-02 19:09	-	-	-
Scan Completed	test's Case	test		2023-11-01 13:32	-	-	-
Scan Completed	G11598492311010319123	test test (test)	21	2023-11-01 12:20	-	-	-
Scan Completed	test's Case	test		2023-10-26 17:25	2023-10-26 17:24	-	-
Scan Completed	wonhoonchoi-test1's Case	wonhoonchoi-test1		2023-10-26 17:22	2023-10-26 17:21	-	-
Ordered	test's Case	test	14, 15, 16, 17	2023-10-24 15:45	-	2023-10-24 ...	Wonhoonchoi_lab2

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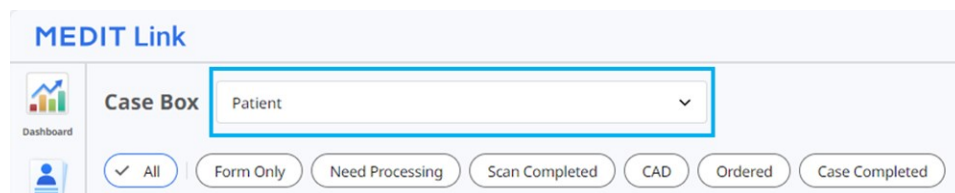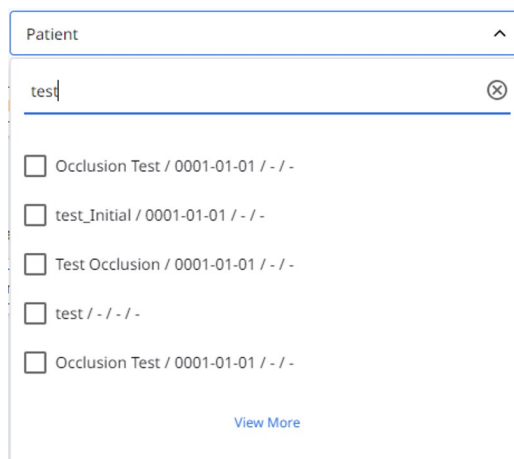
The screenshot shows the search results for 'test' in the Patient dropdown menu. The search bar contains 'test'. Below the search bar, there is a list of search results with checkboxes and a 'View More' link at the bottom.

Search Results
☐ Occlusion Test / 0001-01-01 / - / -
☐ test_initial / 0001-01-01 / - / -
☐ Test Occlusion / 0001-01-01 / - / -
☐ test / - / - / -
☐ Occlusion Test / 0001-01-01 / - / -

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MEDIT Link

Case Box Patient

Dashboard Patient Case Box Order Box Case Table Partners Membership My Info

2023-08-05 ~ 2023-11-03 Case Name Case Name / Contains

Total Cases: 63

Status	Case Name	Patient Name	Form Information	Last Modified Date	Scanned Date	Order Date	Lab Name	Share
Form Only	test's Case	test		2023-11-03 12:50	-	-	-	Share
Scan Completed	wonhoonchoi-test1's Case - Clone	wonhoonchoi-test1	23, 24, 25, 26	2023-11-02 19:37	2023-11-02 09:01	-	-	Share
Scan Completed	wonhoonchoi-test1's Case - Clone	wonhoonchoi-test1	23, 24, 25, 26	2023-11-02 19:21	2023-11-02 09:01	-	-	Share
Ordered	wonhoonchoi-test1's Case	wonhoonchoi-test1	23, 24, 25, 26	2023-11-02 19:21	2023-11-02 09:01	2023-11-02 ...	Wonhoonchoi Lab2	Share
Scan Completed	G11598492311021004395	test test (12345)	21	2023-11-02 19:09	-	-	-	Share
Scan Completed	test's Case	test		2023-11-01 13:32	-	-	-	Share
Scan Completed	G11598492311010319123	test test (test)	21	2023-11-01 12:20	-	-	-	Share
Scan Completed	test's Case	test		2023-10-26 17:25	2023-10-26 17:24	-	-	Share
Scan Completed	wonhoonchoi-test1's Case	wonhoonchoi-test1		2023-10-26 17:22	2023-10-26 17:21	-	-	Share
Ordered	test's Case	test	14, 15, 16, 17	2023-10-24 15:45	-	2023-10-24 ...	Wonhoonchoi Lab2	Share

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All ☒ Form Only ☐ Need Processing ☒ Scan Completed ☐ CAD ☐ Ordered ☐ Case Completed

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Case Box(診所)	All / Form Only / Need Processing / Scan Completed / CAD / Ordered / Case Completed
--------------	---

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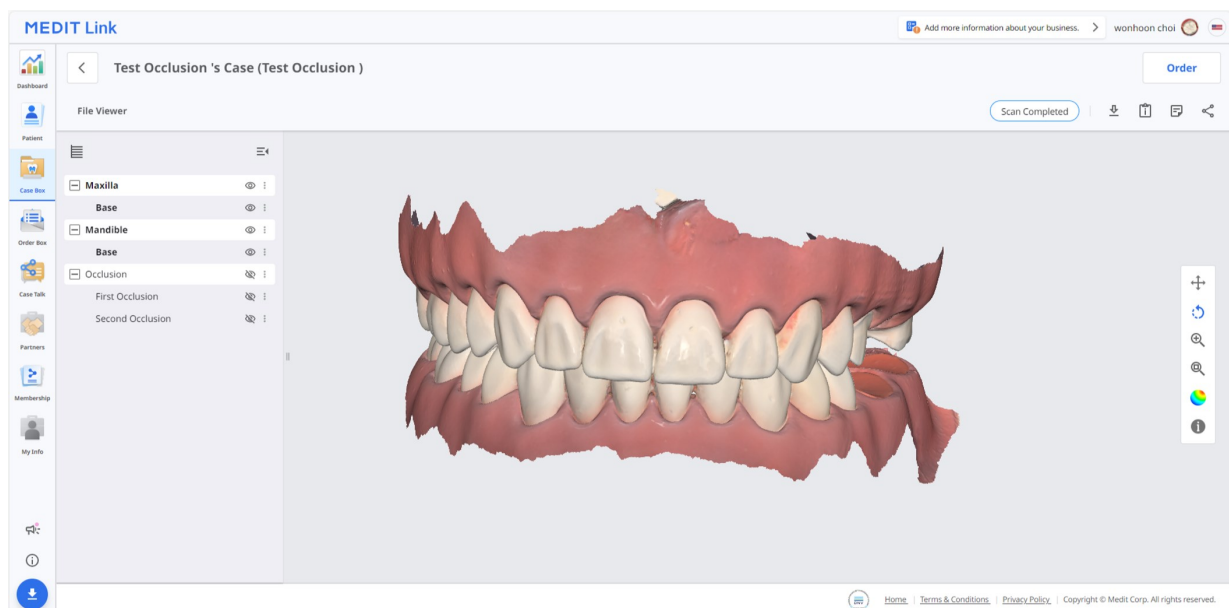
The screenshot displays the Medit Link application interface for viewing order details. The interface includes a sidebar with navigation icons, a top header with the user's name 'wonhoon choi', and a main content area. The main content area is divided into two sections: 'Order' and 'Order List'. The 'Order' section contains fields for Order ID, Patient Name, Occlusion Test, Scanned Date, Order Date, Requested Delivery Date, Expected Delivery Date, AM 08:15, and Memo. The 'Order List' section contains a table with columns: No., Type, Method, Material, and Shade. The table lists four orders (23, 24, 25, 26) all of type 'Inlay' with material 'Multi-Layered ZL...' and shade 'None'. The interface is annotated with numbered circles 1 through 6, corresponding to the table below.

No.	Type	Method	Material	Shade
23	Inlay	-	Multi-Layered ZL...	None
24	Inlay	-	Multi-Layered ZL...	None
25	Inlay	-	Multi-Layered ZL...	None
26	Inlay	-	Multi-Layered ZL...	None

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Standard

Free

- ✓ Free plan for every member
- ✓ 1 GB of cloud storage
- ✓ Efficient for new users

If you have any difficulties with payment, select "Click" on the right. [Click](#)

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☐ Total Selected 0.00 KB [Delete](#)

☐ Raw Data 147.52 GB

Remove raw data.
Raw data does not take up space in cloud storage.

☐ Result Files & Attachments 13.74 GB

Remove result files and attachments such as mesh, image, or video files.

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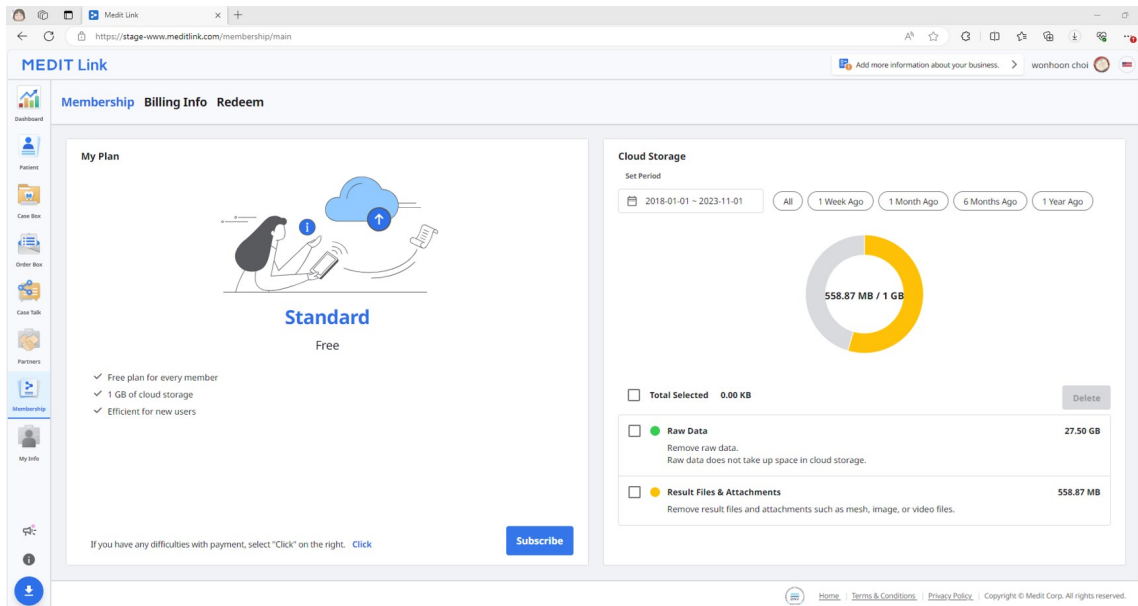
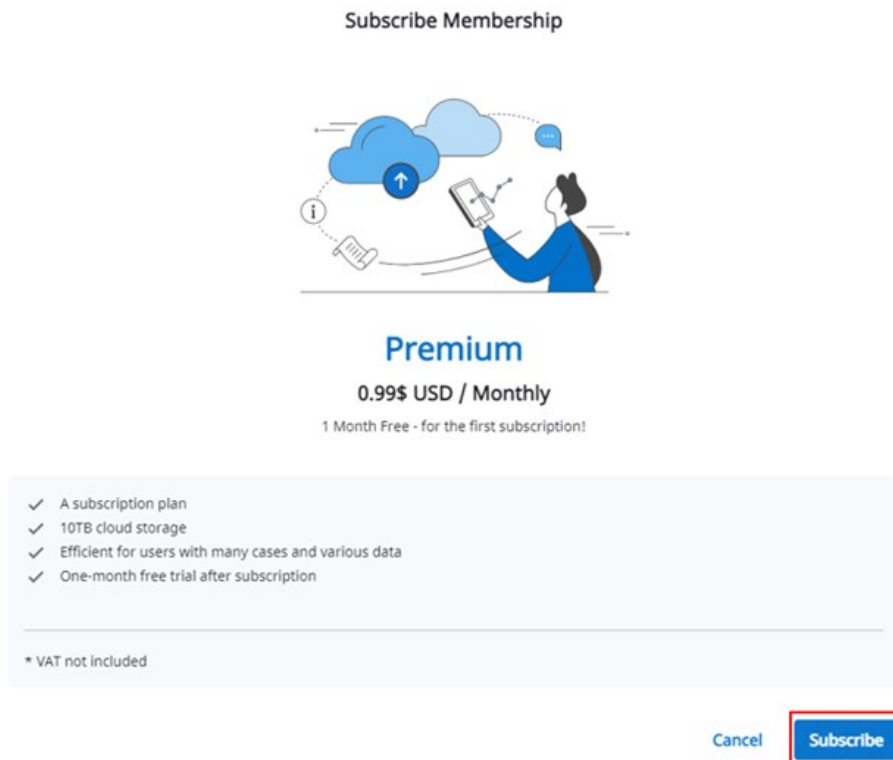
← MEDIT Link TEST MODE

Subscribe to Premium membership for Clinic

\$0.00

Then \$0.99 per day

Medit Link Membership

Premium membership for Clinic	\$0.99
Billed daily	
Subtotal	\$0.99
1 MONTH FREE – for first subscription 100% off for a day	
Total due today	\$0.00

Powered by stripe | Terms Privacy

Pay with card

Email

Card information

1234 1234 1234 1234

VISA

MM / YY

CVC

Name on card

Country or region


South Korea

☐ Save information to pay faster next time

Subscribe

By confirming your subscription, you allow MEDIT SG PTE. LTD. to charge your card for this payment and future payments in accordance with their terms.

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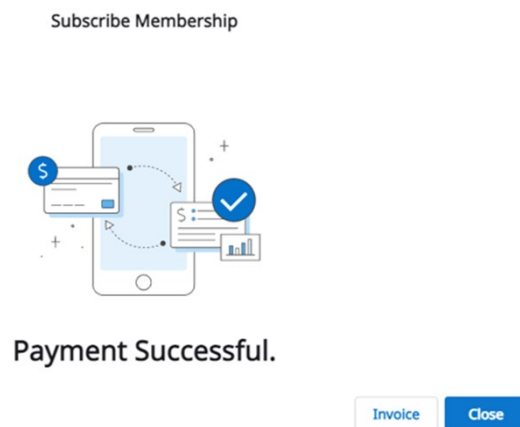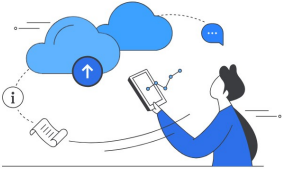
Premium

Free

The first month is free if you join the membership.

- ✓ A paid plan
- ✓ 10 TB of cloud storage
- ✓ Efficient for users with many cases and various data
- ✓ One-month free trial included

The following payment date is 2023-11-30. [View More](#)
If you have any difficulties with payment, select "Click" on the right. [Click](#)

[Unsubscribe](#)

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Cancel Subscription

Premium

In Use

Standard

Free

Are you sure you want to unsubscribe?

- ✓ A free plan for every member
- ✓ 1GB cloud storage
- ✓ Efficient for new users

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Confirm Unsubscribe

Your subscription will be canceled at the end of your billing period on 2021-09-10

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Confirm Unsubscribe

Your subscription has been canceled.

Your current plan can use up to 2021-09-10. After that, the cloud capacity is limited to 1GB.

If the cloud capacity is exceeded, the files will not be backed up, and the order and share functions will not be available.

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Membership **Billing Info** Redeem

My Plan and Billing Information
Plan Premium
Amount 0.99\$ USD / Monthly
Next Payment Date 2021-09-09
Payment Method Visa *****4242 (11/21) Stripe [Edit Payment Method](#)
If you haven't subscribed yet?
Experience a differentiated premium.

Billing History

Payment Date	Membership	From	To	Payment Amount	Status
2021-09-09	Premium	2021-09-09	2021-09-10	0.99 \$ (USD)	Payment Completed
2021-09-09	Premium	2021-09-09	2021-09-09	0 \$ (USD)	Payment Completed

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Membership **Billing Info** Redeem

My Plan and Billing Information
Plan Premium
Amount 0.99\$ USD / Monthly
Next Payment Date 2021-09-09
Payment Method Visa *****4242 (11/21) Stripe [Edit Payment Method](#)
If you haven't subscribed yet?
Experience a differentiated premium.

Billing History

Payment Date	Membership	From	To	Payment Amount	Status
2021-09-09	Premium	2021-09-09	2021-09-10	0.99 \$ (USD)	Payment Completed
2021-09-09	Premium	2021-09-09	2021-09-09	0 \$ (USD)	Payment Completed

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Add

Go back

By adding a payment method, you agree to MEDIT SG PTE. LTD.'s Terms of Service and Privacy Policy.

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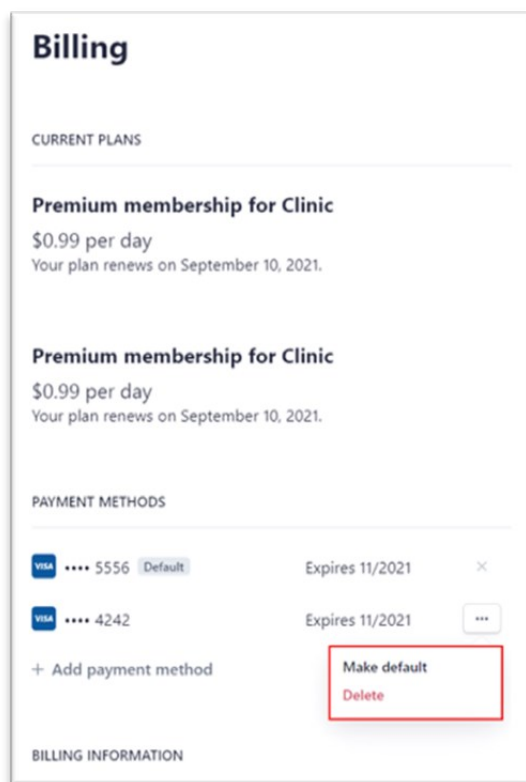
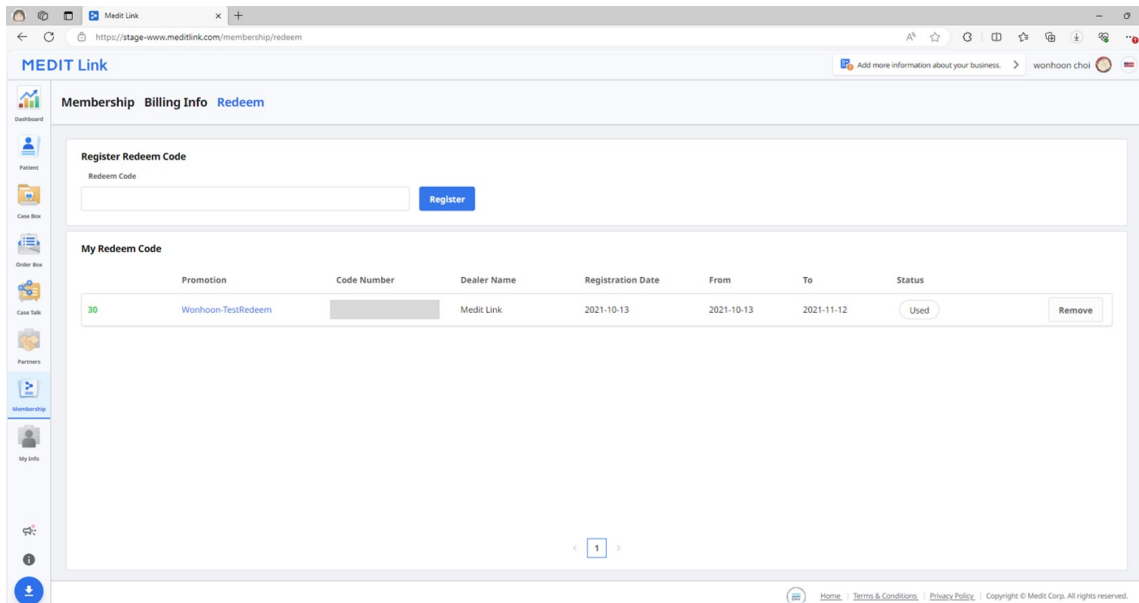
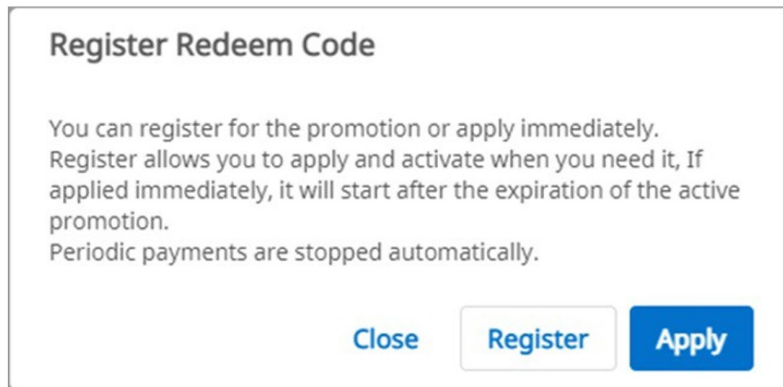
Country

Korea, Republic of

Your chosen language is applied to the email sent by Medit Link and the notification.

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Edit

Cover and Logo

Looking for digitalized dental laboratories?

Find the most ordered Labs around your clinic with 'Featured Part

123

Services

Select

Specialties (choose up to 3)

Select

ⓘ Select the services your business specializes in from the list.

About Us

e.g., Welcome to OOO Dental Clinic.
We are a licensed dental clinic situated in California specializing mainly in aesthetic restorations and implantology.

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-

Buy

ⓘ The scanner information is automatically registered once your scanner is connected to a Medit scan program.

Home

Terms & Conditions

Privacy Policy

Copyright © Medit Corp. All rights reserved.

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Notice
NFT Service Termination Notice
2023-04-19

To our users,

We want to thank all Medit Link users who have supported and used our NFT service.

Unfortunately, we must inform you that this service will be terminated following the upcoming program update.

Please see the details below for more information:

- **Feature in Question**
 - Share via NFT (in Medit Link Case Talk)
- **Effective Date**
 - Starting from the update in the first half of 2023 (scheduled for April 27)

※ The update schedule is subject to change.

We'll continuously work to provide better service.

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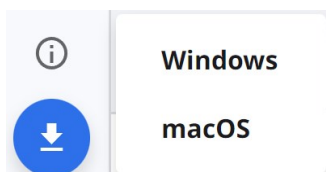
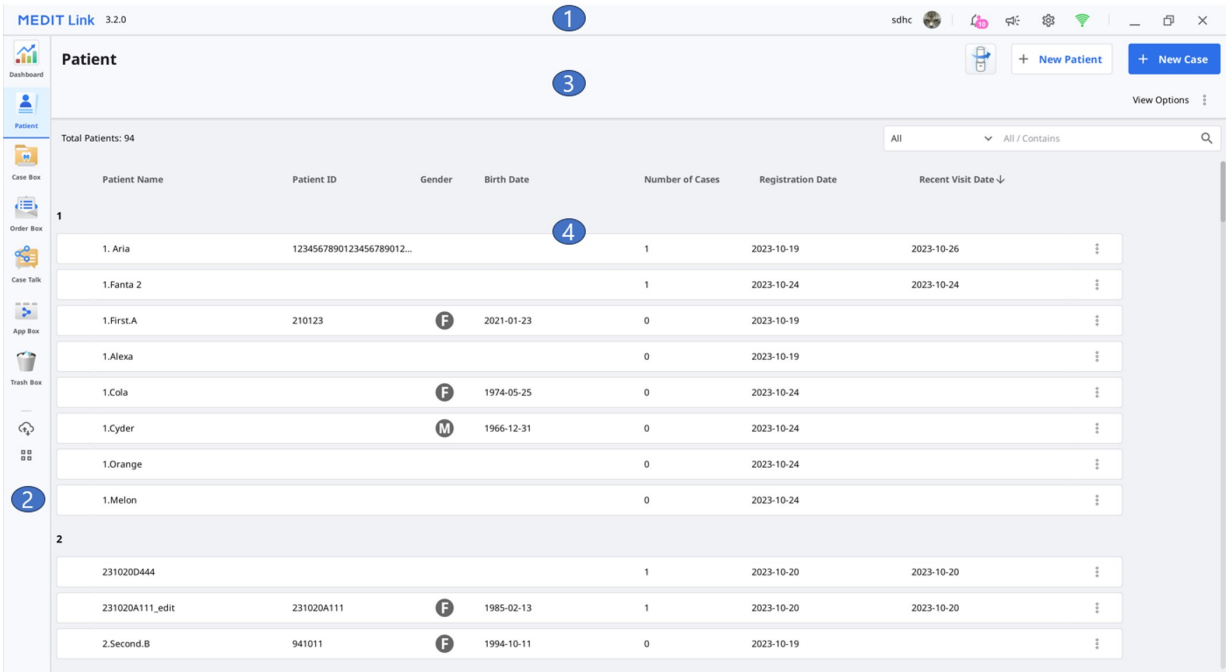
Notification

1☐ Display Unread Messages Only [Mark All As Read](#)

2

Delivery 2020-08-20 오후 5:20 
Wonhoonchoi_Lab_Dev has dispatched the ordered package.

Order 2020-08-20 오후 5:18
Wonhoonchoi_Lab_Dev has accepted the order.

Order 2020-08-17 오후 4:05
Wonhoonchoi_Lab_Dev has accepted the order.

Order 2020-08-17 오후 4:05 
Wonhoonchoi_Lab_Dev has accepted the order.

Order 2020-08-17 오후 4:05 
Wonhoonchoi_Lab_Dev has accepted the order.

Order 2020-08-17 오후 4:05 
Wonhoonchoi_Lab_Dev has accepted the order.

Work 2020-08-14 오후 6:33 
Wonhoonchoi_Clinice_Dev_Member registered a new case test's Case.

Work 2020-08-14 오전 10:15 
Wonhoonchoi_Clinice_Dev_Member registered a new case test's Case.

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The screenshot shows the 'Settings' window with the 'General' tab selected. The left sidebar contains links for General, Storage Management, Set Connected Program, Export Options, Case Converting Tool, and Proxy Option. The main area displays various settings: 'Save Folder' is set to 'C:/ProgramData/Medit/Medit Link/Data/'; 'Language Settings' is set to 'English'; 'Teeth Labeling' has 'FDI' selected; 'Notification' has 'Allow Notifications' checked; 'Upload Raw Data' has 'Upload Automatically' checked; 'Quick Scan' has 'Pre-Op for Maxilla, Pre-Op for Mandible, Face Scan' selected; 'Viewing Options' has 'Advanced Rendering' checked and 'Adjust Color Texture' is set to a medium level. At the bottom, there is an 'Apply' button and a 'Close' link.

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Question

Exporting successfully completed.
Do you want to open the folder containing the exported file?

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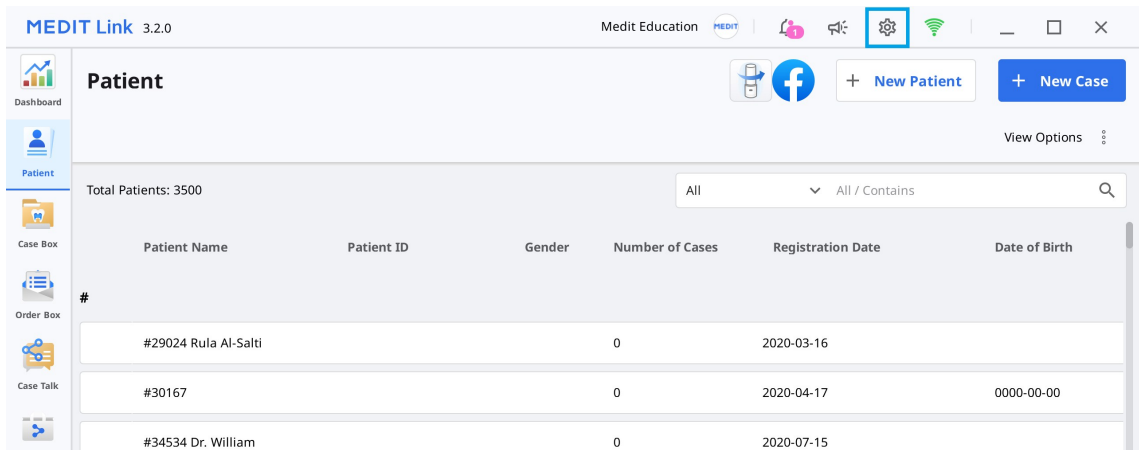
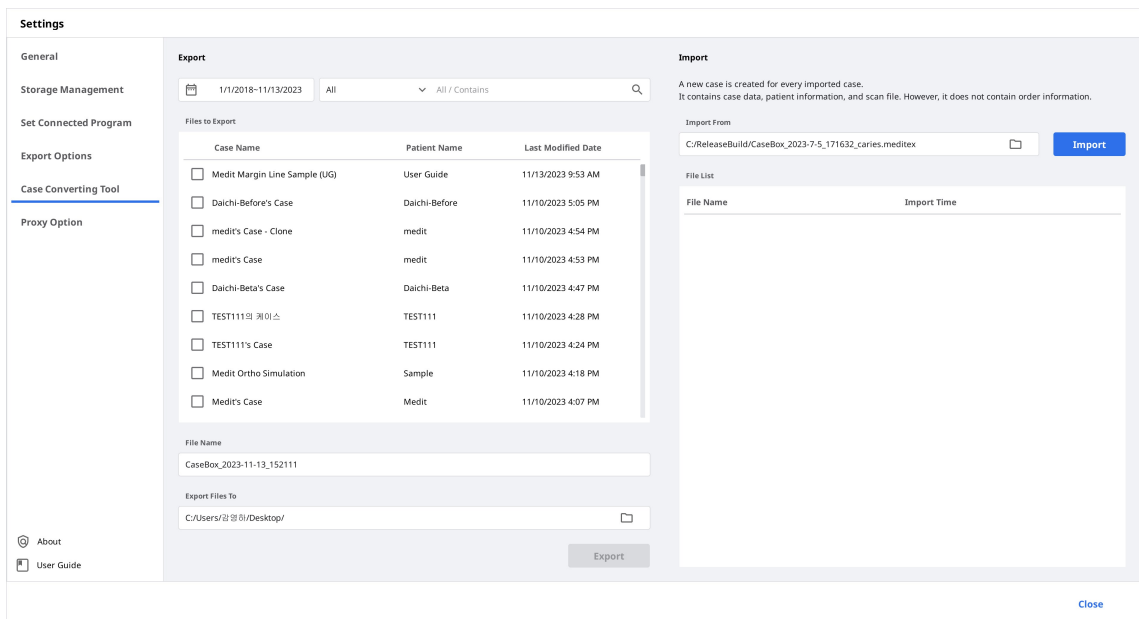
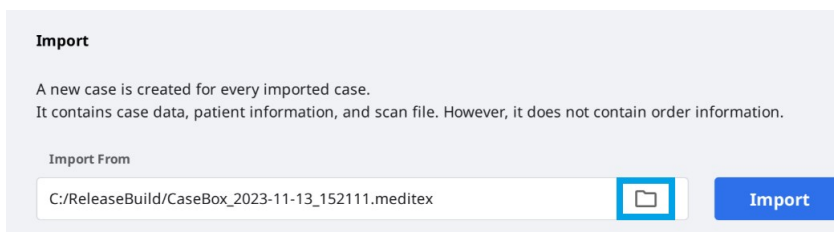
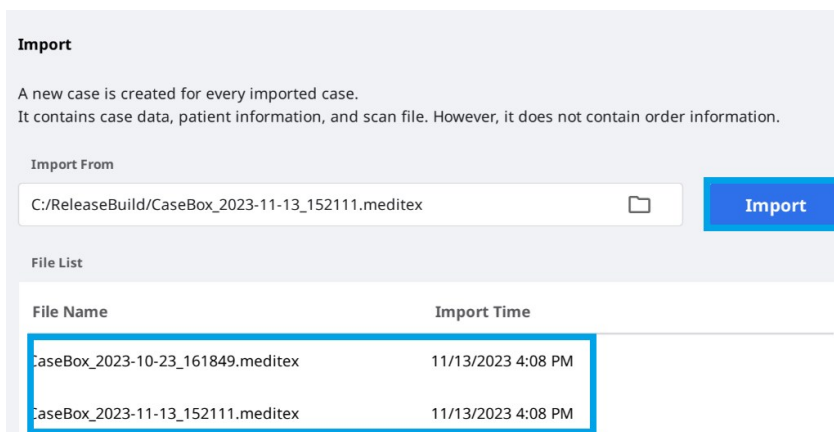
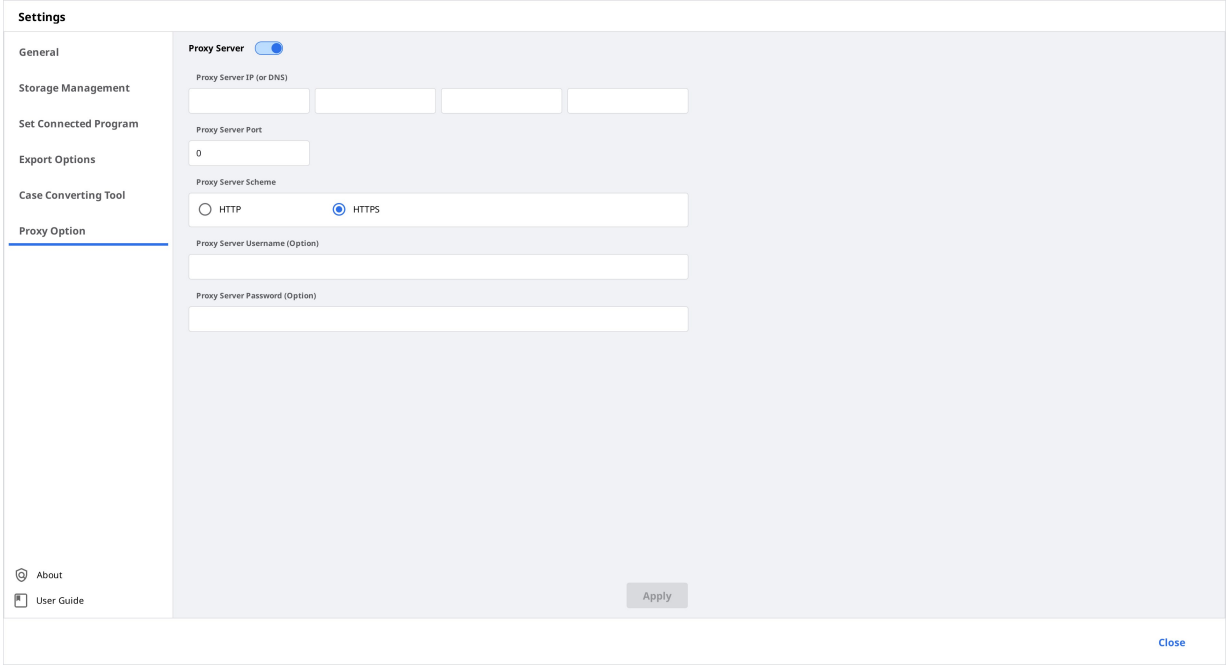
The screenshot shows the 'Settings' window with the 'Proxy Option' tab selected. The 'Proxy Server' toggle is turned on. Below it, there are input fields for 'Proxy Server IP (or DNS)', 'Proxy Server Port' (set to 0), and 'Proxy Server Scheme' (set to HTTPS). There are also optional fields for 'Proxy Server Username' and 'Proxy Server Password'. An 'Apply' button is at the bottom right of the settings area. The left sidebar contains links to 'About' and 'User Guide'.

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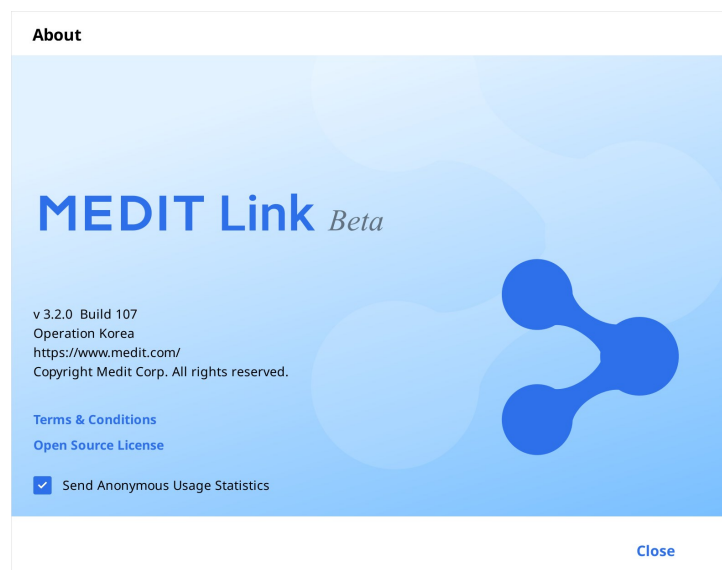
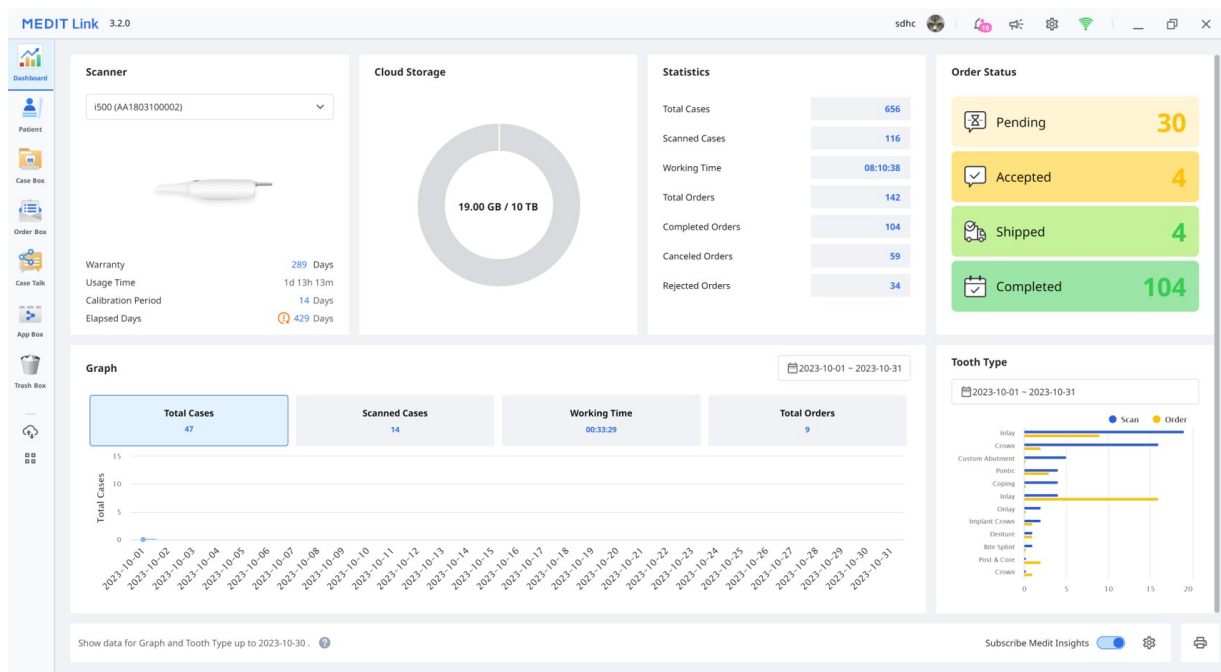
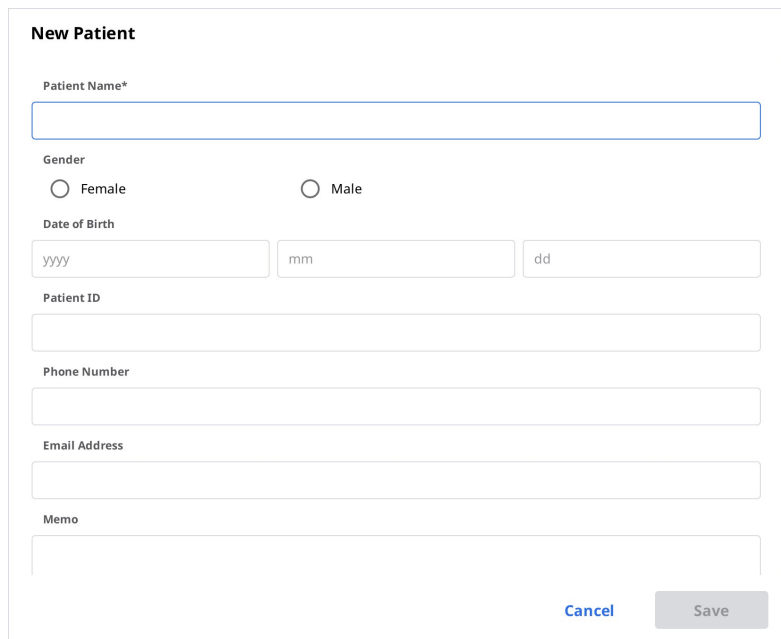
The 'New Patient' form contains the following fields and options:

- Patient Name* (text input)
- Gender: ☐ Female, ☐ Male
- Date of Birth: Three separate input fields for yy, yy, mm, dd.
- Patient ID (text input)
- Phone Number (text input)
- Email Address (text input)
- Memo (text area)
- Buttons: Cancel, Save

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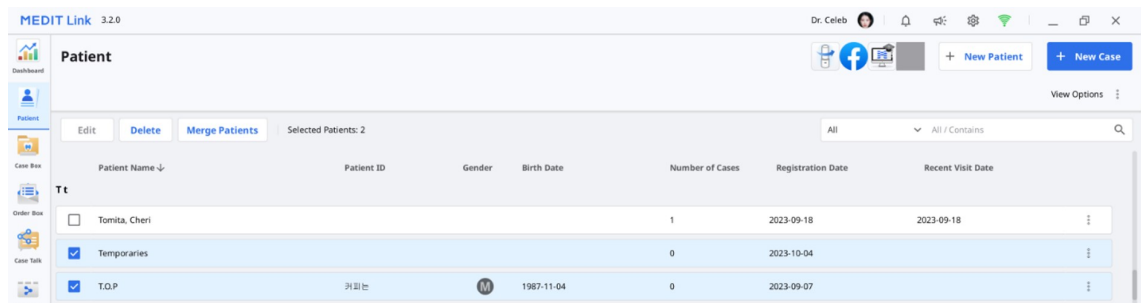
The screenshot shows the 'Patient' management interface in MEDIT Link 3.2.0. It includes a table with the following columns: Patient Name, Patient ID, Gender, Birth Date, Number of Cases, Registration Date, and Recent Visit Date. The table lists three patients: Tomita, Cheri; Temporaries; and T.O.P. The 'Temporaries' and 'T.O.P.' rows are selected with checkboxes.


Patient Name	Patient ID	Gender	Birth Date	Number of Cases	Registration Date	Recent Visit Date
Tomita, Cheri				1	2023-09-18	2023-09-18
Temporaries				0	2023-10-04	
T.O.P.	카미타	♂	1987-11-04	0	2023-09-07	

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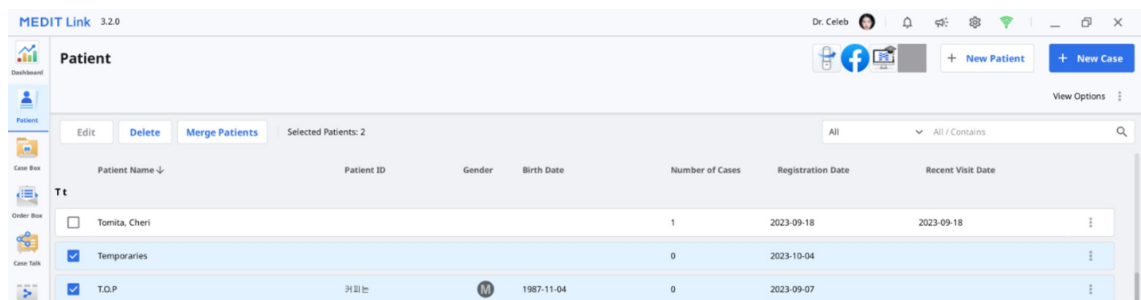This is a duplicate of the screenshot above, showing the 'Patient' management interface with the same table of patients.

Patient Name	Patient ID	Gender	Birth Date	Number of Cases	Registration Date	Recent Visit Date
Tomita, Cheri				1	2023-09-18	2023-09-18
Temporaries				0	2023-10-04	
T.O.P.	카미타	♂	1987-11-04	0	2023-09-07	

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Merge all patient information to

Please select the patient you'd like to merge all data to.
After merging all data together, all information in other patient files will be deleted.

	Patient Name	Patient ID	Gender	Birth Date	Number of Cases	Date Registered	Recent Visit Date
☒	5-9_iScan_CAMOrd...				1	2020-09-18	2020-09-18
☐	test				2	2020-09-17	2020-09-17
☐	test				22	2019-07-26	2020-09-17

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
Order Box/In Box

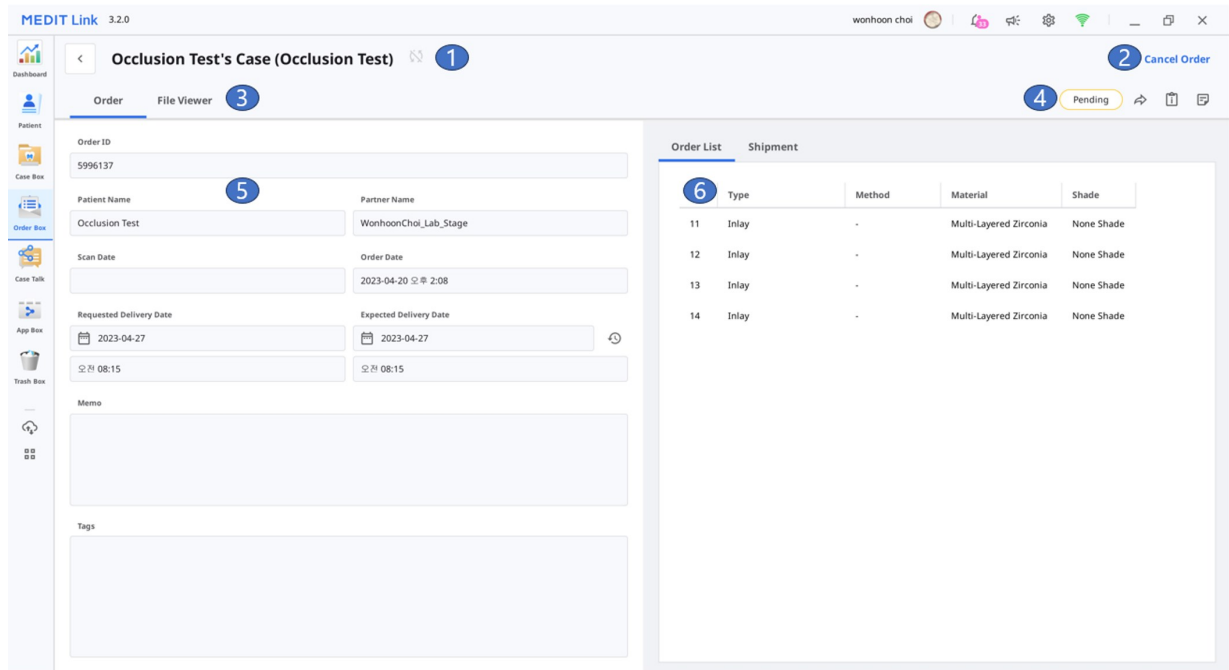The screenshot shows the 'Occlusion Test's Case (Occlusion Test)' screen in the Medit Link 3.2.0 application. The interface is divided into several sections:

- Top Bar:** Includes a 'Cancel Order' button (labeled 2) and a 'Pending' status indicator (labeled 4).
- Left Sidebar:** Contains navigation icons for Dashboard, Patient, Case Box, Order Box, Case Talk, App Box, and Trash Box.
- Main Content Area:**
 - Order ID:** 5996137 (labeled 5).
 - Patient Name:** Occlusion Test (labeled 5).
 - Partner Name:** WonhoonChoi_Lab_Stage.
 - Scan Date:** 2023-04-20 오후 2:08.
 - Order Date:** 2023-04-20 오후 2:08.
 - Requested Delivery Date:** 2023-04-27.
 - Expected Delivery Date:** 2023-04-27.
 - Memo:** A text area for notes.
 - Tags:** A list of tags.
- Order List Table:**

Type	Method	Material	Shade
11 Inlay	-	Multi-Layered Zirconia	None Shade
12 Inlay	-	Multi-Layered Zirconia	None Shade
13 Inlay	-	Multi-Layered Zirconia	None Shade
14 Inlay	-	Multi-Layered Zirconia	None Shade

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Case Talk Name	Patient Name	Form Information	Scanned Date	Shared Date	Expiration Date
BS_4_Scan for Clinics's Case - Clone_2023-10-31 15:04:08	B5***** fo* C*****	Maxilla, Mandible	2020-09-25 14:12	2023-10-31 15:04	2024-01-29 23:59
BS_4_Scan for Clinics's Case - Clone_2023-10-20 11:59:15	BS_4_Scan for Clinics	Maxilla, Mandible	2020-09-25 14:12	2023-10-20 11:59	2023-11-20 23:59
1111	Email_Push	24, 25	2022-09-30 13:48	2022-09-30 16:35	unlimited

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App Name	Category	Version	Status
Medit ClinicCAD	CAD	v 0.9.2.237	Install
Medit Caries Detection	Tool	v 0.9.0.47	Install
Medit Occlusion Analyzer	Analysis Tool	v 1.0.1.72	Installed
Medit Margin Lines	Tool CAD	v 1.0.0.23	Installed
MIPACS Integration	PACS/PMS	v 0.9.0.0	Install
Medit Splints	Tool	v 1.1.0.9	Update
Medit Design (Formerly Known as Medit Compare)	Tool	v 2.1.3.95	Install
Medit Model Builder	Tool CAD	v 1.3.2.65	Install
Medit Ortho Simulation	Tool	v 1.3.0.22	Install

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
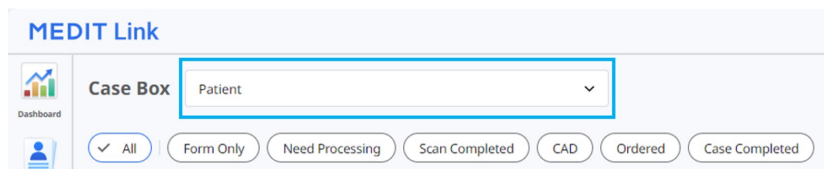
Case Box Patient

Dashboard

✓ All Form Only Need Processing Scan Completed CAD Ordered Case Completed

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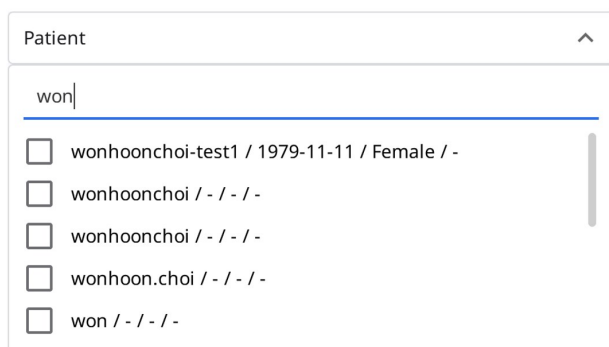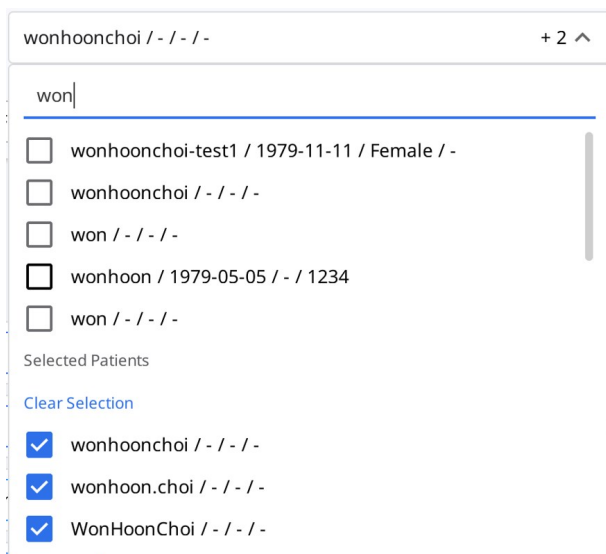
Status	Patient Name	Case Name	Form Information	Last Modified Date	Scan Date	Lab Name	Order Date
Form Only	wonhoonchoi	originalcase	-	2023-10-05 02:10			
Form Only	wonhoonchoi	wonhoonchoi's Case	-	2023-10-05 02:09			
Scan Completed	wonhoonchoi	FileViewerTest_1	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-06-15 02:14	2023-02-09 02:35		
Scan Completed	wonhoonchoi	FileViewerTest_1 - Clone	2-Inlay / 3-Inlay / 4-Inlay / 5-Inlay / 6-Inlay / 8...	2023-05-25 02:05	2023-02-09 02:35		
Form Only	wonhoonchoi	FileViewerTest_1 - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-03-15 02:06			
Scan Completed	wonhoonchoi	FileViewerTest_1 - Clone - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-03-09 02:41	2023-02-09 02:35		
Need Processing	wonhoonchoi	FileViewerTest_1 - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-02-14 02:13	2023-02-09 02:35		
Scan Completed	wonhoonchoi	FileViewerTest_1 - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-02-10 02:11	2023-02-09 02:35		
Form Only	wonhoonchoi	testsamplecase	-	2023-02-10 02:39			
Scan Completed	wonhoonchoi	Multi Occlusion	-	2023-02-09 02:26	2023-02-09 02:01		
Scan Completed	wonhoonchoi	FileViewerTest_1 - Clone - Clone	8-Inlay / 10-Inlay / 12-Inlay / 13-Inlay / 14-Inl...	2023-02-09 02:50	2023-02-09 02:35		
Case Completed	wonhoonchoi	wonhoonchoi's Case	12-Inlay / 13-Inlay / 14-Inlay	2022-10-13 02:30	2022-10-11 02:56		
Case Completed	wonhoonchoi	wonhoonchoi's Case	-	2022-10-13 02:10	2022-10-13 02:10		

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
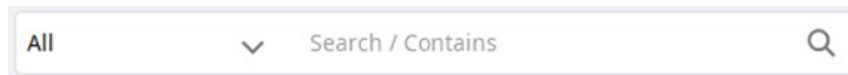
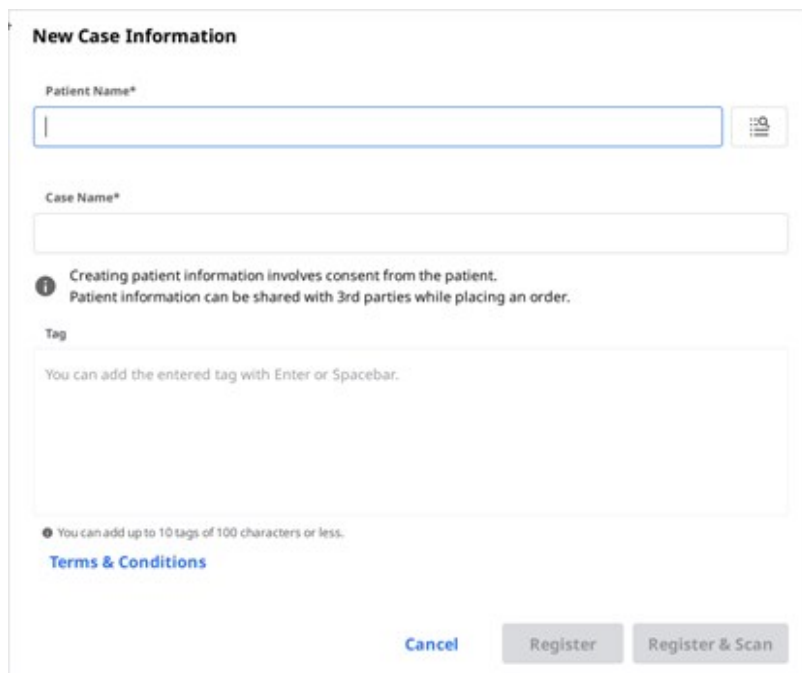
The form titled "New Case Information" contains the following fields and elements:

- Patient Name***: A text input field with a dropdown icon on the right.
- Case Name***: A text input field.
- Information**: A note stating "Creating patient information involves consent from the patient. Patient information can be shared with 3rd parties while placing an order."
- Tag**: A text area with the instruction "You can add the entered tag with Enter or Spacebar."
- Footer**: Includes a note "You can add up to 10 tags of 100 characters or less.", a link for "Terms & Conditions", and three buttons: "Cancel", "Register", and "Register & Scan".

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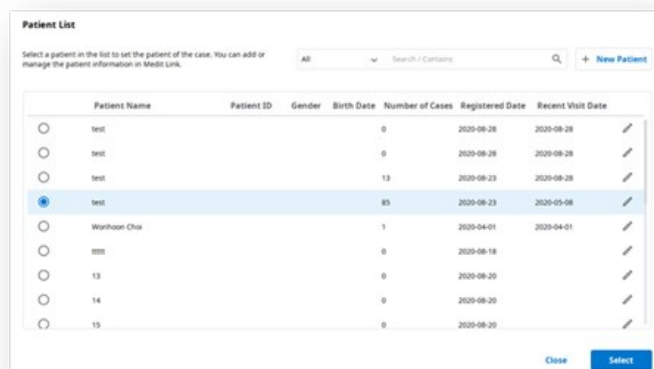
The "Patient List" dialog box displays a table of existing patients. At the top, there is a search bar and a "New Patient" button. The table has the following columns: Patient Name, Patient ID, Gender, Birth Date, Number of Cases, Registered Date, and Recent Visit Date. The first four rows have "test" as the patient name. The fifth row is "Wenhsien Chai". The sixth row is "1000". The seventh row is "13". The eighth row is "14". The ninth row is "15". The "Number of Cases" column shows values 0, 0, 13, 85, 1, 0, 0, 0, 0. The "Registered Date" column shows values 2020-08-28, 2020-08-28, 2020-08-23, 2020-08-23, 2020-04-01, 2020-08-18, 2020-08-20, 2020-08-20, 2020-08-20. The "Recent Visit Date" column shows values 2020-08-28, 2020-08-28, 2020-08-28, 2020-05-08, 2020-04-01, and empty for the last three rows. Each row has a radio button in the "Patient Name" column and an edit icon in the "Recent Visit Date" column. The fourth row (with 85 cases) is selected. At the bottom, there are "Close" and "Select" buttons.

Patient Name	Patient ID	Gender	Birth Date	Number of Cases	Registered Date	Recent Visit Date
☐ test				0	2020-08-28	2020-08-28
☐ test				0	2020-08-28	2020-08-28
☐ test				13	2020-08-23	2020-08-28
☒ test				85	2020-08-23	2020-05-08
☐ Wenhsien Chai				1	2020-04-01	2020-04-01
☐ 1000				0	2020-08-18	
☐ 13				0	2020-08-20	
☐ 14				0	2020-08-20	
☐ 15				0	2020-08-20	

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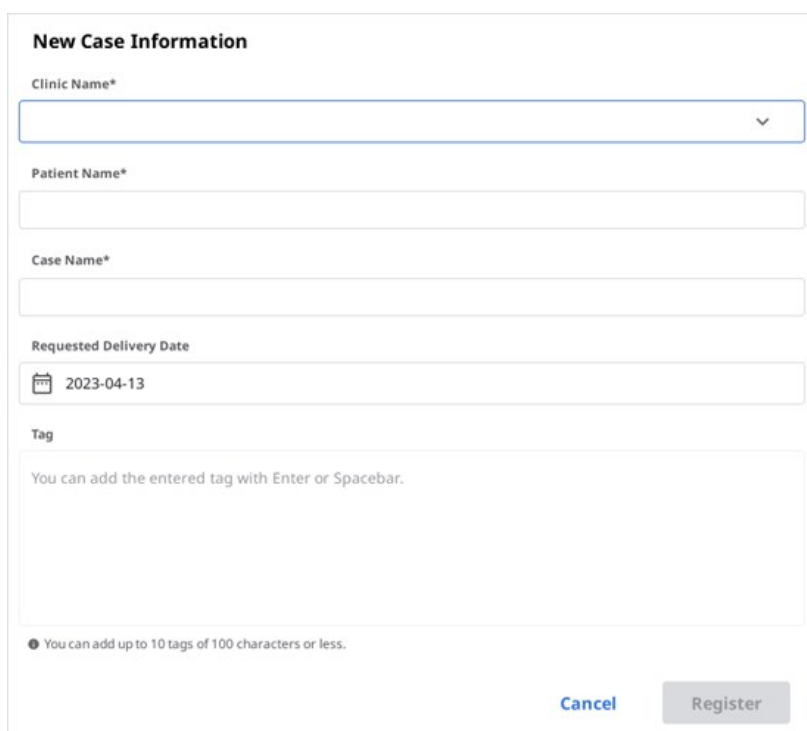The form is titled "New Case Information". It contains the following fields:

- Clinic Name***: A dropdown menu.
- Patient Name***: A text input field.
- Case Name***: A text input field.
- Requested Delivery Date**: A date picker showing "2023-04-13".
- Tag**: A text area with the instruction "You can add the entered tag with Enter or Spacebar." and a note "You can add up to 10 tags of 100 characters or less."


At the bottom right, there are two buttons: "Cancel" and "Register".

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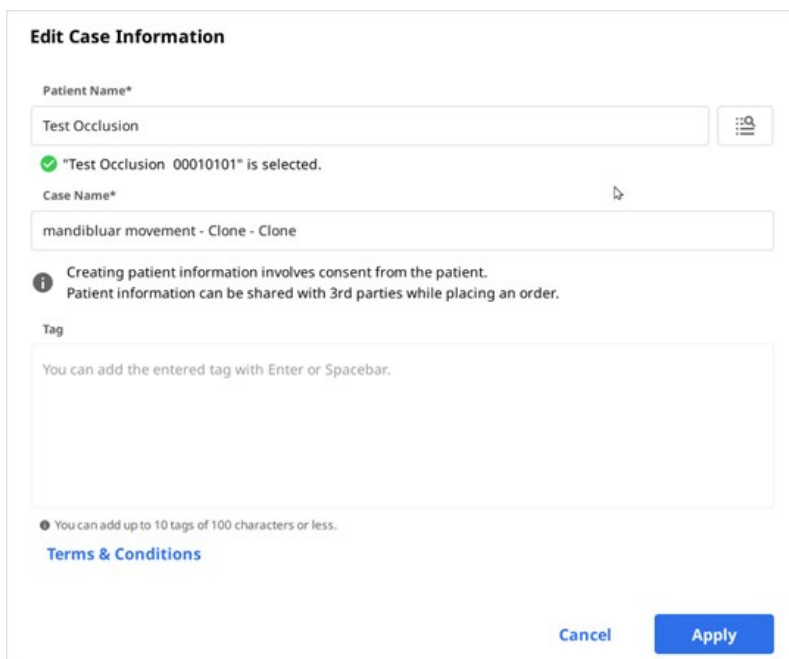The form is titled "Edit Case Information". It contains the following fields:

- Patient Name***: A text input field containing "Test Occlusion" and a search icon.
- Case Name***: A text input field containing "mandibular movement - Clone - Clone".
- Tag**: A text area with the instruction "You can add the entered tag with Enter or Spacebar." and a note "You can add up to 10 tags of 100 characters or less."

Below the Case Name field, there is an information icon and a message: "Creating patient information involves consent from the patient. Patient information can be shared with 3rd parties while placing an order."

At the bottom right, there are two buttons: "Cancel" and "Apply".

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Patient List

Select a patient in the list to set the patient of the case. You can add or manage the patient information in Medit Link.

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Options File List

Export Files To
C:/Users/Medit-Wonhoonchoi/Desktop/test/

Folder Name
Multi Occlusion - Scan - Before Scanning

File Name
Multi Occlusion - Scan - Before Scanning

CAD Type
Maestro

Mesh Format
Convert mesh files to the below formats
☐ meditMesh ☐ OBJ ☒ PLY ☐ STL
☒ Do not convert the attached files

Image Format
☐ Convert images to the below format
☐ JPG ☐ JPEG ☒ PNG ☐ BMP


☐ Compress Files
☐ Combine Individual Mesh
☐ Export All Occlusal Relationship
☐ Save Case Information as PDF [Preview](#)

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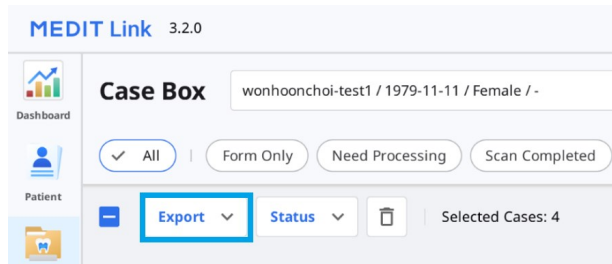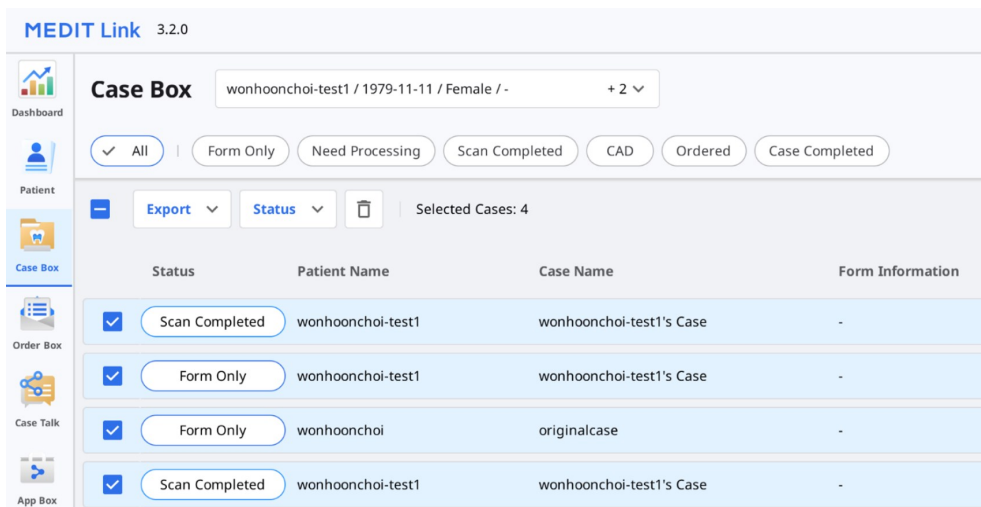
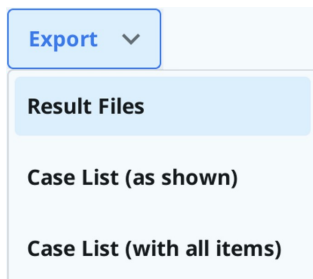
Options

Export Files To
C:/Users/Wonhoon Choi Laptop/Desktop/ExportTest/3/

Folder Name
Case Creation Date - Case Name

File Name
Case Creation Date - Case Name

Type
Medit

Mesh Format
Convert mesh files to the below formats
☐ meditMesh ☐ OBJ ☒ PLY ☐ STL
☐ Do not convert the attached files

Image Format
Convert images to the below format
☐ JPG ☐ JPEG ☒ PNG ☐ BMP

☐ Compress Files
☐ Combine Individual Mesh
☐ Export All Occlusal Relationship
☐ Save Case Information as PDF

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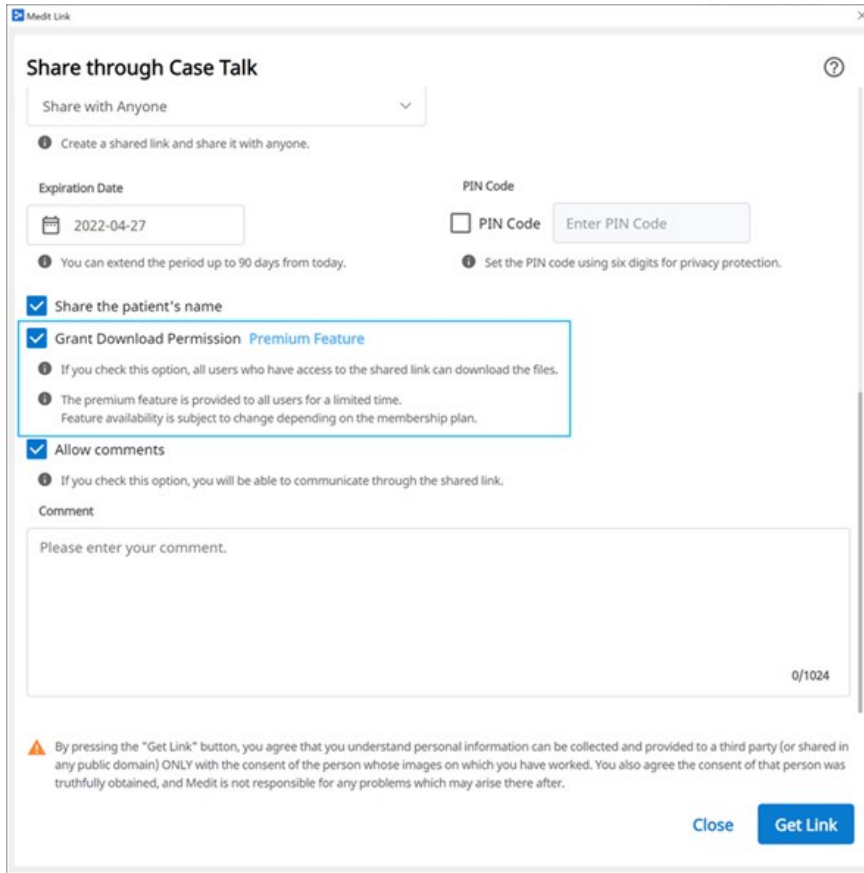
Share through Case Talk

Share with Anyone

Create a shared link and share it with anyone.

Expiration Date: 2022-04-27

PIN Code: ☐ PIN Code Enter PIN Code

You can extend the period up to 90 days from today. Set the PIN code using six digits for privacy protection.

☒ Share the patient's name

☒ Grant Download Permission **Premium Feature**

If you check this option, all users who have access to the shared link can download the files.

The premium feature is provided to all users for a limited time. Feature availability is subject to change depending on the membership plan.

☒ Allow comments


If you check this option, you will be able to communicate through the shared link.

Comment: Please enter your comment. 0/1024

By pressing the "Get Link" button, you agree that you understand personal information can be collected and provided to a third party (or shared in any public domain) ONLY with the consent of the person whose images on which you have worked. You also agree the consent of that person was truthfully obtained, and Medit is not responsible for any problems which may arise there after.

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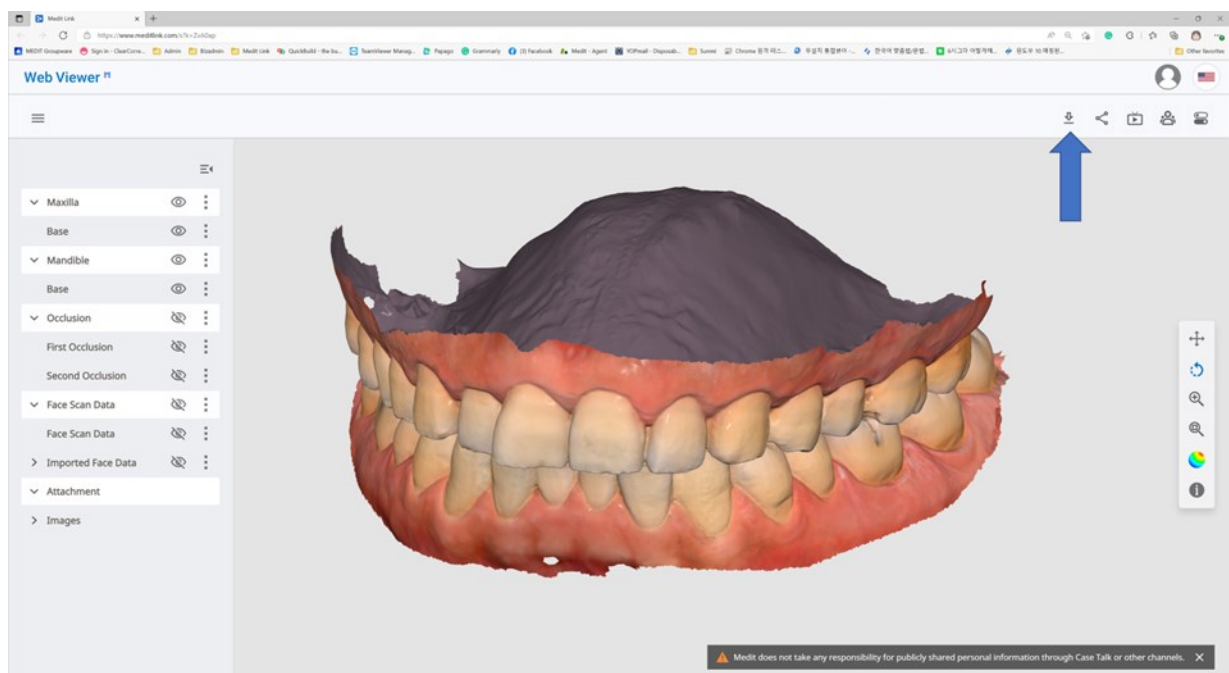Download Premium Feature

The premium feature is provided to all users for a limited time.
Feature availability is subject to change depending on the membership plan.

Mesh Format

Convert mesh files to the below formats

☐ meditMesh ☐ OBJ ☒ PLY ☐ STL

i The selected format will be applied to the result files created in the Medit Link application.

i The attachment files will not be converted but will keep their original format.

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Case Name: mandibular movement - Clone - Clone

Patient Name: Test Occlusion

☒ Share patient information

Requested Delivery Date: 2023-04-13

Delivery Time: AM 8:00

Memo

Tags

You can add the entered tag with Enter or Spacebar.

● You can add up to 10 tags of 100 characters or less.

☒ Include all occlusion relationships in the order

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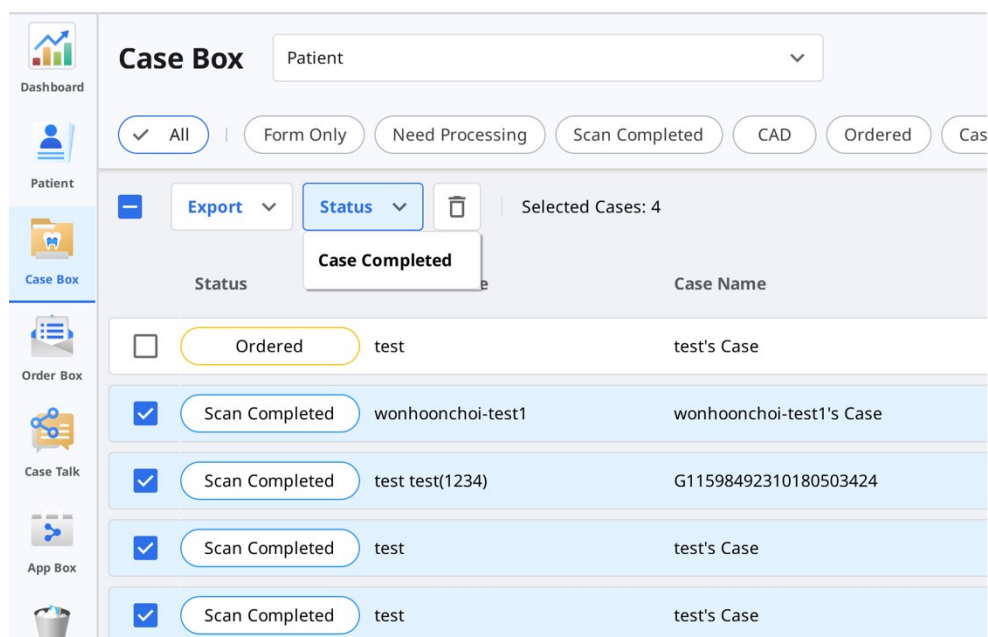
The screenshot displays the 'Case Box' interface. On the left is a sidebar with navigation icons for Dashboard, Patient, Case Box (highlighted), Order Box, Case Talk, and App Box. The main area is titled 'Case Box' and includes a 'Patient' dropdown menu. Below this are filter buttons: 'All' (selected), 'Form Only', 'Need Processing', 'Scan Completed', 'CAD', 'Ordered', and 'Cas'. A toolbar contains an 'Export' button, a 'Status' dropdown menu (which is open, showing 'Case Completed' as an option), and a trash icon. To the right of the toolbar, it says 'Selected Cases: 4'. Below the toolbar is a table with columns 'Status' and 'Case Name'. The table contains five rows of case data.

Status	Case Name
☐ Ordered	test
☒ Scan Completed	wonhoonchoi-test1
☒ Scan Completed	test test(1234)
☒ Scan Completed	test
☒ Scan Completed	test

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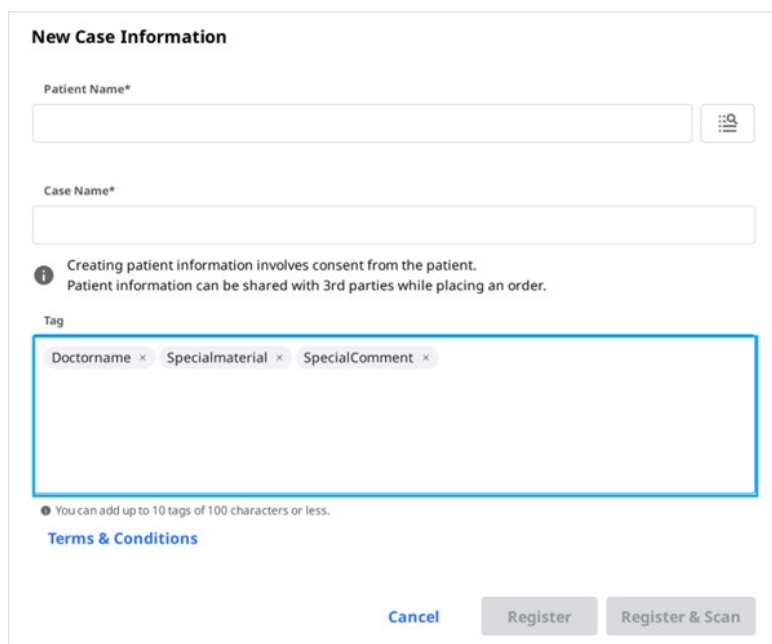
Patient Name*

Case Name*

*Creating patient information involves consent from the patient.
Patient information can be shared with 3rd parties while placing an order.*

Tag

Doctorname × Specialmaterial × SpecialComment ×

You can add up to 10 tags of 100 characters or less.

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